

GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

Data Sheet



# Table of Contents

System.....3

SIGNA™Works..... 13

SIGNA™Works Features ..... 31

Image Acquisition..... 40

RF Coils Suite ..... 52

eXpress Workflow..... 67

Visualization..... 75

Siting..... 80



GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

System



# Magnet

## The Foundation of Quality

When it comes to delivering on the promise of 3.0T image quality while enhancing the openness and patient experience, no other MR component has greater impact than the magnet. The Discovery™ MR750w GEM 3.0T system features a compact, lightweight, superconducting magnet designed to provide excellent homogeneity ensuring uniform signal and fat-suppression over a larger FOV. While improving the patient experience with a 70 cm bore size, the Discovery™ MR750w GEM 3.0T magnet supports a large 50 cm FOV and may reduce exam time since fewer acquisitions are needed to cover a large anatomy.

### Magnet Specifications

Operating field strength	3.0 Tesla
Magnet shielding	Active
EMI shielding factor	97.5% 0.6 hertz excitation and 94.5% dc step
Size (without enclosures) (L x W x H)	1.74 x 2.12 x 2.30 meters
Size (with enclosures) (L x W x H)	1.96 x 2.13 x 2.40 meters
Magnet weight with cryogenics	14,060 lbs (6,378 kg)
Magnet cooling	Cryogenic
Long-term stability	< 0.1 ppm per hour over a 24 hour period
Cryogen refill period	Zero-Boil Off†
Boil-off rate	Zero-Boil Off†
Fringe field – (axial x radial)	7.8 m x 4.8 m at 1 Gauss 5.2 m x 2.8 m at 5 Gauss
Manufacturer	GE Healthcare

### Patient Focused Design

Patient Bore (L x W x H)	130 cm x 70 cm x 70 cm
Patient Aperture	78 cm at magnet flair 70 cm at isocenter
Patient comfort module	Head- or feet-first imaging Dual-flared patient bore 2 way in-bore intercom system Adjustable in-bore lighting system Adjustable in-bore patient ventilation system

### LV-vrms Homogeneity Specifications

Diameter of Spherical Volume – DSV	Guaranteed ppm	Typical ppm
10 cm		0.02
20 cm	< 0.050	0.03
30 cm	< 0.150	0.08
40 cm	< 0.500	0.27
45 cm	< 1.500	0.7
40 (z) x 50 cm	< 3.000	1.8
50 (z) x 50 cm	< 4.000	2.5

†Normal Operating Conditions

# Gradients

The gradients of an MR system play a crucial role when it comes to imaging performance, throughput, and consistency during clinical practice. Gradient speed, accuracy, and reproducibility often determine the success of demanding acquisitions.

## **ART (Acoustic Reduction Technology) Quiet Technology**

State-of-the-art clinical imaging demands the routine use of ultra-fast imaging techniques. At 3.0T, the strong gradients interact with the magnetic field to create mechanical forces resulting in acoustic noise. GE has implemented Quiet Technology on many components of the system to reduce acoustic noise and improve the patient environment.

## **Gradient Coil Isolation and Acoustic Damping**

The full performance of the Extreme Gradient Driver is used while helping to maintain a safe environment for the patient. Clear separation between the gradient coil, RF body coil, and patient support structures ensures minimal component interactions. In addition, mass-damped acoustic barriers are used under the system enclosures to further reduce acoustic noise for the patient.

## **RF Coil Isolation**

During gradient pulses, the RF body coil acts as a secondary source of noise. To further reduce the noise heard by the patient, the RF body coil mounting has been optimally designed with features to reduce acoustic noise.

## **Vibro-Acoustic Isolation**

To isolate the magnet from the building and reduce the transmission of acoustic noise in the structure, GE has designed a vibro-acoustic dampening pad that sits under the feet of the magnet. The dampening characteristics of the pad are optimized based on the magnet geometry and weight.

## **Gradient Waveform Optimization**

User selectable mode to further reduce acoustic noise.

### **Gradient Performance\***

Peak amplitude	44 mT/m
Slew-rate	200 T/m/s
Maximum FOV (x,y,z)	50 cm x 50 cm x 50 cm
Duty Cycle	100%

\*Based on the product of the gradient coil gain, inductance, and peak gradient amplifier current & voltage.

### **Gradient Amplifier & Coil (water-cooled)**

Peak amplifier current and voltage	830A/1650V
Control	Full-digital control
Frequency dependent feed-forward model to match amplifier output to gradient coil	
Dedicated active feedback control loop to regulate current errors	
Gradient current accuracy	300 uAs
Shot-to-shot repeatability <sup>†</sup>	100 uAs
Symmetry <sup>†</sup>	200 uAs

<sup>†</sup>Typical gradient fit expressed in terms of the absolute integrated errors in micro-Amperes-second (uAs). Gradient integral precision is the maximum integrated current error over a full-scale, echo-planar gradient waveform. Shot-to-shot repeatability is the largest difference between integrated errors across waveforms. Symmetry is the largest difference in integrated current error when comparing positive and negative gradient waveforms.






The RF acquisition technology of the Discovery™ MR750w GEM 3.0T enables greater clinical performance and higher image quality especially for data-intensive applications and provides an improvement in SNR versus previous generation systems.

## OpTix Optical RF Architecture

Simultaneous RF Receivers (A/D Converters)	32
Coil input ports	144 with 40 ports embedded in table
Receiver sampling per channel	80 Mhz
Quadrature demodulation	Digital
Receiver dynamic range at 1 Hz BW	> 165 dB
Receiver resolution	Up to 32 bits

## MultiDrive

At 3.0T, precise control over the RF environment in a 70 cm patient bore has been challenging until now. The Discovery™ MR750w GEM 3.0T RF transmit architecture technology consists of a liquid-cooled 30 kW solid-state RF power amplifier with multiple output channels. By optimizing the phase and amplitude of each RF amplifier output channel that is applied to GE's 70 cm 4-port drive whole-body RF transmit coil, 4-port drive whole-body RF transmit coil improves RF uniformity and signal homogeneity regardless of patient shape, size, and/or body habitus.

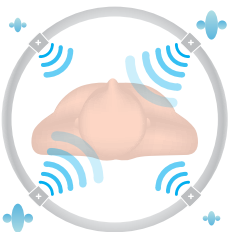

## RF Transmit Architecture

RF amplifier	Multiple output Small footprint Water cooled
Maximum output power	15 kW body per channel (30 kW peak total) 4.5 kW Head
Maximum B1 field with whole body RF coil	16 uT at 75 kg (> 25 uT at 20 kg)
Transmit gain	40 db coarse, > 84 dB instantaneous
RF exciter frequency range	127.72 +/- 0.625 MHz
Receiver resolution	< 0.6 Hz/step
Frequency stability	14 parts per billion (0 to 50 C)
Phase resolution	0.005 deg/step
Amplitude control	16 bit with 12.5 ns resolution
Amplitude stability	< 0.1 dB over one minute at rated power
Digital RF pulse control	2 amplitude modulators 2 frequency/phase modulators
Transmit / Receive Body Coil	Fully integrated 4 port drive, 16 rung quadrature birdcage 70 cm inner diameter 50 cm FOV



# Volume Reconstruction Engine & Host Computer

Reconstruction performance today is challenged by explosive growth in data, and increased computational complexity. The amount of data to be stored and processed continues to increase with the advances in MR system technology. The Discovery™ MR750w GEM 3.0T meets that challenge head-on with innovations in reconstruction to take full advantage of computing power and by leveraging both hardware and software technology.

Reconstruction System Gen6		
	PERFORMANCE	ADVANCED
Operating system	Scientific Linux	Scientific Linux
Processor	Intel® Xeon® E5-2680 v3 (12 core) x 2	Intel® Xeon® E5-2967 (2 x 14 core)
Clock rate	2.5 GHz	2.6 GHz
Memory	160 GB	256 GB
Network	1 GbE	10 GbE
Hard disk storage	2 x 400 GB SSD	3 x 400 GB SSD
2D FFT/second (256 x 256 full FOV)	62,000 2D FFTs/second	75,000 2D FFTs/second
GPU	NA	NA



# Volume Reconstruction Engine & Host Computer



## Host Computer

Operating system	Scientific Linux (RT)
Processor	Intel® Xeon® E5-1620 v3 (4 core, 8 threads)
Clock rate	3.5 GHz
Memory	32 GB
Network	Gigabit (10/100/1000) Ethernet
Hard disk storage	1024 GB SSD
Graphics subsystem	PCI-Express x16 2GB Single DVI-I 3D
Media drives	CD/DVD drive
Cabinets	Single, tower configuration

## Orchestra Reconstruction Platform

Orchestra is a high performance computing software library toolbox that enables new possibilities for integration of advanced reconstruction elements. Delivering enhanced productivity gains by increased image reconstruction speed and minimizing workflow disruptions. A powerful platform not only built to support the most demanding applications such as HyperSense, but also to provide our collaborators with easy access to the product reconstruction algorithms.





# Computing Platform

## Operator Console

The Discovery™ MR750w GEM 3.0T system comes equipped with a scan control keyboard assembly that contains intercom speaker, microphone and volume controls, and an emergency stop switch. Start-scan, pause-scan, stop-scan, and table advance to isocenter hot keys are also included.

## Display and DICOM Data

The Discovery™ MR750w GEM 3.0T system generates MR Image, Secondary Capture and Grayscale Softcopy Presentation State (GSPS) DICOM objects. The DICOM networking supports both send and query retrieve as well as send with storage commit to integrate with the site's PACS archive. DICOM filming support includes both Basic Grayscale and Basic Color Print Service Classes. Additionally, the Discovery™ MR750w GEM 3.0T system supports the CT and PET image objects for display allowing the user to refer to cross-modality studies.

Display	
AutoView	560 x 560 Image Window (standard)
Window/Level (W/L)	7 user-programmable keys on scan control keyboard plus one key for returning to prior setting 6 user-programmable buttons in image viewer Arrow keys on scan control keyboard On-image through middle mouse button Save State stores user-selected image orientation, user annotation and window level
Image display	Zoom/Roam/Flip/Rotate/Scroll/Explicit Magnify and Magnifying Glass Image Measurement Tools Grid On/Off Cross Reference/User Annotation Exam/Series Page Hide Graphics/Erase Annotation/Screen Save Accelerator Command Bar Compare Mode/Reference Image Scoutview Cine Paging (up to 4 windows and 128 images/window) Add/Subtract/Edit Patient Data
Image display performance	256 Image buffer (256 x 256) at 30 fps
Image annotation	Shadowed to permit ease in reading Two graphic/text planes overlay the entire screen Grid placement with anatomical reference on an image Drawing and annotation may be added to and removed from images



# Computing Platform (continued)



## Filming

### Filming

- Drag and Drop filming
- One-button Print Series
- One-button Print Page
- Multi-image formats – from 1 to 24 images displayed simultaneously in various layouts
- DICOM Basic Grayscale Print Service Class
- DICOM Basic Color Print Service Class

## Wide-screen display monitor

### Display monitor

- 24" Widescreen LCD Flat Panel
- 1920 x 1200 dot resolution
- Non-interlaced, flicker-free presentation
- Contrast ratio 1000:1
- Digital DVI Interface



# Scan Parameters

Sequences	Parameters	Matrix 64	Matrix 128	Matrix 256	Matrix 512
2D Spin Echo	Min. TR (ms)	N/A	2.8 ms	3.5 ms	5.0 ms
	Min. TE (ms)	N/A	1.624 ms	1.984 ms	2.816 ms
2D Fast Spin Echo	Min. TR (ms)	N/A	3.4 ms	3.9 ms	8.1 ms
	Min. TE (ms)	N/A	1.608 ms	1.952 ms	2.792 ms
	Min. slice thickness	0.2 mm			
	Min. ESP (ms)	N/A	1.608 ms	1.952 ms	2.792 ms
	Max. ETL	480			
3D Fast Spin Echo	Min. TR (ms)	N/A	41 ms	53 ms	74 ms
	Min. TE (ms)	N/A	5.0 ms	7.0 ms	10.0 ms
	Min. slice thickness	0.3 mm			
	Min. ESP (ms)	N/A	1.672 ms	2.304 ms	3.664 ms
	Max. ETL	395			
2D Fast Gradient Echo	Min. TR (ms)	0.532 ms	0.66 ms	0.924 ms	1.324 ms
	Min. TE (ms)	0.184 ms	0.184 ms	0.188 ms	0.192 ms
3D Fast Gradient Echo	Min. TR (ms)	0.540 ms	0.660 ms	0.900 ms	1.270 ms
	Min. TE (ms)	0.184 ms	0.184 ms	0.18 ms	0.184 ms
Inversion Recovery	Min. TR (ms)	N/A	58.5 ms	59.4 ms	61.5 ms
	Min. TE (ms)	N/A	1.608 ms	1.952 ms	2.792 ms
	Min. TI (ms)	N/A	50 ms	50 ms	50 ms



# Scan Parameters (continued)

Sequences	Parameters	Matrix 64	Matrix 128	Matrix 256	Matrix 512
3D FIESTA	Min. TR (ms)	0.940 ms	1.248 ms	1.930 ms	3.36
	Min. TE (ms)	0.248 ms	0.324 ms	0.452 ms	0.628
Echo Planar Imaging	Min. TR (ms)	4.0 ms	5.0 ms	6.0 ms	N/A
	Min. TE (ms)	1.1 ms	1.4 ms	1.8 ms	N/A
	Min. FOV [cm]	4 cm			
	ESP at 25 cm	0.456 ms	0.656 ms	1.056 ms	N/A
	ESP at 48 cm	0.328 ms	0.460 ms	0.672 ms	N/A
	ESP at 99 cm	0.228 ms	0.320 ms	0.556 ms	N/A
	Images per second	138	71	36	N/A
	b value	Maximum(s/mm2): 10.000 Max # for ADC: 40			
	Diffusion Tensor directions	Max: 150			
	SLICE THICKNESS and FOV				
Minimum slice thickness in 2D					0.2 mm
Minimum slice thickness in 3D					0.1 mm
Min/Max FOV					10 mm / 500 mm
Min/Max Matrix					32-1024



GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

SIGNA™ Works




# SIGNA™Works

The latest software platform provided by GE, it includes the base pulse sequences, workflow enhancements and visualization tools to enable high productivity with exceptional quality and outcomes. SIGNA™Works, starting with the acquisition provides the tools needed to enable superb results in the various clinical fields. With 6 optimized Works categories, GE delivers preset protocols for the most demanding Neuro, Muskuloskeletal, CardioVascular, Body, Oncology and Paediatric areas. In addition to enabling the routine imaging, SIGNA™Works provides the user with a streamlined and efficient operating environment with in-line processing through single-click outcomes for even the most demanding processes.

## **SIGNA™Works provides:**

- Software platform with a wide range of assets for image acquisition, display and post processing.
- Strategically packaged to deliver speed, high quality diagnostic images and reliable post processing to each clinical area.
- An intelligent combination of MR pulse sequences and advanced techniques, designed to bring solutions for enhanced care and productivity.
- From SE, FSE, frFSE, Inversion Recovery, SSFSE, SSFSE-IR, GRE, FGRE, SPGR, FSPGR to Volumetric imaging, Motion Correction, Diffusion Weighted, Vascular imaging and beyond.

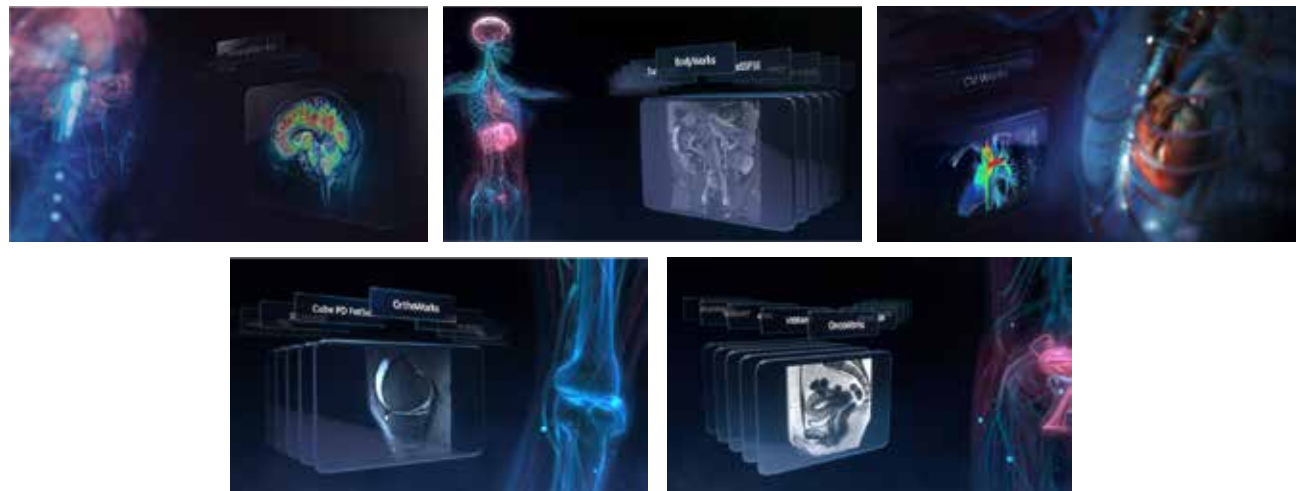





NeuroWorks includes the basic imaging acquisitions and processing to the latest in motion correction, functional and volumetrics. Supporting both simple reconstruction to real-time perfusion results with Brainstat AIF.


## Volumetric Imaging

3D Cube	PD, T1, T2, T1 FLAIR and T2 FLAIR Isotropic high resolution volumetric One sequence, reformat in all planes
3D Cube DIR	DIR, typically but not limited to CSF and white matter suppression
BRAVO T1	BRAVO, 1 mm isotropic sequence of choice for functional data overlay
Visualization	3D reformat MPR Volume segmentation Volume rendering Auto-contour

## Motion Correction

PROPELLER MB	Multiple contrasts – T1, PD, T2, T1 FLAIR, T2 FLAIR and DWI Motion reduction Magnetic susceptibility effects reduction
Visualization	Registration Motion correction



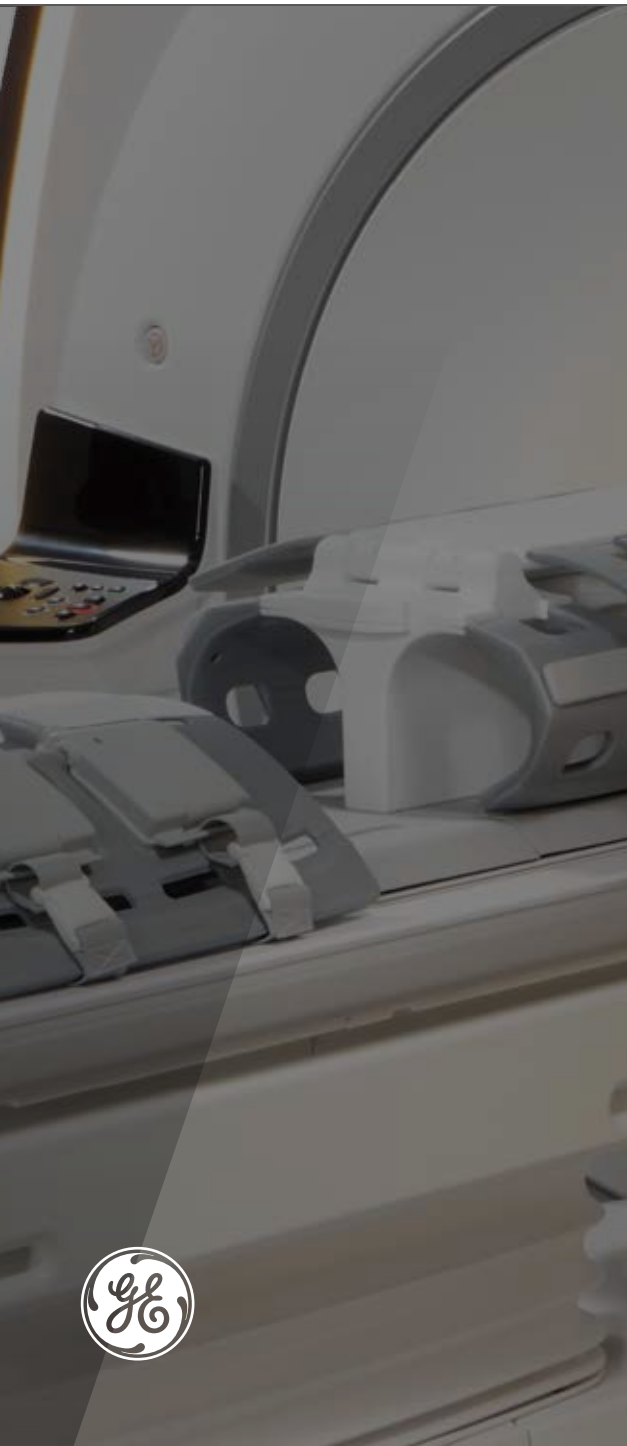

## Enhanced Diffusion Weighted

eDWI	Multi b-value 3:1, Tetrahedral Smart NEX Inversion recovery for robust FatSat RTFA: Increases SNR by 50% and distortion reduction for accurate post processing when compared to dual spin echo
Visualization	ADC and eADC


## One Touch Protocol

READYBrain	Automated multi series, multi plane prescription Combine with Auto Scan for one touch protocol In-line for Auto Post processing
------------	---------------------------------------------------------------------------------------------------------------------------------------

## Dynamic Brain Function

BrainSTAT Perfusion and Analysis	Blood flow Blood volume Mean transit time Time to peak parametric Fusion
BrainSTAT Arterial Input Function (AIF)	Manage tracer arrival differences due to patient flow dynamics Automatically or manually specify the AIF to normalize maps
Visualization	Brain STAT



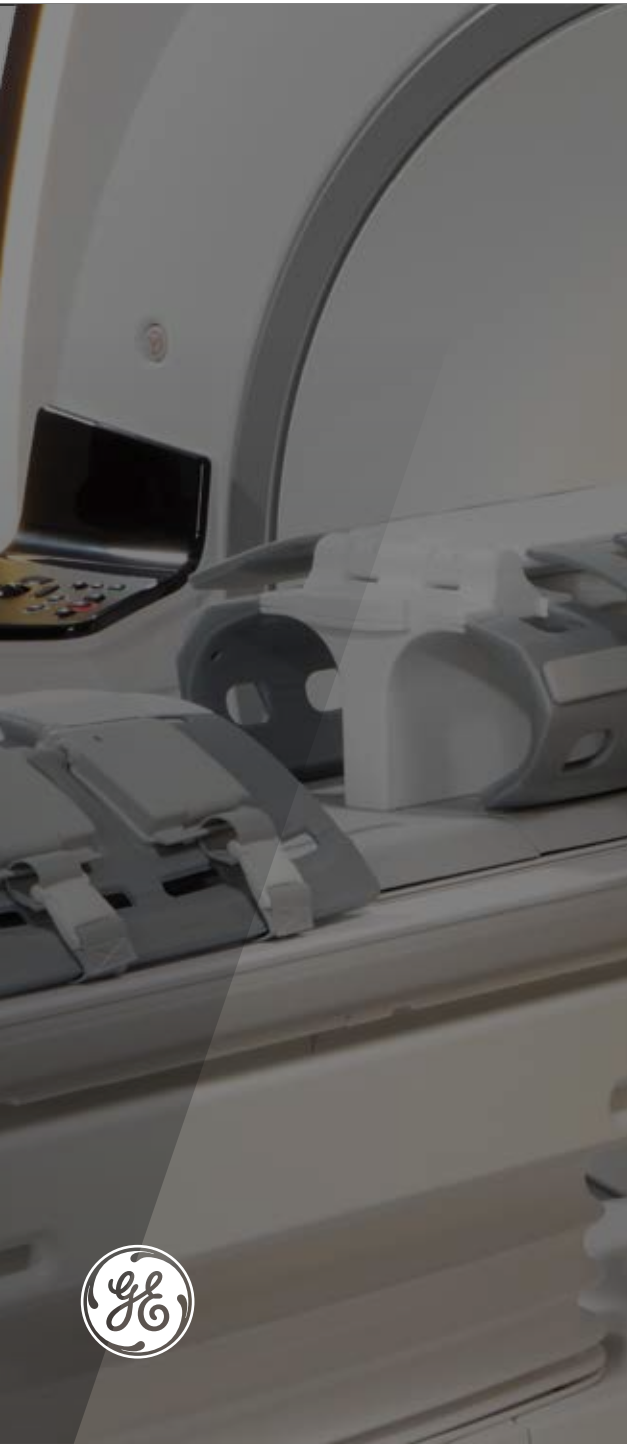

## Spectroscopy

PROBE PRESS	Concentrations of in-vivo metabolites evaluation
	Acquisition and display
	Reduced flip angles for lower min TE values
	Up to twice the SNR when compared to PROBE STEAM
Visualization	Brain Spectroscopy

## Spine Imaging

2D/3D MERGE	High SNR T2* contrast
	Gray/white matter differentiation
	Foraminal detail
3D COSMIC	SSFP to emphasize T2 signal for improved contrast
	Nerve root and disc detail
Visualization	3D reformat MPR
	Volume segmentation
	Volume rendering





The latest in torso imaging is delivered with volumetric imaging supporting advanced parallel imaging standard. Including, Snapshot imaging with optimized Single Shot FSE, 3D isotropic imaging for MRCP, Dynamic Imaging and Routine Volumetric imaging enabled with Motion Free navigation for post-contrast uses with high temporal resolution results. Motion correction is further enhanced with both the PB navigators as well as PROPELLER including T1-weighted results. Turbo class of acquisitions, streamlines the speed and enables higher quality results. Advanced processing is made one-touch with the new READYView on Console capabilities.


## Volumetric Imaging

3D Cube	Isotropic high resolution volumetric One sequence, reformat in all planes
3D Dual Echo	In- and out-of-phase Used to help identifying fatty infiltration, focal fatty sparing, liver lesions, and other conditions High spatial resolution
Visualization	3D reformat MPR Volume segmentation Volume rendering Auto-contour

## Motion Correction

PROPELLER MB	Motion reduction
PB Navigators	Free breathing tracker
Respiratory Trigger	Free breathing bellows
Visualization	Registration Motion correction



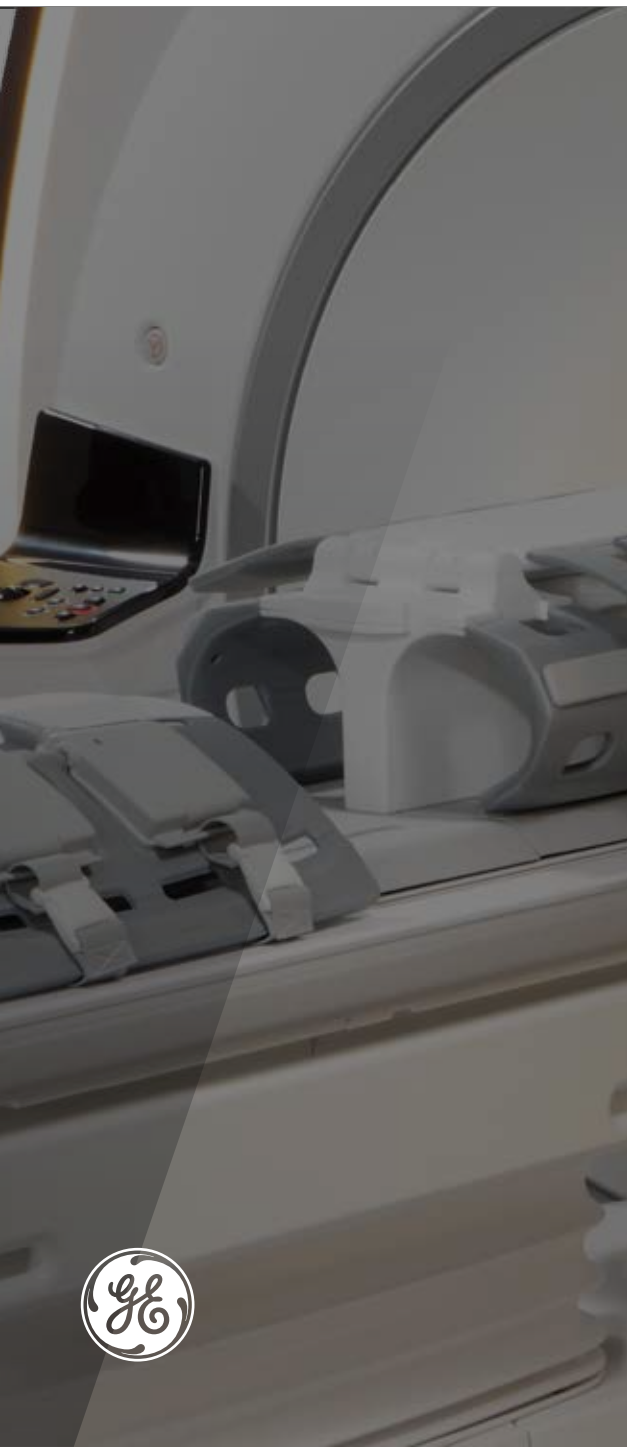


## Enhanced Diffusion Imaging

eDWI	Multi b-value, 3:1, Tetrahedral Smart NEX Inversion recovery for robust fat sat RTFA: Increases SNR by 50% and distortion reduction for accurate post processing when compared to dual spin echo
Visualization	ADC and eADC Fusion

## Dynamic Body Imaging

LAVA	SPGR Fast Liver Acquisition SPECIAL for robust fat suppression
Turbo LAVA	ARC acceleration for full organ coverage Shorter breath holds
Multi-phase Dynaplan	Customizable phase delay for dynamic studies Series per phase Auto subtraction Pause after mask
Visualization	MR standard SER



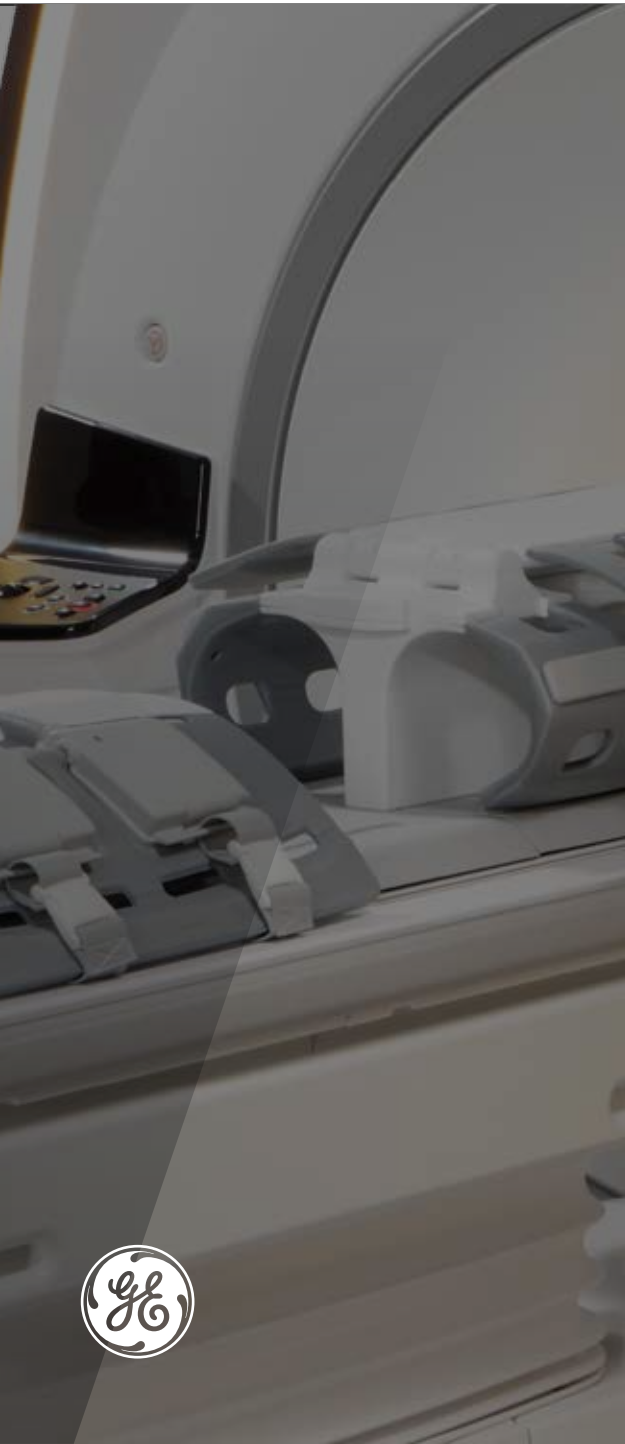

## Non-Invasive Non Contrast Biliary System – MRCP

3D frFSE MRCP	T2 Prep for background suppression Breath-hold and PB navigator
2D SSFSE	T2-weighted, with sub second single slice acquisition High signal from fluids Good suppression of other tissues Snapshot acquisition, motion artifacts virtually eliminated Thin slices and thick slab protocols Single breath-hold acquisition MIP post processing
2D FatSat FIESTA	Excellent contrast between ducts and gallbladder with surrounding anatomy Fat sat for increased conspicuity
2D frFSE	T2-weighted High resolution Supplementary information for assessment of extra ductal masses
Visualization	3D Reformat MPR MIP & HD MIP







CVWorks provides GE's extensive coverage for the latest techniques enabling high performance CardioVascular imaging outcomes. Single Breath-hold imaging for whole heart coverage are available from Morphology to Delayed enhancement. Enabling simplified generation of superb results including head-to-toe MRA support to single acquisition Time of Flight and additional non-contrast imaging for flow. With SmartPrep and Fluoro triggering enable for first-time right contrast injections.

## Myocardium Delayed Enhancement

### MDE PLUS

Single-Shot Myocardial Delayed Enhancement (SSH MDE)	<ul style="list-style-type: none"> <li>Shorten breath-holds or free breathing for better patient tolerance</li> <li>Potential for reduced scan time</li> <li>Imaging arrhythmic patients</li> <li>Snapshot imaging for motion reduction</li> </ul>
------------------------------------------------------	----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------


Adiabatic IR Pulse	<ul style="list-style-type: none"> <li>Robust Myocardial Suppression</li> <li>Fat Suppression</li> <li>Adiabatic fat suppression pulse</li> <li>Improved characterization of enhancing tissue</li> </ul>
--------------------	----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

MDE Plus: Phase Sensitive MDE (PSMDE)	<ul style="list-style-type: none"> <li>Inversion Recovery FGRE sequence</li> <li>Phase-sensitive image reconstruction</li> <li>Consistent myocardial suppression, even with sub-optimal TI</li> <li>Improved contrast for myocardial</li> <li>Potential to shorten overall exam time</li> </ul>
---------------------------------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

## Single Breath Hold Whole Heart

Black Blood SSFSE	<ul style="list-style-type: none"> <li>Difficult patients with irregular heartbeats or limited breath-hold capacity</li> <li>Potential to shorten exam times</li> <li>Shorten breath-holds for better patient tolerance</li> <li>Whole chest survey</li> </ul>
-------------------	----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------



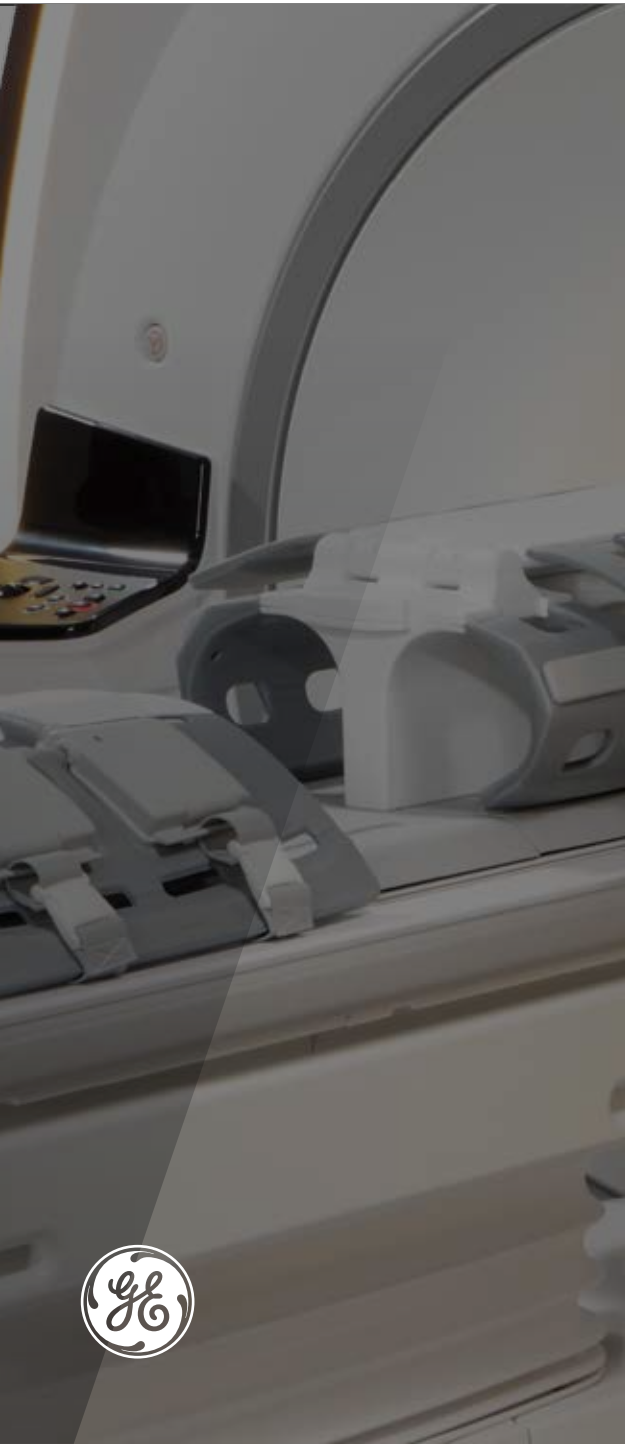

## Viability Imaging

CINE IR	<p>Multiphase FGRE Cine acquisition...quick assessment of optimal TI time for MDE</p> <p>Captures image contrast evolution at different TI times</p> <p>Adiabatic Inversion Recovery for uniform myocardial suppression</p> <p>Support both 1 RR and 2 RR mode</p>
---------	--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

## T2\* Mapping

StarMap	<p>T2* mapping compatible with gating for cardiac evaluation</p> <p>Non-invasive evaluation of the entire organ</p>
READYView	R2 Star

## Navigator Free Breathing Acquisition

Navigators	<p>Used with 3D IR Prepared FGRE or 3D Fat Sat FIESTA</p> <p>Free breathing navigator diaphragm tracking</p>
------------	--------------------------------------------------------------------------------------------------------------


## Flow Imaging

Flow Analysis	<p>Flow velocity and volume flow quantification</p> <p>Peak and average flow charts and graphics</p> <p>Automated contour detection</p> <p>Brain, chest and abdominal clinical applications</p>
---------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

## Contrast Enhancement Tracking

SmartPrep	Automated bolus tracking
Fluoro triggered	Real Time bolus tracking
Visualization	MIP & HD MIP



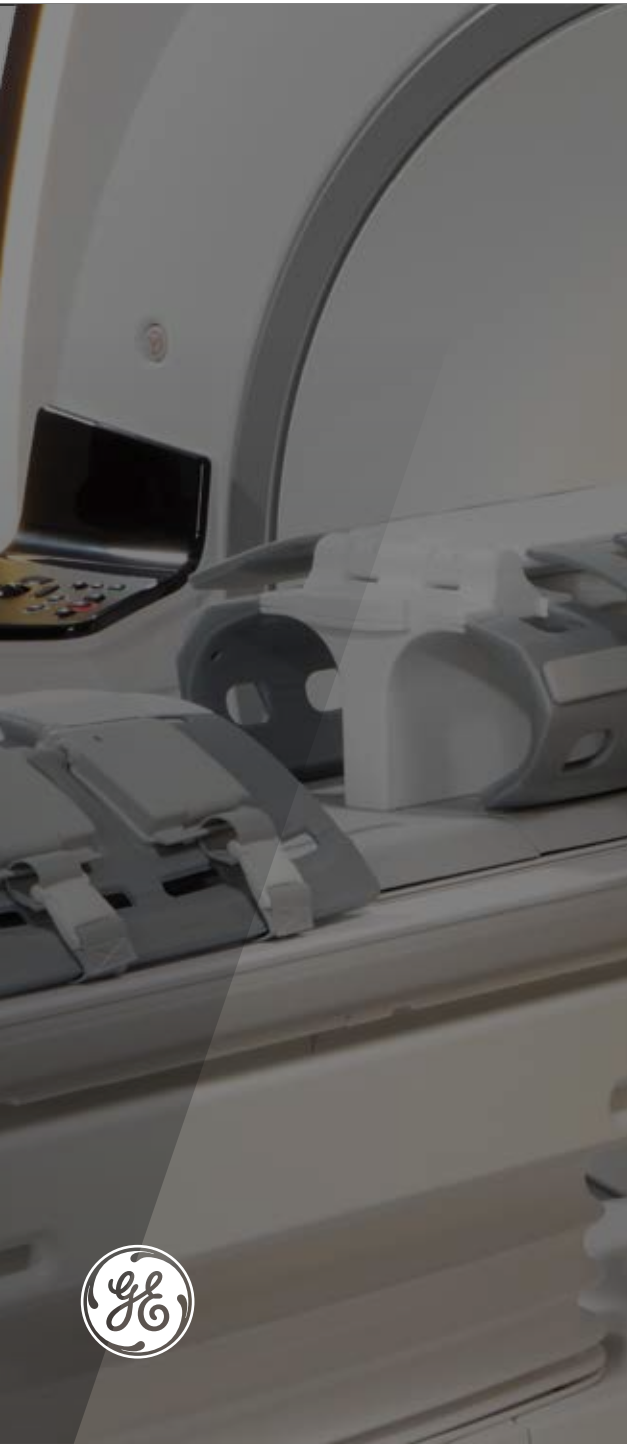

## Peripheral Vascular Runoff

QuickStep	Multi-station, multi phase acquisition
	Automatically prescribes, acquires, and combines images from multiple stations
	Entire exam complete with no user intervention in as little as 7 minutes
	Auto subtraction

## Non Contrast Vascular Imaging

2D TOF	Carotid bifurcation, venous anatomy, aortic arch, peripheral vesselsr
3D TOF	Circle of willis, intracranial vasculature, abdominal vasculature
3D TOF Multi Slab	Intracranial vasculature, carotid bifurcation, aortic arch, peripheral vessels, venous anatomy
2D Phase Contrast	Localizer, flow direction and velocity for intracranial and extracranial vasculature, portal or hepatic vein, quantitative measurement of flow velocity
3D Phase Contrast	Intracranial vasculature, renal arteries
Visualization	MIP & HD MIP





OrthoWorks delivers routine imaging that is not always a given. From motion correction to advanced volumetric imaging, GE's latest MSK techniques provide you with the contrasts you need for the basic imaging to enhanced cartilage imaging. And with multiple tissue suppression methods available, OrthoWorks enables the best of what can be achieved in a standard configuration.

## High Resolution Imaging

FSE & frFSE	Intermediate PD, T1, T2-weighted imaging
	Compatible with FatSat, STIR, SPECIAL and ASPIR
	Gold standard for articular cartilage, cartilage ligaments, menisci and subcondral bone


## Volumetric Imaging

3D Cube	PD, T1, T2
	Isotropic high resolution volumetric
Visualization	One sequence, reformat in all planes
	3D reformat MPR
	Volume segmentation
	Volume rendering

## Motion Correction

PROPELLER MB	Multiple contrasts – T1, PD, T2
	Motion reduction
Visualization	Registration
	Motion correction



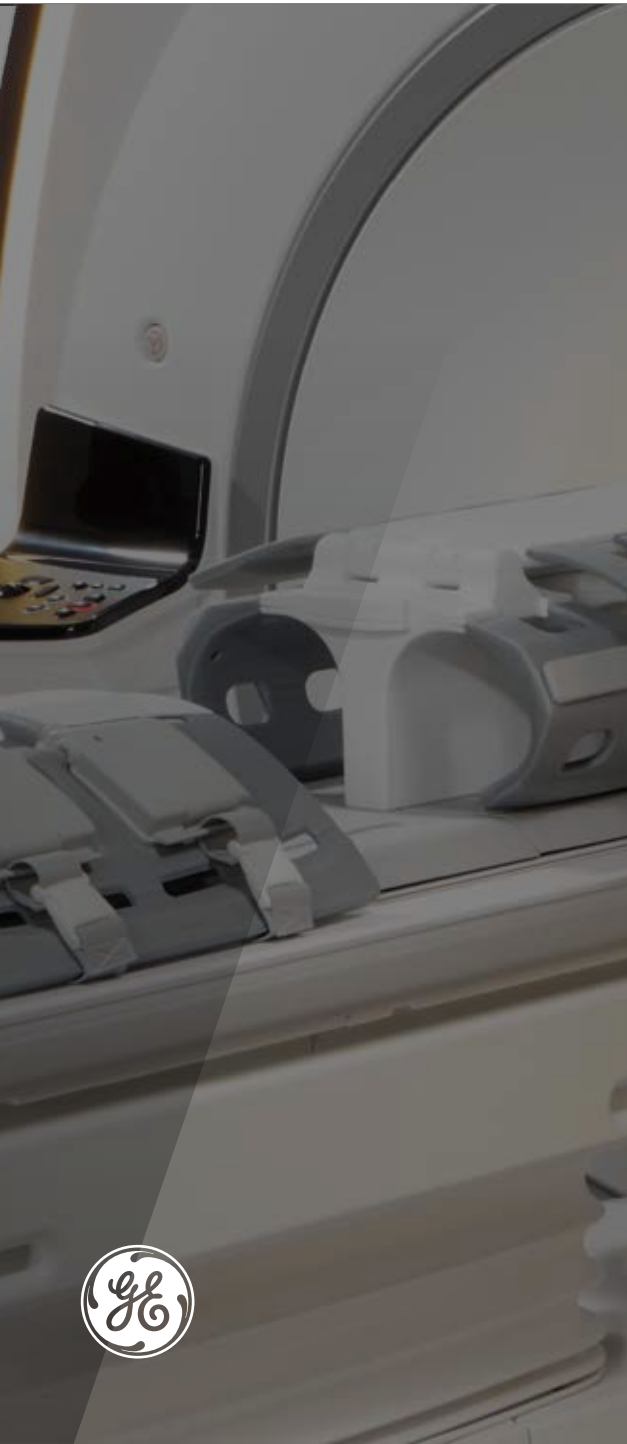

## T2\*-weighted Imaging

3D MERGE	High SNR T2* contrast Visualization of ligaments while adding soft tissue contrast Reduced chemical shift
3D COSMIC	Fast, high resolution volumetric imaging SSFP to emphasize T2 signal for improved contrast
Visualization	3D reformat MPR Volume segmentation Volume rendering

## Artifact Reduction Standard Sequence

MARS	FSE High bandwidth protocols High resolution, small FOV imaging
------	--------------------------------------------------------------------

## Fat Suppression

Chemical FatSat	Frequency selective fat saturation
STIR	Inversion recovery fat null point method
ASPIR	Solution for poor fat suppression due to B <sub>1</sub> inhomogeneity
SPECIAL	Hybrid method between chemical FatSat and STIR
Spectral Spatial	Water excitation only



# OncoWorks

OncoWorks delivers a complete platform for your needs in prostate, breast and radiation therapy planning. From the basic routine acquisitions to whole body imaging including volumetric and enhanced diffusion capabilities, GE enables superb linearity from the gradient platform and hardware performance. GE provides the necessary preset protocols to supply you with optimal imaging for your oncology needs that is further enhanced visualization capabilities so that your results can be a single click away.

## Volumetric Imaging

3D Cube	PD, T1, T2, T1 FLAIR and T2 FLAIR Isotropic high resolution volumetric One sequence, reformat in all planes
3D Cube DIR	DIR, typically but not limited to CSF and white matter suppression
BRAVO T1	BRAVO, 1 mm isotropic sequence of choice for functional data overlay
Visualization	3D reformat MPR Volume segmentation Volume rendering Auto-contour

## Enhanced Diffusion Weighted

eDWI	Multi b-value 3:1, Tetrahedral Smart NEX Inversion recovery for robust fat sat RTFA: Increases SNR by 50% and distortion reduction for accurate post processing when compared to dual spin echo
Visualization	ADC and eADC





# OncoWorks (continued)



## Dynamic Imaging

Multi-phase SPGR	SPGR dynamic fast acquisition SPECIAL for robust fat suppression
Visualization	MR standard SER

## Whole Body Scanning

FSE-IR/3D SPGR/DWI	Whole body imaging Multiple stations with large FOV Metastasis screening
--------------------	--------------------------------------------------------------------------------




# PaedWorks

PaedWorks is the GE solution to address your specific needs in Paediatric imaging, from standard sequences supported with the latest in motion control for brain to toes. GE delivers standard acoustic reduction technologies and further addresses clinical needs for volumetric imaging, whole body imaging and enhanced diffusion results. The streamline processing enables simplified one-click processing and visualization of complex results. PaedWorks covers your needs for all anatomies and provides optimized protocols and preset procedures.

Volumetric Imaging	
3D Cube	PD, T1, T2, T1 FLAIR and T2 FLAIR Isotropic high resolution volumetric One sequence, reformat in all planes
3D Cube DIR	DIR, typically but not limited to CSF and white matter suppression
BRAVO T1	BRAVO, 1 mm isotropic sequence of choice for functional data overlay
3D Dual Echo	In- and out-of-phase used to help identifying fatty infiltration, focal fatty sparing, liver lesions, and other conditions High spatial resolution
Visualization	3D reformat MPR Volume segmentation Volume rendering



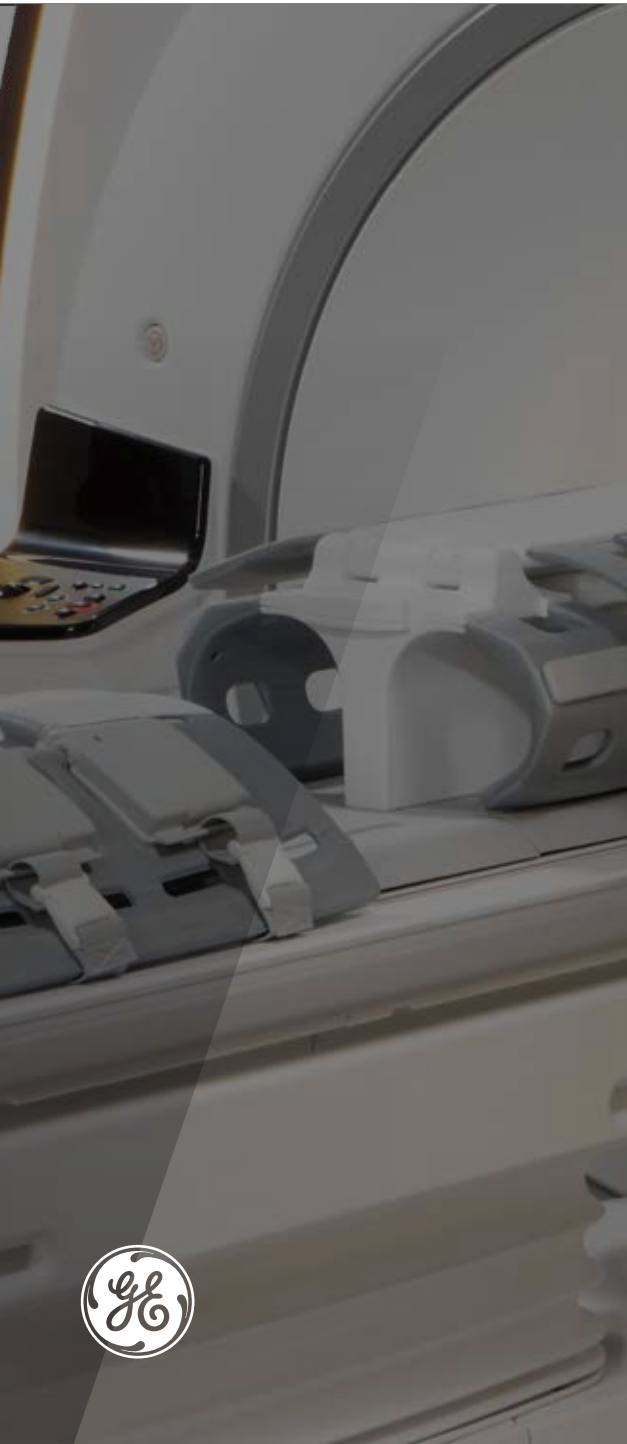

## Motion Correction

PROPELLER MB	Motion Reduction
PB Navigators	Free breathing tracker
Visualization	Registration Motion correction


## One Touch Protocol

READYBrain	Automated multi series, multi plane prescription
Not recommended for under 1 year of age	Combine with auto scan for one touch protocol In-line for auto post processing

## Dynamic Brain Function

BrainSTAT Perfusion and Analysis	Blood flow
	Blood volume
	Mean transit time
	Time to peak parametric
	Fusion
BrainSTAT Arterial Input Function (AIF)	Manage tracer arrival differences due to patient flow dynamics Automatically or manually specify the AIF to normalize maps
Visualization	BrainSTAT



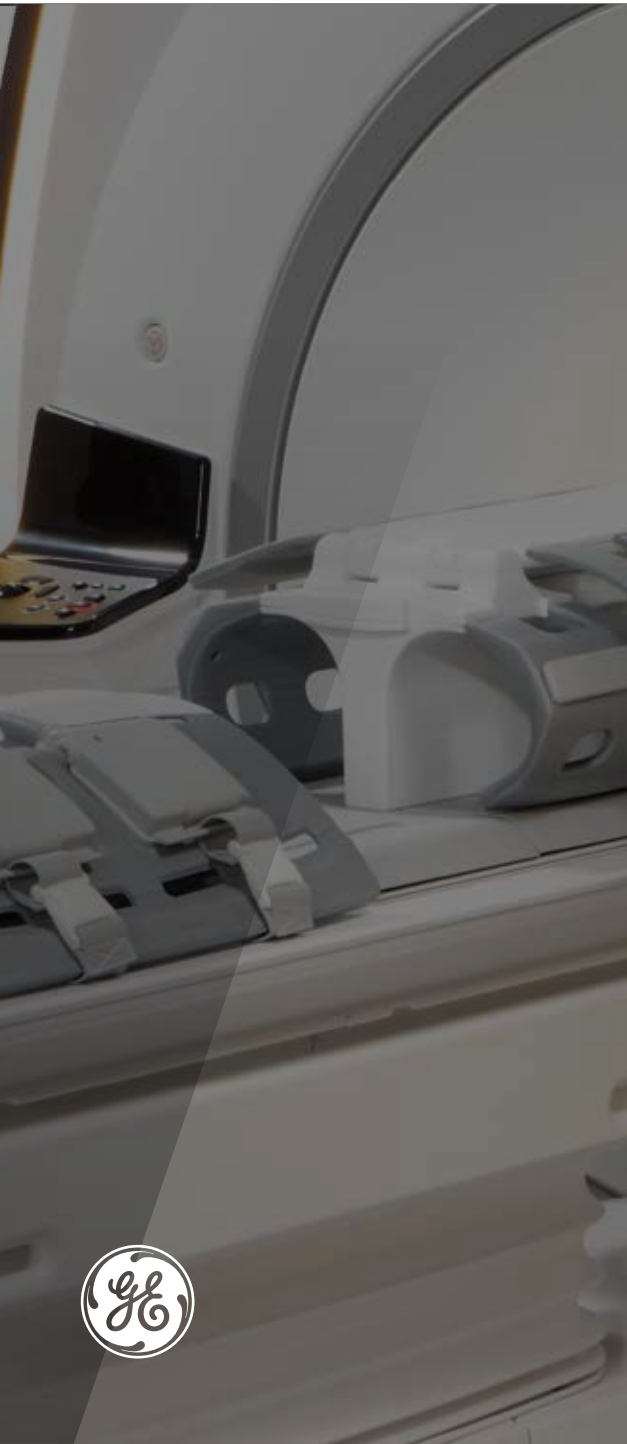

## Spectroscopy

PROBE PRESS	Concentrations of in-vivo metabolites evaluation
	Acquisition and display
	Reduced flip angles for lower min TE values
	Up to Twice the SNR when compared to PROBE STEAM
Visualization	Brain spectroscopy

## Spine Imaging

2D/3D MERGE	High SNR T2* contrast
	Gray/white matter differentiation
	Foraminal detail
3D COSMIC	SSFP to emphasize T2 signal for improved contrast
	Nerve root and disc detail
Visualization	3D reformat MPR
	Volume segmentation
	Volume rendering



# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

SIGNA™ Works  
Features






## Benefits

- Increase productivity by reduced scan times
- Faster 3D imaging acquisitions
- Combined with ARC for higher acceleration factors

## Going further than common sense

HyperSense is an acceleration technique based on sparse data sampling enabling faster imaging without the penalties commonly found with conventional parallel imaging.

HyperSense is intended to be used with volumetric acquisitions, it is combined with (ARC) parallel imaging delivering optimal signal to noise ratio with shorter acquisition times.

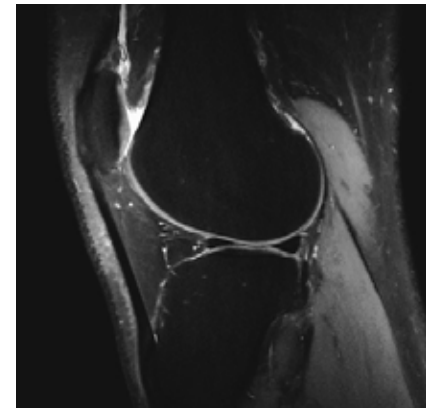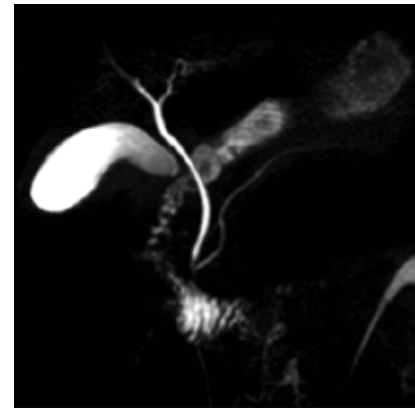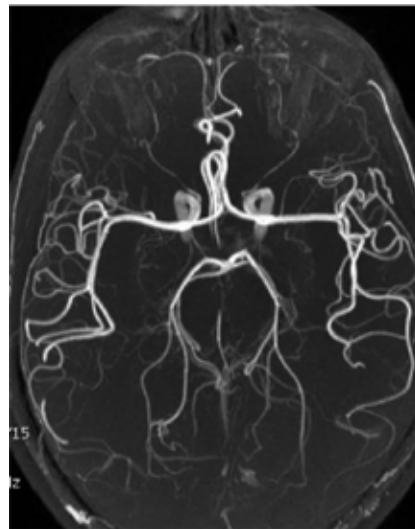






## Benefits


- Simultaneous excitation: multiple slices at multiple locations
- Acquisition time reduction without compromising post processing metrics
- More diffusion directions, number of slices or higher temporal resolution without extra scan time
- Shorter breath-holds
- Combine with ARC for higher acceleration factor

## Quality and Speed Synchronized


HyperBand provides a reduction in scan time by simultaneously exciting multiple slices at multiple locations. It can lead to higher acceleration reduction factors when combined to other methods of parallel imaging. The benefits of HyperBand acceleration include enhancements on productivity and patient experience by shortened breath-holds, increased anatomy coverage and higher resolution image acquisition.

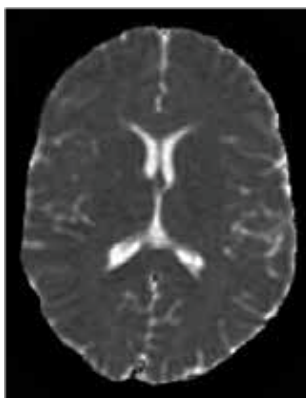

ADC Map

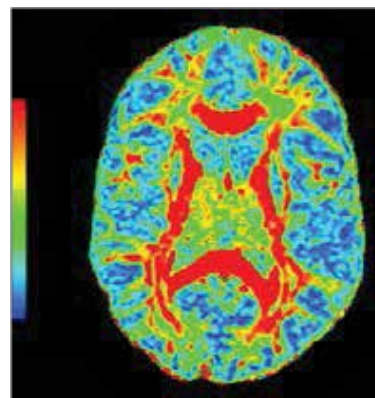

FA Map

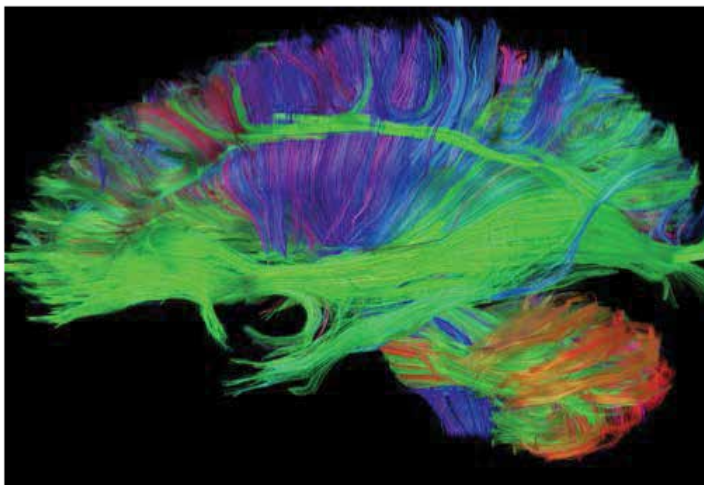

FiberTracks




## Benefits


- Significant scan time reduction while maintaining SNR efficiency
- High resolution small FOV isotropic volumetric imaging
- FLEX for large FOV robust fat suppression

## Tailored 3D imaging that fits to perfection


Delivers small field of view organ specific volumetric imaging acquisition that can reduce artifacts originating from outside of the prescribed FOV. HyperCube can be applied with or without fat suppression using Flex or chemical saturation methods. Provides significant savings of imaging time without sacrificing contrast quality and it can be used across the entire body.

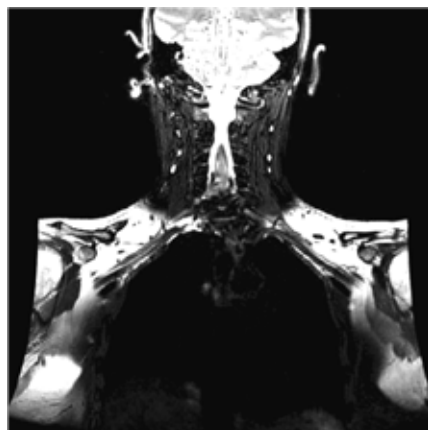


FatSat Failure  
Chem FatSat

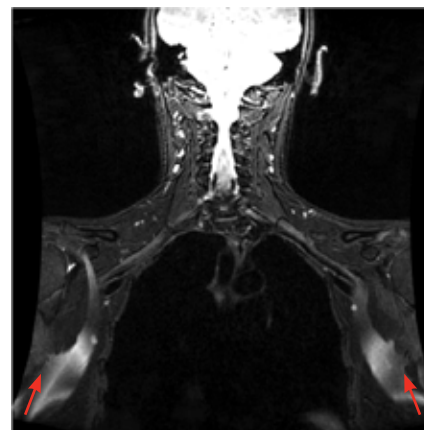

Arm wrap causing swaps  
Cube with Flex

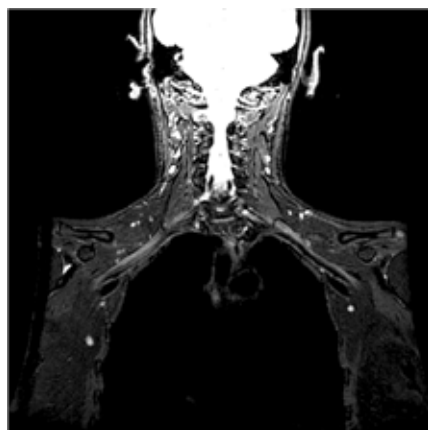

Time 3:37  
HyperCube with Flex

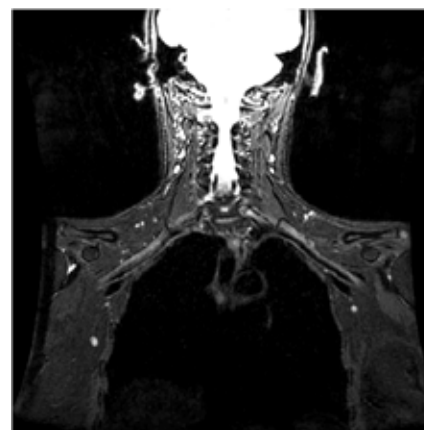

Time 2:37  
HyperCube with Flex and HyperSense

# Flex for Cube and FSE




## Benefits

- 2D and 3D dual echo fat-water separation technique
- Uniform fat suppression for large FOV challenging off-center anatomies
- Dixon based, less sensitive to  $B_0$  inhomogeneity
- Choice of single pass acquisition for significant scan time reduction
- Water, Fat, in-phase and out-of-phase images

## Unlimited solutions, consistent results

Flex uses a dual echo fat-water separation technology to provide robust and homogeneous fat suppressed images. Flex is compatible with ARC acceleration and can be used with a fast triple echo selection for significant scan time reduction. Enhanced uniformity and control of fat water swaps allow large field of view and off-center imaging where uniformity is a challenge. Delivering fast 2D and 3D acquisitions with reconstructed in-phase, out-of-phase, water and fat images, Flex represents productivity gains in all clinical areas.

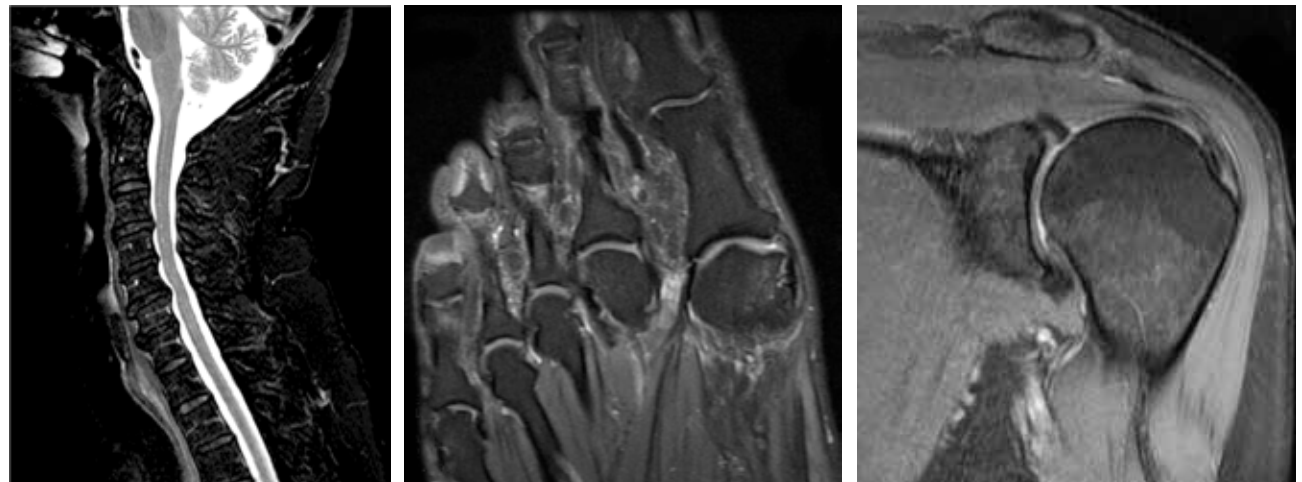




## Benefits

- Multiple contrasts a single scan
- Up to 50% faster than acquiring all contrasts separately\*
- Ability to change the contrast after acquisition by modifying TR, TE and/or TI
- Enhanced image slice registration owing to the absence of inter-acquisition patient motion
- Parametric Maps: T1, T2, R1, R2, PD

## One and Done

MAGiC (MAGnetic resonance image Compilation), enables one and done imaging capability by delivering multiple contrasts in a single scan. MAGiC utilizes a multi-delay, multi-echo acquisition. The data acquired is processed using a technique to generate T1, T2, PD and Inversion Recovery (IR) weighted images (including: T1-FLAIR, T2-FLAIR, STIR, Dual IR and PSIR weighted images), all at once, reducing scan time by up to 50% compared to acquiring all contrasts separately.\* MAGiC generates all the different contrasts from the same acquisition, leading to enhanced image slice registration, owing to the absence of inter-acquisition patient movement. Because of the efficiency of MAGiC, the user has the flexibility to explore more advanced imaging, such as Spectroscopy\*\*, Susceptibility Weighted Imaging\*\* etc., in the same time required to perform the routine exam without MAGiC. MAGiC provides the user the ability to change the contrast of the images after acquisition. This is performed by adjusting the TR, TE, and/or TI parameters post-acquisition, to generate the specific contrast desired. MAGiC also enables users to generate parametric T1, T2, R1, R2, PD maps for further analysis of MRI acquisition data.

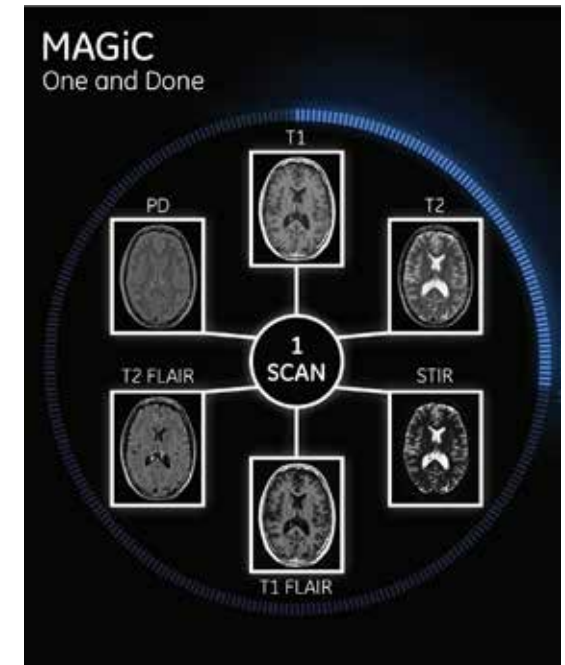

One MAGiC scan delivers six contrasts

\*Based on MAGiC clinical study of 109 patients from 6 separate institutions.

\*\*Optional package (MAGiC in itself does not deliver advanced imaging)

This feature may not be commercially available in all countries. Due to regulatory reasons, future availability cannot be guaranteed. Please contact your local GE MR organization for further details. General Electric Company reserves the right to make changes in specifications and features shown herein, or discontinue any products described at any time without notice or obligation. Please contact your GE representative for the most current information.

© 2016 General Electric Company-All rights reserved. GE, GE Monogram and imagination at work are trademarks of General Electric Company. GE Healthcare, a division of General Electric Company. DOC1710687






# MAGiC DWI




## Benefits

- Multiple synthetic b-values from a single DWI scan
- High b-values in shorter scan times
- Compatible with FOCUS Diffusion


MAGiC DWI generates multiple synthetic b-values from a single DWI scanned series allowing the user to view diffusion contrasts changes in real time after the acquisition. It delivers high b-values without stressing protocol parameters and resulting in shorter scan times without sacrificing contrast or anatomy coverage. Synthetic Diffusion is not limited to diffusion directionality or coil type.

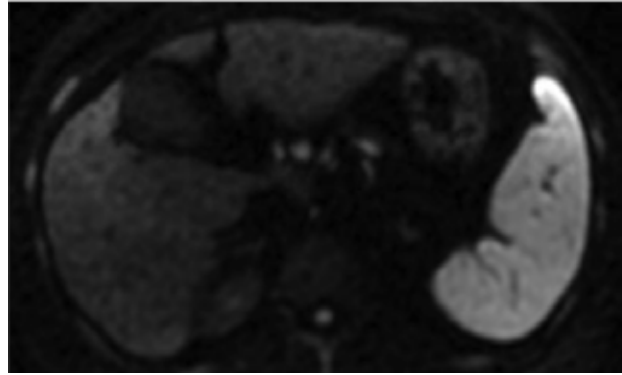


Synthetic b value = 1000

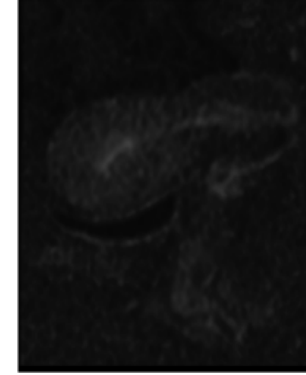


Synthetic b value = 1400

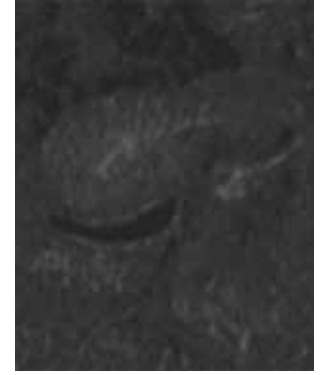

Synthetic b value = 2000

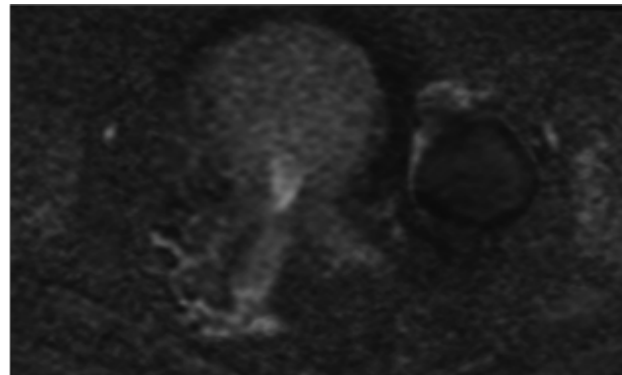

Synthetic b value = 1400

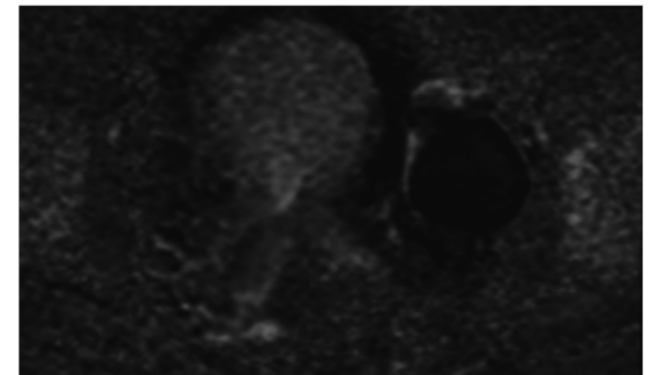

Synthetic b value = 2000






## Benefits

- 3D cine acquisition in any dimension
- Free breathing whole chest coverage
- Allows velocity encoding in all directions
- Single and view sharing frames for higher temporal resolution
- Effortless workflow

## Confident Functional Accuracy

ViosWorks is a 3D cine-based acquisition that can be planned in any dimension and allows for velocity encoding in all directions to assess vascular flow. The acquisition delivers fast imaging with the use of Hyperkat acceleration including both, single and view sharing frames for higher temporal results. Provides high spatial resolution to enable visualization of flow through complex structures.

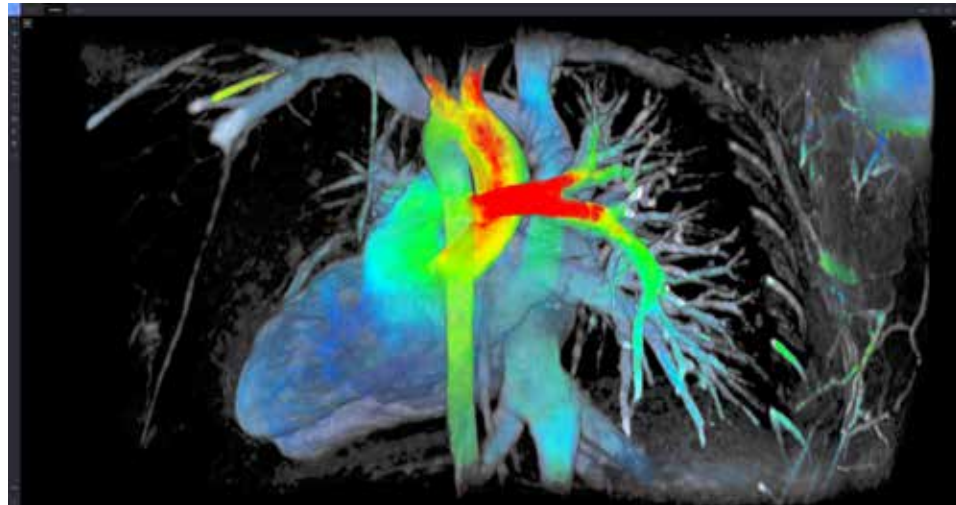






## Benefits


- Quantitative measurement of T1 and T2 relaxation times
- Automatic motion correction for T1 Mapping
- Two methods of acquisition for T1\* or true T1 measurements

## Achieving measurable benefits


CardioMaps is a powerful diagnostic technique that supports detection of cardiac pathologies by quantitative measurement of T1 and T2 relaxation times. The T1 Mapping acquisition includes automatic motion correction that compensates for cardiac and/or respiratory motion, providing reliable results. T1 Mapping offers two methods of acquisition: Inversion-recovery Look-Locker with FIESTA readout (MOLLI) for apparent T1 (T1\*) measurements or saturation-recovery SMART1Map for true T1 measurements.

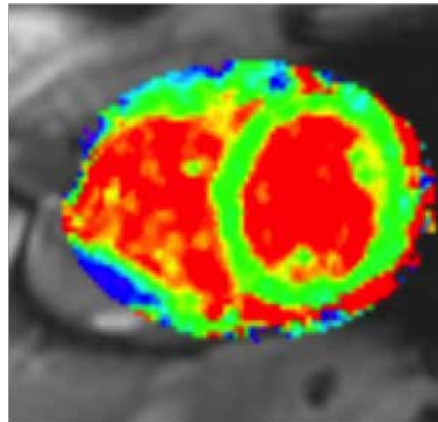


T1 CardioMap

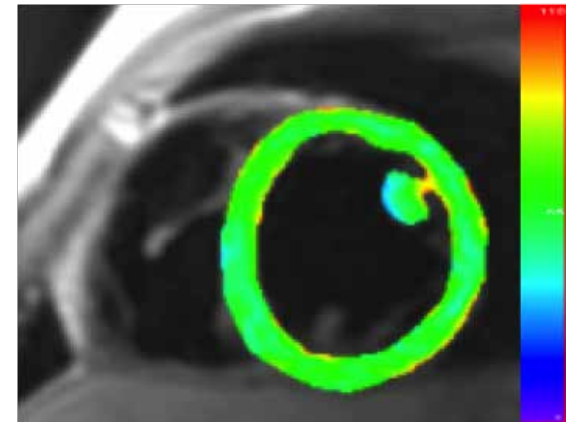

T2 CardioMap

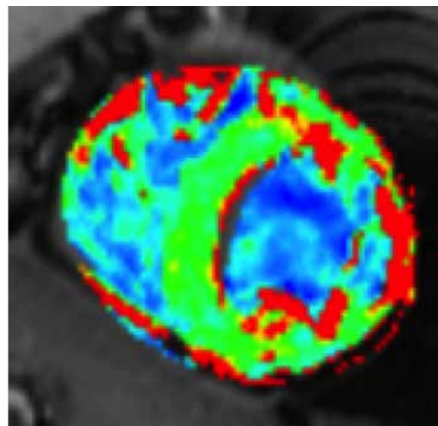

T1 CardioMap  
Without Motion Correction

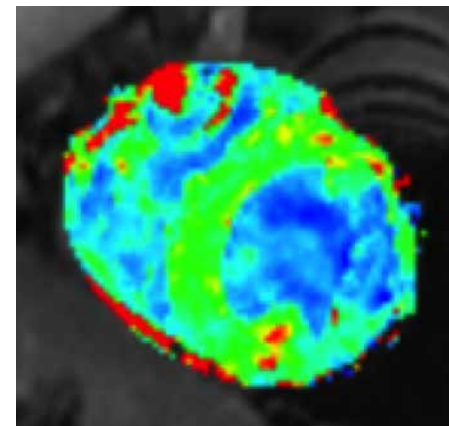

T1 CardioMap  
With Motion Correction




# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

Image  
Acquisition



# Pulse Sequences

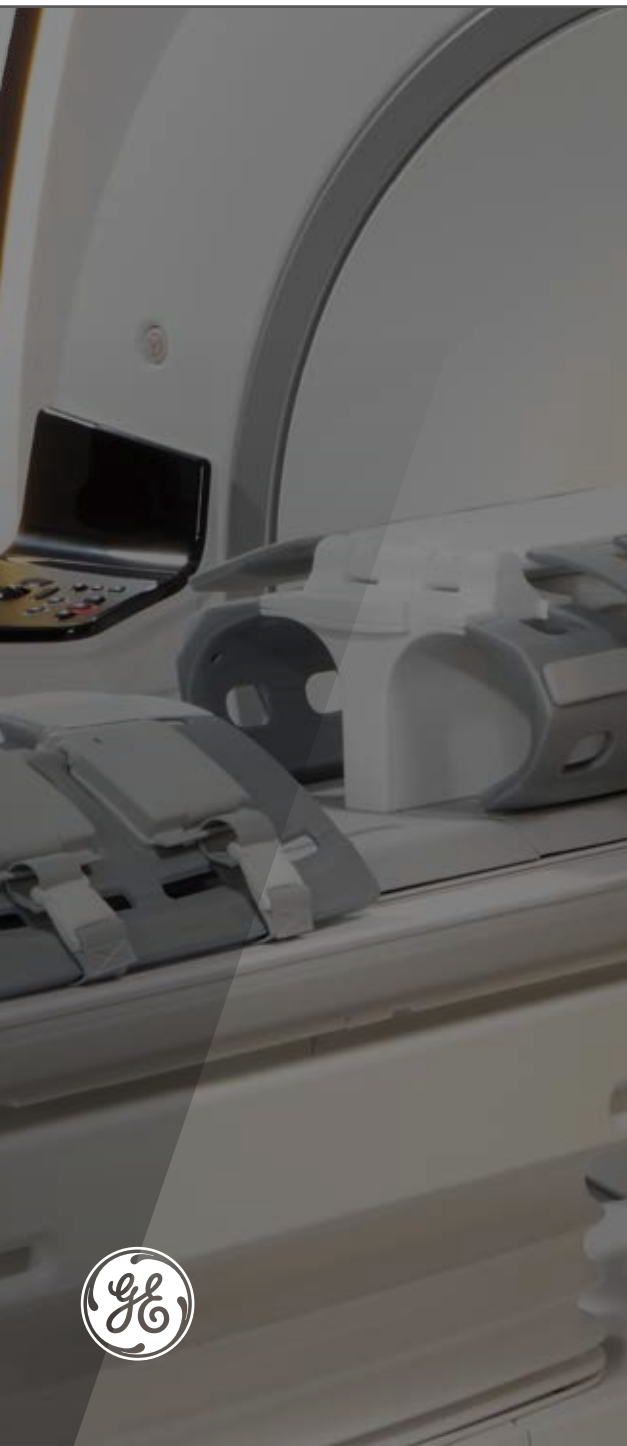


SPIN Echo	
SE FSE frFSE	Standard pulse sequences that are used to generate T1, Proton Density and T2 contrasts. The FSE technique enables long TR and long TE choices in reduced scan times. frFSE produces images with more T2 contribution allowing shorter TR values and resulting in shorter scan times when compared to FSE.
IR FSE-IR	IR techniques provide uniform suppression of tissues by applying an inversion pulse to null signal. FSE-IR reduces scan time while still achieving efficient tissue suppression.
3D FSE 3D frFSE	Three-dimensional imaging acquisitions mostly used for T2-weighted contrast.
T1 FLAIR T2 FLAIR	T1 and T2 Fluid Attenuated Inversion Recovery (FLAIR) pulse sequences allow the suppression of signal from cerebrospinal fluid (CSF). This sequence provides contrast to differentiate white and gray matter to T1- and T2-weighted brain and spine imaging.
Double IR/Triple IR (Black Blood)	These pulse sequences are included to allow Black Blood imaging for studies of cardiac morphology. Triple IR adds fat suppression to Black Blood imaging. It also can be combined with Single Shot.
Double IR/Triple IR Single Shot	Single Shot Black Blood acquisitions allow larger volume acquisitions in fewer breath-holds.
SSFSE SSFSE-IR	Single Shot Fast Spin Echo is a technique that permits single slice data acquisition in less than one second. It is frequently used for MRCP studies in a single breath hold and myelograms.
3D MRCP	3D frFSE sequence that combined with the T2 Prep option provides improved background tissue suppression for MRCP exams.
T2 MAP*	T2 MAP is a multiple acquisition; multiple echoes FSE based method to obtain images that represent different T2-weighting values. The acquired data is processed to produce T2 color maps that are used for cartilage evaluation.

\*Optional



# Pulse Sequences (continued)

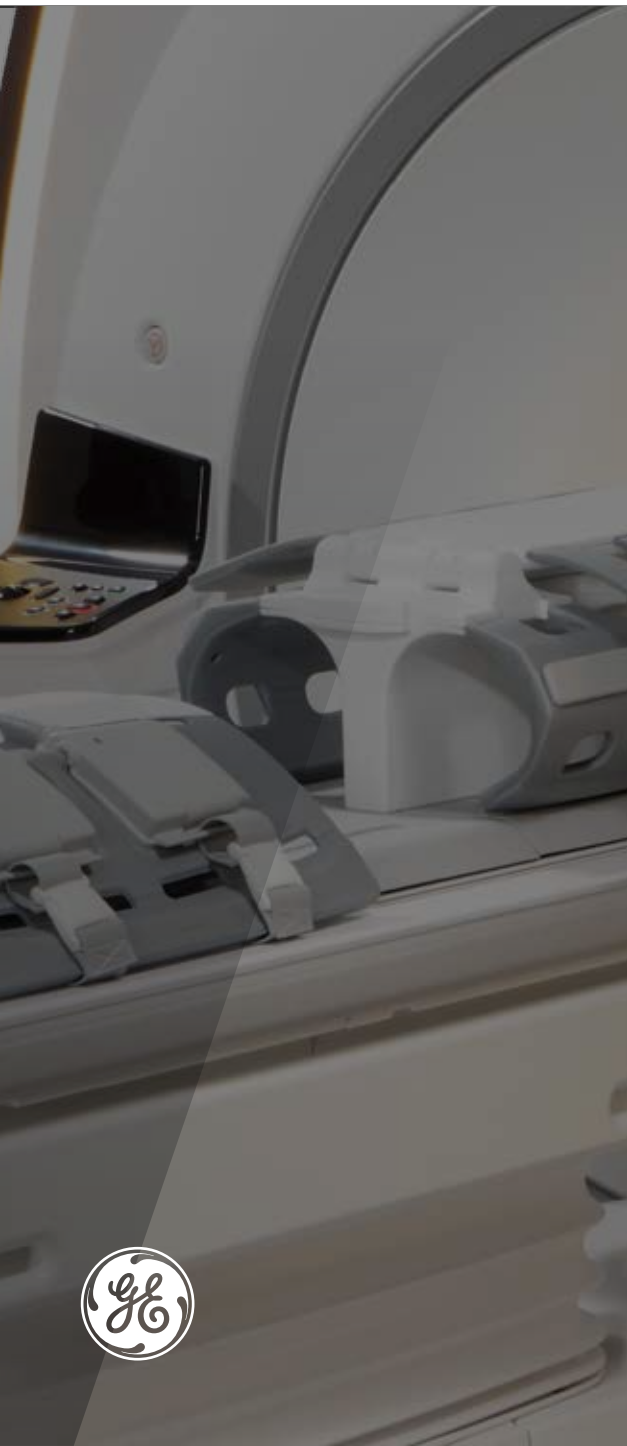

## SPIN Echo (Continued)

Cube Cube FLAIR	3D FSE technique that applies modified refocusing pulses for increased SNR. It is used to acquire isotropic data that can be reformatted in any plane.
Cube DIR	Cube DIR, double inversion recovery, is designed to achieve signal suppression from either gray or white matter and CSF.
Cube PROMO*	Prospective Motion correction is a real time 3D navigator based motion correction technique compatible with Cube T2, Cube DIR and Cube T2 FLAIR.
2D IDEAL*	2D FSE 3-point Dixon Water Fat Separation method that acquires 4 contrasts in one acquisition: Water, Fat, in-phase and out-of phase.
MAVRIC SL*	Multi-Spectral imaging technique designed to reduce metal artifact near MR conditional implants. The sequence provides T1, PD or STIR contrast.
3D ASL*	3D FSE based technique that uses a "labeling" pulse to quantify cerebral blood flow.

\*Optional

## Gradient Echo


2D and 3D GRE/SPGR 3D GRE Dual Echo 2D and 3D FGRE/FSPGR 2D MFGRE (Multi Echo) 2D CINE GRE / SPGR	Gradient Echo basic techniques offer a variety of possibilities to support imaging of all anatomies and can be acquired in 2D, 3D and Cine modes. The sequences generate T1 or T2 contrasts and support single, dual and multi echo acquisitions.
2D and 3D MDE	Myocardial Delayed Enhancement is a technique used for tissue characterization to provide the assessment of myocardial perfusion.
PSMDE	Phase Sensitive MDE increases the contrast between enhanced and normal tissue even with non-optimal inversion delay times.
SSMDE and SSPSMDE	MDE and PSMDE Single Shot based sequence that provides multi slice coverage with reduced breath hold times.
2D and 3D FIESTA 2D FIESTA CINE 2D FatSat FIESTA 3D FIESTA-C	Fast Imaging Employing Steady-state Acquisition generates great contrast differentiation between tissues of low T2/T1 ratios and high T2/T1 ratios. Provides high SNR images in short acquisition times. FIESTA sequences offer benefits for Neuro, Cardiac and Abdominal imaging.

\*Optional





# Pulse Sequences (continued)

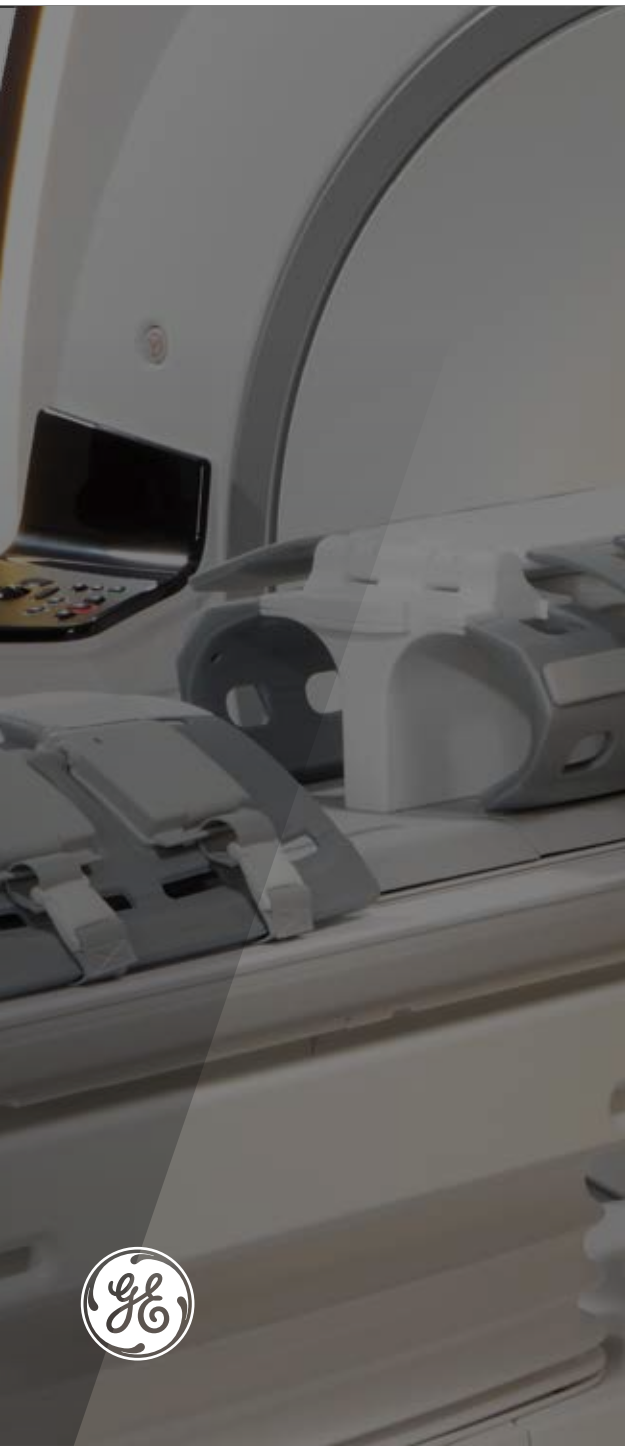


## Gradient Echo (Continued)

2D and 3D MERGE FGRE	T2* contrast technique that acquires multiple echoes at several different TE values.
2D Fastcard GRE/SPGR	Prospective gating sequence designed for breath hold, aortic arch gated imaging.
2D FastCINE GRE/SPGR	Retrospective gating sequence, beneficial to cardiac wall motion studies, assessment of valve function and visualization of regurgitation and stenosis.
2D FGRE-ET* 2D FGRE-ET Realtime*	Fast Gradient Echo sequence combined with an EPI echo train for acquiring multiple phase encoding steps per TR. Used for first pass myocardial perfusion studies. Compatible with Real time for cardiac planning and imaging uncooperative patients.
2D FGRE TC*	Fast Gradient Time Course used for myocardium tissue evaluation on first pass studies. Allows multiple planes radial acquisitions.
2D CINE-IR	FAST-CINE GRE IR Prep sequence is designed for myocardial viability studies. Supports TI time selection for consistent results.
2D Real time FGRE/FIESTA	Free breathing, real time planning sequence for whole heart coverage.
2D Tagging*	Fast Cine GRE based sequence for visualization of cardiac contractile function.
3D Heart*	3D FGRE/FIESTA navigated sequence for free breathing coronary artery imaging.
3D COSMIC	Coherent Oscillatory State Acquisition for the Manipulation of Imaging Contrast is a modified FGRE sequence with steady-state free precession segmented acquisition for high SNR, high contrast spine imaging.
3D LAVA	Liver Acquisition with Volume Acceleration is a 3D SPGR technique designed to image the liver. SPECIAL is the fat suppression method applied and parallel imaging provides shorter scan times.
3D LAVA Flex*	3D FSPGR technique that acquires in-phase, out-of-phase, water only and fat only images in one acquisition. LAVA Flex uses ARC; a self-calibrated 2D parallel imaging technique that allows acceleration in phase and slice direction.
3D Turbo LAVA 3D Turbo LAVA Flex*	LAVA Turbo provides a reduction of breath-hold timing for both LAVA and LAVA Flex acquisitions by as much as 20% reduction compared to conventional LAVA and LAVA Flex acquisitions.

\*Optional

# Pulse Sequences (continued)

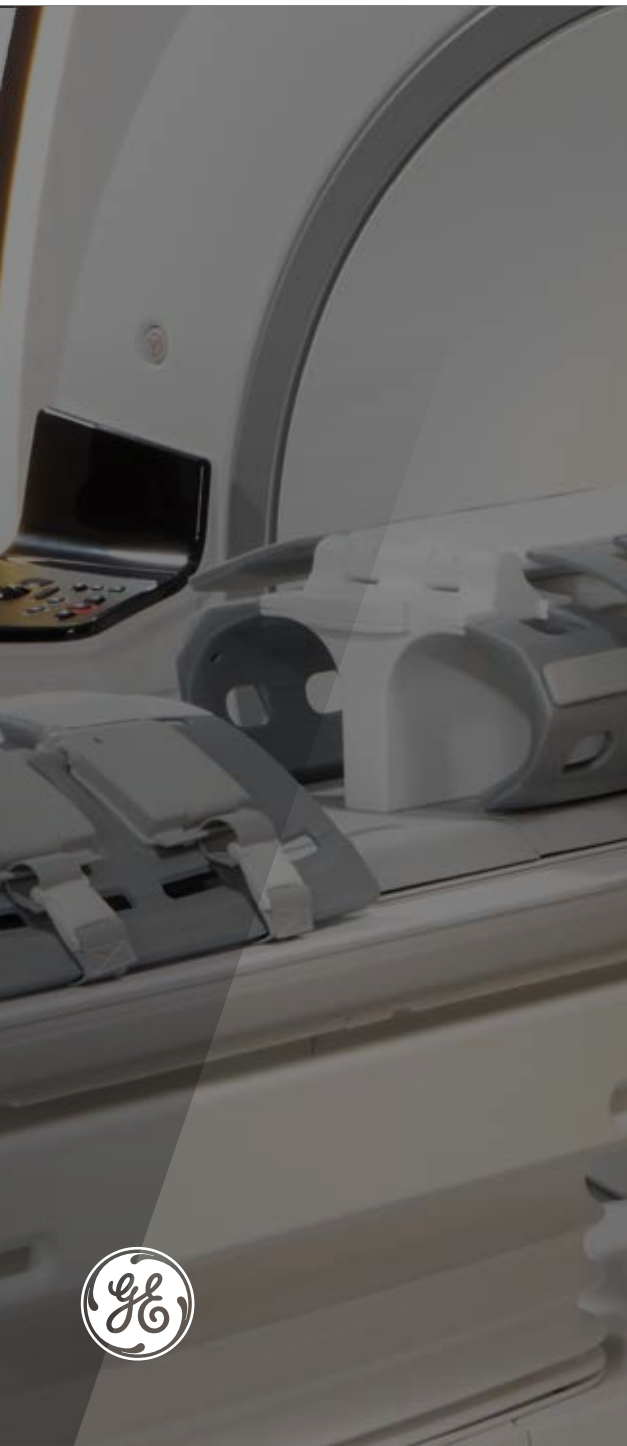


## Gradient Echo (Continued)

3D VIBRANT*	Simultaneous bilateral breast imaging technique in the axial and sagittal plane. SPECIAL and dual-shim volume capabilities provide homogeneous fat suppression.
3D VIBRANT Flex*	Acquires in-phase, out-of-phase, water only and fat only images in a single scan. It provides robust fat saturation and applies ARC, 2D Self-calibrated acceleration method for high spatial and temporal resolution images.
3D QuickSTEP	QuickStep is an automated multi-station run-off acquisition. This application automatically prescribes, acquires, and combines images from multiple stations for fast acquisition and simplified workflow.
3D TRICKS*	The Time Resolved Imaging of Contrast KineticS (TRICKS) is a fast 3D dynamic acquisition for high temporal and spatial resolution MR angiography imaging. Combined with elliptical-centric data sampling for consistent results.
3D SWAN*	High-resolution susceptibility weighting 3D multi echo gradient acquisition designed for small vessels visualization, as well as large vascular structures and iron or calcium deposits in the brain.
3D IDEAL*	IDEAL is a 3-point Dixon water fat separation method that generates in-phase, out-of-phase, water images and fat images in one single scan. Provides homogeneous fat saturation for imaging for challenging anatomies as such as neck and spine.
3D IDEAL-IQ*	Whole liver 3D coverage in a single breath-hold, IDEAL IQ provides a non-invasive, quantitative assessment of triglyceride fat content in the liver that can aid in diagnosing steatosis.
StarMap*	StarMap is an acquisition and post processing technique that helps evaluate iron content in the heart and liver. Multiple echoes are acquired at different TE times for each pixel resulting in images that represent variations of T2*-weighting. After the acquisition the images are post processed to generate color and grayscale T2* and R2* Maps.
DISCO*	Differential Sub-Sampling with Cartesian Ordering, combine TRICKS and LAVA Flex technologies to acquire high temporal resolution 3D dynamic images with robust fat suppression and without compromising spatial resolution.

\*Optional




# Pulse Sequences (continued)

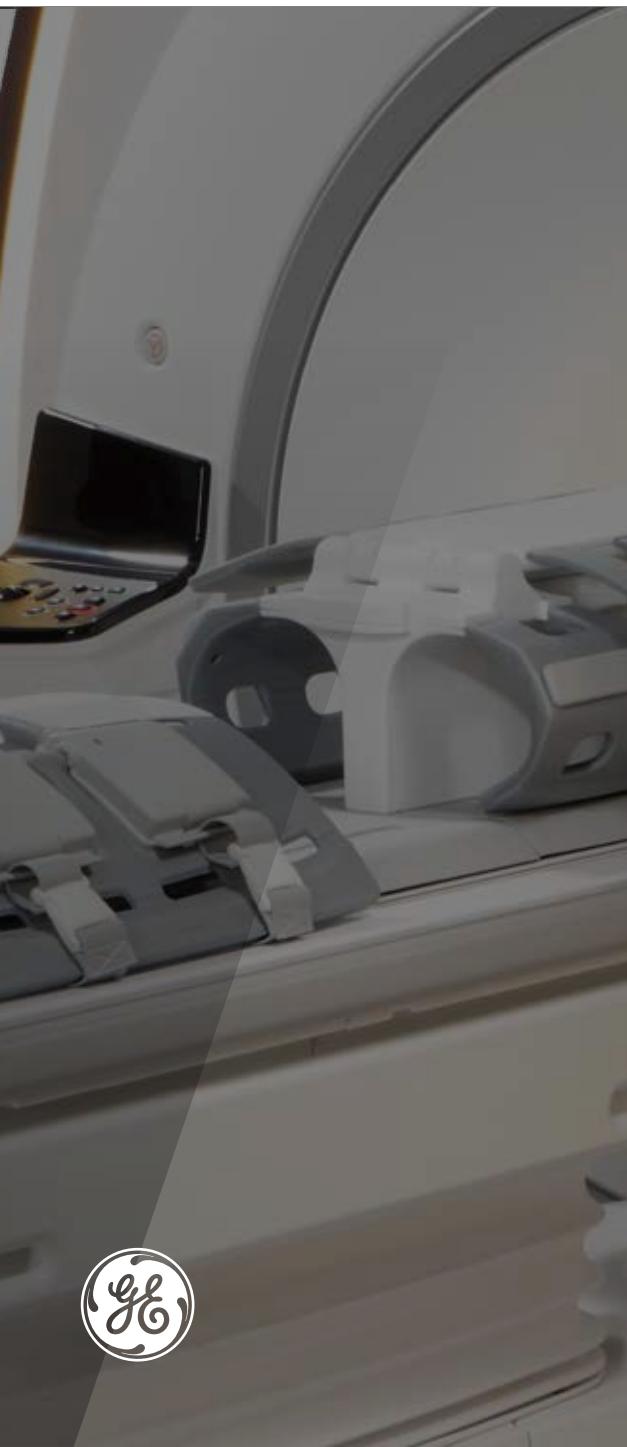

Vascular	
Inhance Inflow IR*	3D Fiesta based non-contrast-enhanced MR angiography technique that provides static background tissue and venous flow suppression for imaging arteries. It uses SPECIAL for uniform fat suppression and respiratory gating compatibility reduces respiratory motion artifacts during free-breathing renal exams.
Inhance 3D Velocity*	3D Phase Contrast based technique designed to acquire angiographic images in brain and renal arteries with robust background suppression in a short scan time. Respiratory triggering compatibility enabling abdominal angiography.
Inhance 2D Inflow*	Designed for imaging arteries that follow almost a straight path (i.e. femoral, popliteal, and carotid arteries) Inhance 2D Inflow acquires data during the systolic phase only. Compatible with Cardiac or Peripheral Gating and ASSET.
Inhance 3D Delta Flow*	3D FSE cardiac gated based non-contrast-enhanced MRA application designed for peripheral arterial imaging. This technique uses the differences between systolic and diastolic flow to help generate arterial signal contrast with robust background and venous suppression. ASSET compatibility provides shorter scan times.
2D TOF 2D Gated TOF 2D Fast TOF FGRE/SPGR 3D TOF 3D Fast TOF FGRE/SPGR	2D TOF Imaging, 2D Gated TOF Imaging, 3D TOF Imaging and Enhanced 3D TOF Imaging are used for MR angiography imaging. Based on conventional gradient echo scanning, TOF imaging techniques rely primarily on flow-related enhancements to distinguish moving from stationary spins.
2D CINE Phase Contrast 2D FastCINE Phase Contrast	This pulse sequence is included specifically for studies of cardiac function. Through the use of retrospective gating, it allows full R-R coverage.
2D Phase Contrast 3D Phase Contrast	These techniques demonstrate flow velocities and directional properties in vessels and other moving fluids such as CSF and aortic flow.

\*Optional

# Pulse Sequences (continued)

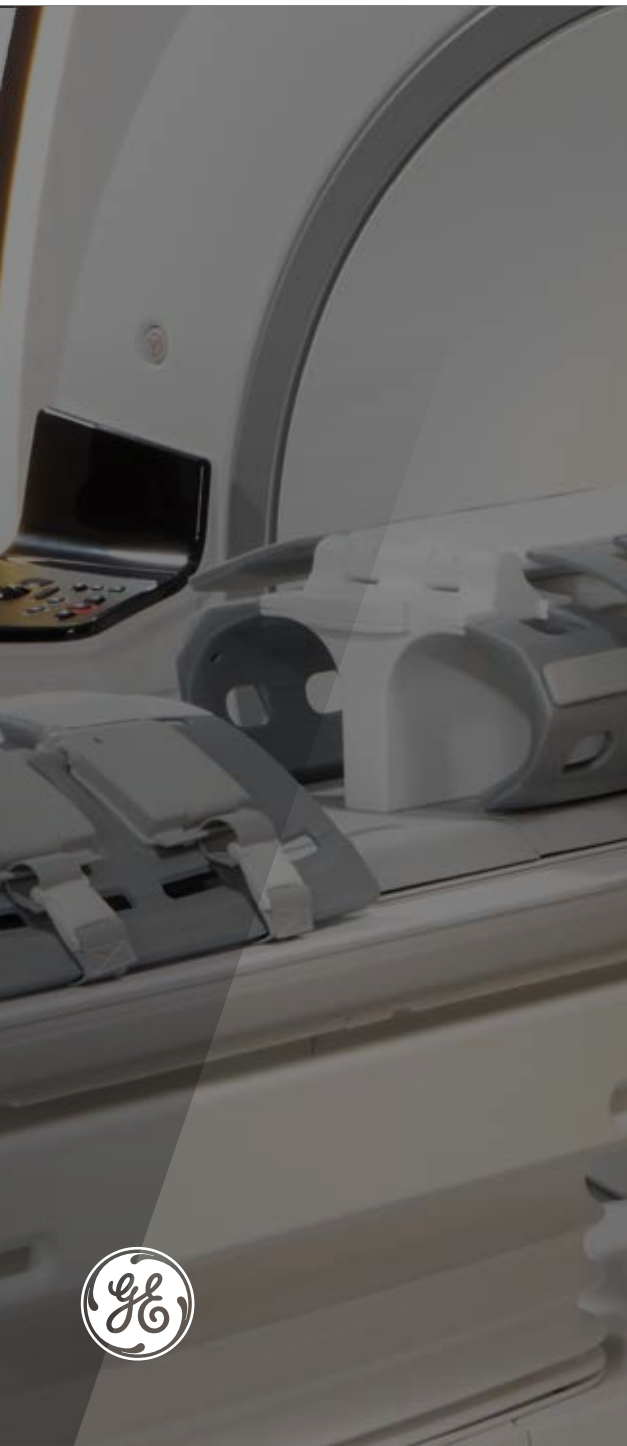


EPI	
fMRI – BrainWave*	Real time acquisition, processing and display of functional imaging.
GRE-EPI SE-EPI FLAIR-EPI DW-EPI	Standard on all systems are gradient echo, spin echo, flair, and diffusion weighted echo planar imaging. The EPI sequence supports single and multishot imaging, multi-phase imaging, as well as cardiac gating. Diffusion EPI produces images that can detect acute and hyper-acute stroke with b-value up to 10,000 s/mm <sup>2</sup> , multi-NEX compatibility and the ability to generate ADC and T2-weighted TRACE images. The FLAIR option suppresses the CSF signal.
DTI*	DTI (Diffusion Tensor Imaging) is an EPI technique that acquires diffusion information in up to 150 different directions. The image contrast is based on the degree of diffusion anisotropy in the tissues. Post processing include Fractional Anisotropy (FA), Apparent Diffusion Coefficient (ADC), 2D directional maps and 3D fiber track models.
eDWI	Enhanced DWI (eDWI) provides high SNR diffusion images with short acquisition times. Supports Multi b-values with SMART NEX for variable NEX selection per B-value, “3 in 1” diffusion weighting to all three gradients simultaneously, tetrahedral selection with four different diffusion weighting combinations for shorter TE values and inversion recovery for fat signal reduction.
RTFA	The RTFA algorithm leads to a reduction in distortion of the diffusion image per diffusion axis. RTFA is designed to reduce image blurring and distortions typically associated with diffusion imaging throughout the body. RTFA also allows for increased utilization of single spin echo DWI which results in an increase in SNR by up to 50% compared to dual spin echo and, when combined with the improved resolution leads to an increase in image quality that can be utilized for image presentation, fusion and ADC map outputs.
RTCF	Real Time Center Frequency (RTCF) option can be applied to DWI & DTI to enable using the optimal center frequency for each slice. This is intended to help improve fat suppression and signal drop off at areas of high B <sub>0</sub> inhomogeneity (off-center, or area with high tissue susceptibility). It is also intended to reduce station-to-station misalignment in whole body diffusion imaging.
FOCUS DWI*	FOV Optimized & Constrained Undistorted Single-shot (FOCUS) DWI utilizes 2D selective excitation pulses to limit the prescribed phase encode FOV eliminating artifacts from motion, imaging back folding or unsuppressed tissue.
MR-Touch*	MR-Touch is software and hardware application designed to measure relative tissue stiffness with MR. The acquisition uses a EPI based sequence that synchronizes induced vibrations to acquire a series of phase-contrast images over time.

\*Optional



# Pulse Sequences (continued)

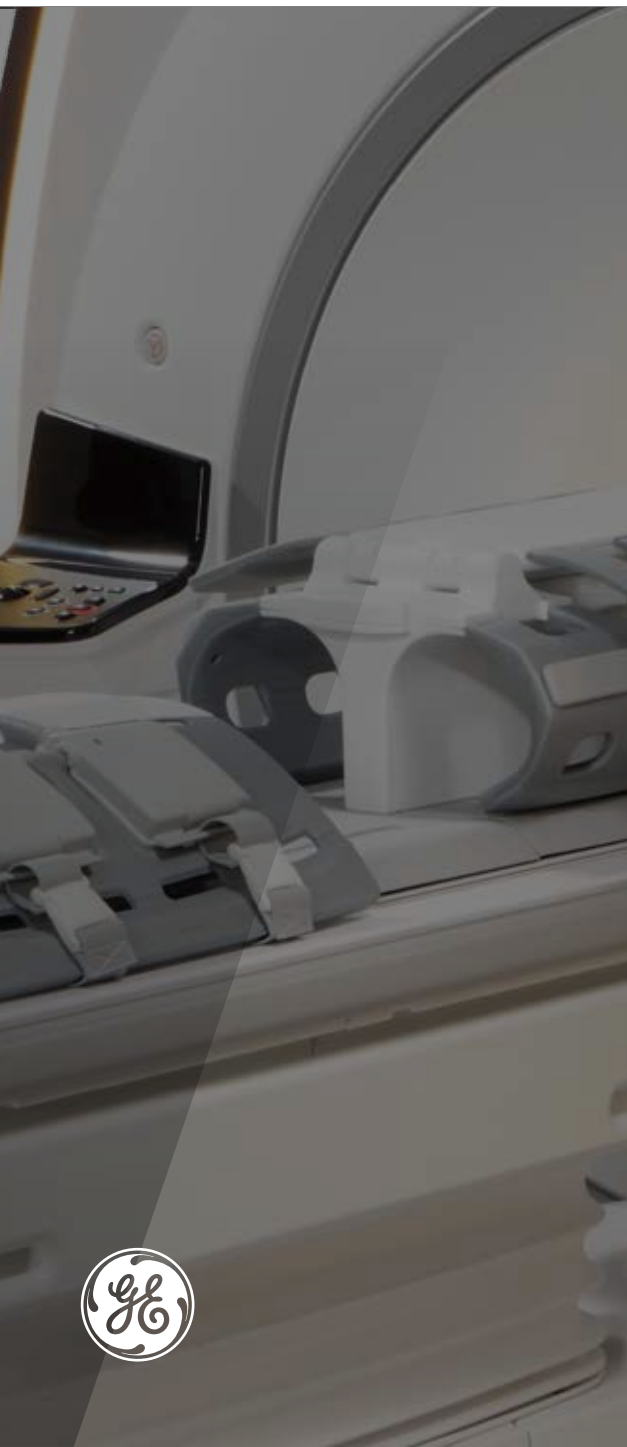

## Spectroscopy

PROBE-PRESS PROBE-STEAM*	PROBE Single-Voxel spectroscopy allows non-invasive evaluation of the relative concentrations of in-vivo metabolites. The sequence provides acquisition and display of volume localized, water-suppressed H1 spectra in single-voxel mode. The sequence consists of three slice selective RF pulses with crusher gradients. PRESS provides up to twice the SNR over STEAM.
PROBE-PRESS CSI 2D & 3D*	PROBE 2D and 3D CSI enable simultaneous multi-voxel spectroscopic acquisitions in the brain. It is available with PRESS excitation to maximize SNR. Post processing includes automatically generated metabolic maps.
BREASE*	A TE-averaged PRESS (Point RESolved Spectroscopy) acquisition that provides the necessary biochemical information to help characterize breast tissue by assessing the presence of choline.
MNS *	The Multi Nuclear spectroscopy option consists of software and hardware necessary to acquire, display and store 31P and 13C signals. The MNS pulse sequences are FID CSI, Echo CSI and Spin Echo (MRS). Sage 7 is the software tool that allows MNS data processing, display and storage. (3.0T only)
TEA-PRESS*	TEA PRESS is a TE-Averaged variant of the PRESS CSI pulse sequence. It collects spectra across a range of TE values and averages the results together to reduce the appearance of signals whose intensity varies as a function of TE. This allows signals whose intensity does not vary with TE to be accentuated in comparison. This is the underlying pulse sequence behind the BREASE application.

\*Optional


## PROPELLER MB

Silent T1, PD, T2, DWI, T1 FLAIR and T2 FLAIR PROPELLER MB*	PROPELLER MB is a multi shot per blade sequence that uses a radial k-space filling pattern acquisition and a post processing correction algorithm to significantly reduce the effects of motion artifacts. PROPELLER MB is compatible with spatial and chemical Sat, ASPIR, STIR T1, PD and T2 Auto TI/TR and Navigator.
T1, PD and T2 PROPELLER MB	
T2 FLAIR PROPELLER MB	
T1 FLAIR PROPELLER MB	
DWI PROPELLER MB	

\*Optional



# Pulse Sequences (continued)

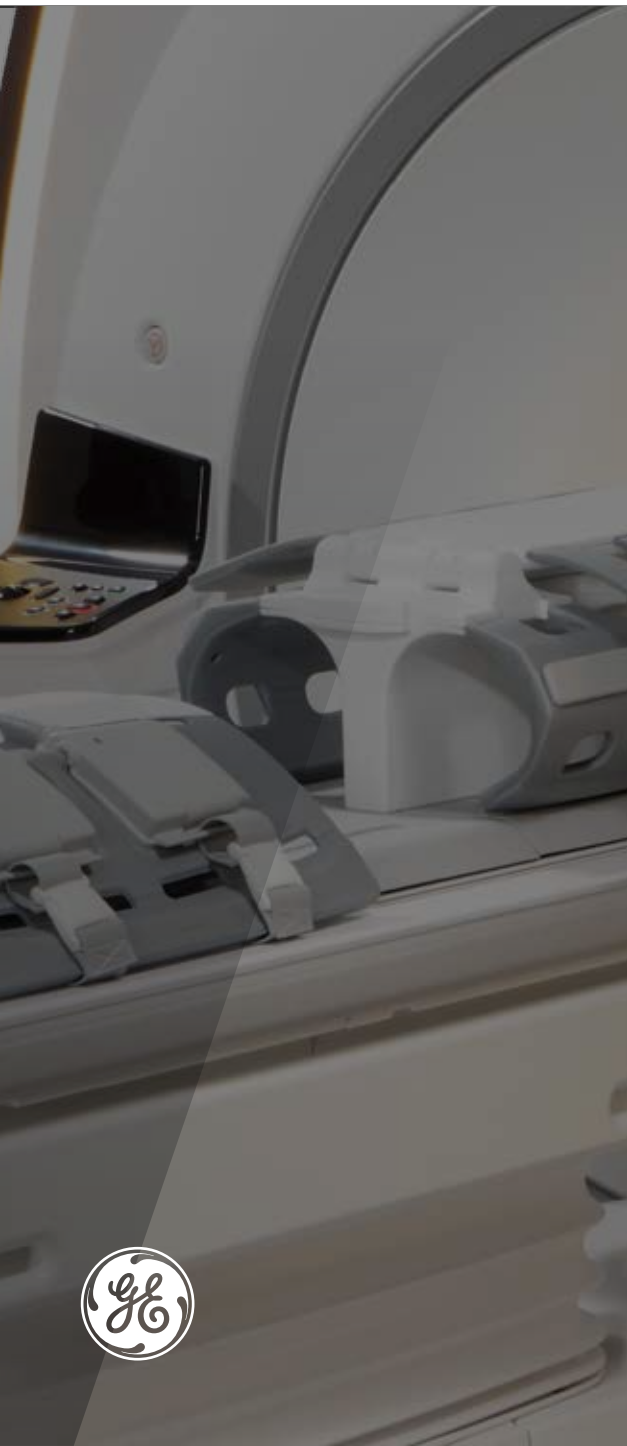

## Silenz\*

Silenz T1  
Silenz PD  
Silenz MRA

Silenz is a 3D Zero-TE sequence comprising high bandwidth excitation and reduced gradient-switching radial acquisition that results in sound levels that are within 3-4dB of ambient level. Silenz has added flexibility in sequence prescription for anisotropic resolution enabling faster scan times and includes axial as well as oblique geometries.

\*Optional

## Fat Suppression Technology

FatSat

Applies a frequency selective saturation pulse at the frequency of fat before the imaging excitation pulse with the result being a signal measurement primarily from water.

STIR

STIR is an inversion recovery method that takes advantage of the T1 difference between water and fat to allow selection of the signal to suppress. In order to eliminate the signal from tissues, the TI time must match exactly the null point of the tissue that needs to be suppressed.

SPECIAL

Hybrid fat suppression technique that incorporates features from both the frequency selective FatSat and the STIR techniques by using a spectrally selective inversion pulse that inverts only the fat magnetization and leaves the only the water peak available for excitation.

Spectral Spatial

Method that applies selective pulses for water excitation only, while fat is left untouched, thereby producing no signal.

ASPIR

ASPIR method is a solution for poor fat suppression due to  $B_1$  inhomogeneity. It is based on the frequency and the relaxation fat behaviors. Applies a spectrally selective adiabatic inversion pulse to excite the fat spins, imaging pulses are then applied after TI null time when longitudinal magnetization of fat crosses zero. The disadvantages include sensitivity to  $B_0$  and longer scan times.

IDEAL\*

IDEAL is a 3-point Dixon technique that acquires three images at slightly different echo times to generate phase shifts between water and fat. The water/fat separation method is very efficient at providing homogeneous image quality. One acquisition provides four contrasts: water, fat, in-phase and out-of-phase images.


Flex\*

Flex is a 2-point dixon technique delivering faster scan times compared to IDEAL 3-point dixon. It is based on the difference between fat and water resonance frequencies using two flexible echo times for further scan time reduction. One acquisition provides four contrasts: Water, Fat, in-phase and out-of-phase images.

\*Optional



# Pulse Sequences (continued)

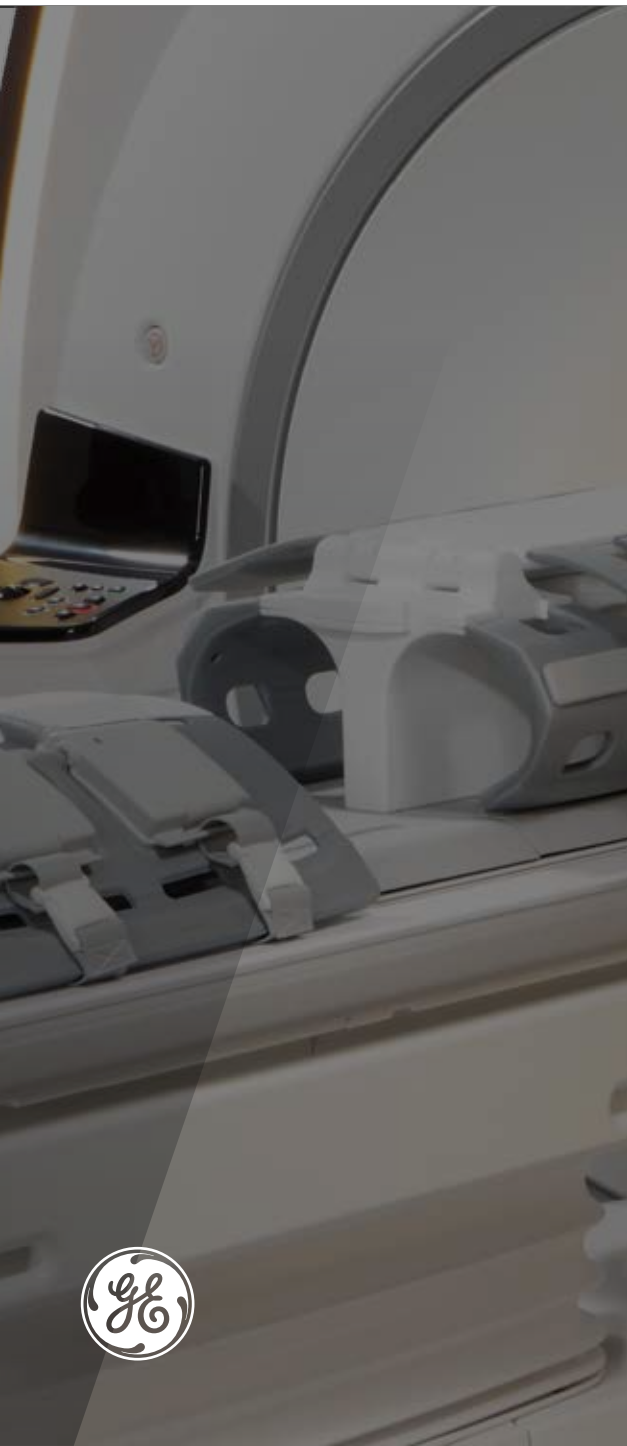


## Motion Correction Technology

PROPELLER MB	PROPELLER MB is a multi shot per blade sequence that uses a radial k-space filling pattern acquisition and a post processing correction algorithm to significantly reduce the effects of motion artifacts. It is compatible with spatial and chemical Sat, ASPIR, STIR Auto TI/TR and navigator.
PROMO*	Prospective motion correction is a real time 3D navigator based motion correction technique compatible with Cube T2, Cube DIR and Cube T2 FLAIR.
PB Navigators	Pencil beam navigators allow free breathing body and cardiac imaging by tracking the motion of the diaphragm. There are two navigator modes: navigator gating, uses a predefined signal acceptable range during the expiration and navigator triggering, uses signal to trigger data collection during the expiration.
Respiratory Trigger	Reduces breathing motion artifacts by synchronizing the acquisition with the respiratory cycle.
VCG	Vector cardiac gating reduces motion artifacts by synchronizing the acquisition with the cardiac cycle.
PG	Peripheral gating reduces motion artifacts caused by pulsating blood.

\*Optional



# Pulse Sequences (continued)

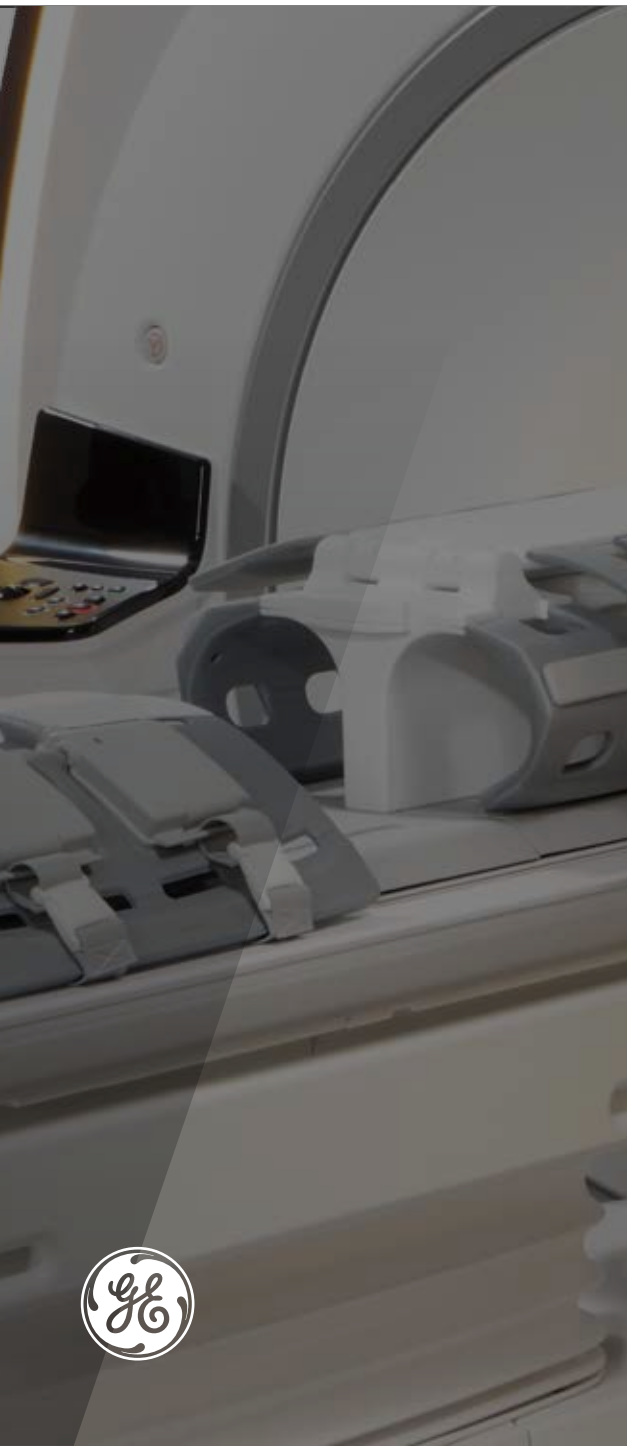


Acceleration Technology	
Fractional Nex	Technique in which only partial k-space data is collected and the remaining data is estimated. It uses the phase conjugate symmetry reconstruction method, which only half of the phase encode steps are acquired for scan time reduction.
ASSET	Array Spatial Sensitivity Encoding Technique acquires under sampled multicoil data generating aliased images. These are post processed with coil sensitivity maps from the calibration scan to unfold the images.
ARC	Auto-Calibrating Reconstruction for Cartesian Imaging is a highly accelerated parallel imaging auto-calibrating method that doesn't require coil sensitivity maps. It enables smaller FOV prescriptions, less sensitivity to motion and prevents artifacts caused by coil calibration inaccuracies.
HyperSense	High performance acceleration based on sparse or compressible images. It can be extended to include inherent compressibility in dimensions besides k-space. While parallel imaging suffers from SNR loss due to scan time reduction and coil spatial encoding, with HyperSense there is no SNR loss caused by the coil geometric factor.
Hyperkat*	HyperKat s an advanced k-t acceleration method that employs time-shifted sampling in data acquisition and exploits both spatial and temporal correlation with motion-adaptive time window selection in image reconstruction.
HyperCube*	Small FOV organ specific volumetric imaging acquisition method that enables outside phase FOV HyperCube* signal suppression. The technique can help to reduce artifacts originated outside of the prescribed field of view.

\*Optional





# Pulse Sequences (continued)

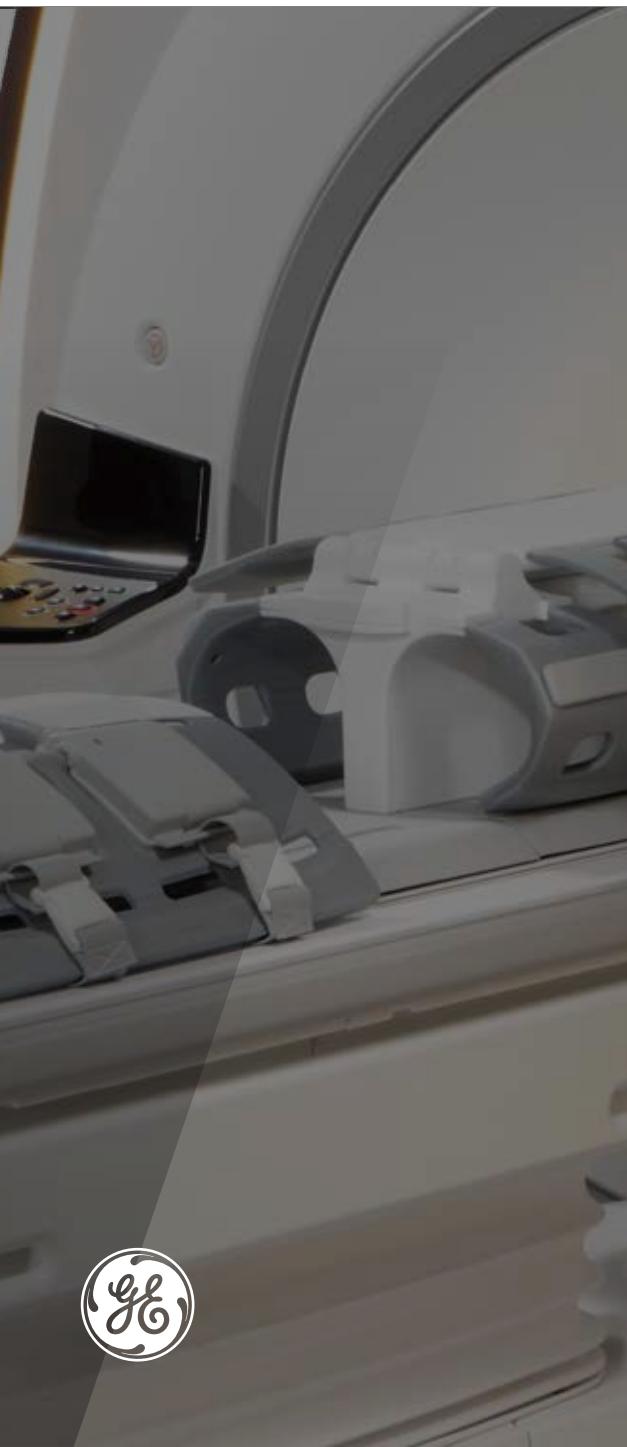

## Uniformity Correction Technology

SCIC	SCIC is a method that uses statistics from the actual image for inhomogeneity correction. It uses an optimized set of parameters for tuning that are defined according to the different anatomies and coils.
PURE	PURE corrects the field inhomogeneity by collecting a calibration scan from the (uniform) body coil and the (non-uniform) surface coil and calculating maps that relate the intensity correction values to the images.
deFINE	deFINE is an integrated in-line imaging processing method that provides edge enhancement and smoothing algorithms allowing the user to customize the image appearance.
reFINE	reFINE is an advanced image uniformity correction algorithm that addresses non-uniformity due to coil sensitivity profiles and dielectric shading effects. It reduces organ-motion induced misregistration artifacts, effects of low signal in dark regions and edge effects at tissue interfaces and borders. reFINE optimizes parameter settings for each application, coil, and body anatomy maximizing image uniformity results.

\*Optional

## Noise Reduction Technology

ART	Acoustic Noise Reduction Technology optimizes the gradient waveform to reduce the gradient noise without compromising performance.
Silenz*	Silenz is a 3D Zero-TE sequence comprising high bandwidth excitation and reduced gradient-switching radial acquisition that results in sound levels that are within 3-4dB of ambient level. Silenz has added flexibility in sequence prescription for anisotropic resolution enabling faster scan times and includes axial as well as oblique geometries.
Silent PROPELLER*	Silent PROPELLER gradient waveform approach reduces the acoustic noise level to less than 11dB above the ambient room noise.

\*Optional



GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite


RF Coils Suite



# eXpress Table & Posterior Array

The GEM Suite, Geometry Embracing Method, consists of a set of receive-only RF arrays designed to reduce multiple physical coil changes within a single exam and between different exams.

- Detachable table with embedded Posterior Array
- 100cm S/I Coverage
- 40 Elements with dedicated spine configurations
- Head-first or Feet-first
- Automatic coil mode selection
- Acceleration in All Directions

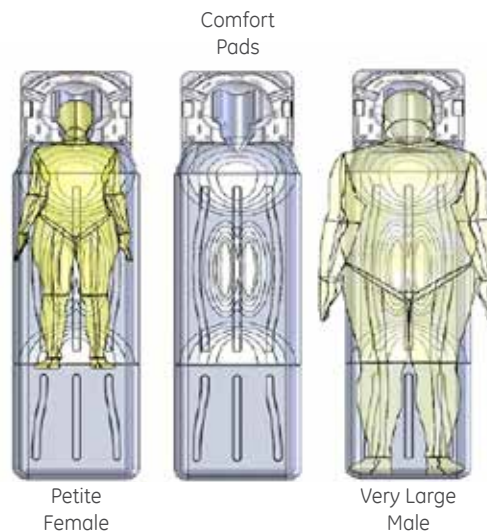

# Head & Neck Unit

The GEM HNU is a standard component of the GEM Suite. The HNU consists of four imaging components: a head base-plate, an anterior neuro-vascular face array, and the open face adapter. The coil may be positioned at either end of the GEM table to support head-first or feet-first imaging.




Head and Neck NV with Comfort Tilt



Head and Neck



Cervical Open Face

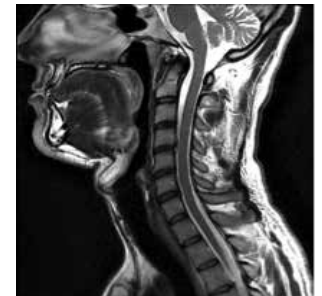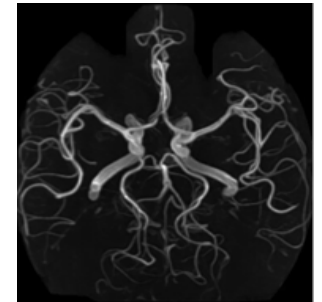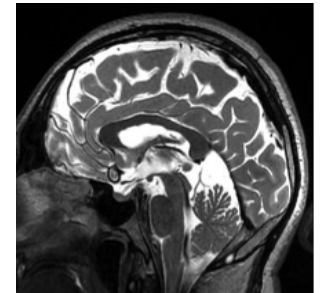

Head Neck Unit NV Specifications

Length	49.5 cm (19.5 in)
Width	38.8 cm (15.3 in)
Height	35.4 cm (13.9 in)
Weight of HNU base	5.2 kg (11.4 lb)
Weight of Anterior Adapter	2.7 kg (5.9 lb)
S/I Coverage	50 cm (19.7 in), when combined with the PA and AA
R/L Coverage in head mode	24 cm (9.4 in)
R/L Coverage for NV	50 cm (19.7 in), when combined with the PA and AA
Head-first or feet-first imaging	
Up to 28 elements in the FOV, when combined with the PA and AA	



# Head & Neck Unit (continued)



## Head Neck Unit Cervical Specifications

Length	49.5 cm (19.5 in)
Width	38.8 cm (15.3 in)
Height	32.6 cm (12.8 in)
Weight of Cervical Adapter	1.7 kg (3.7 lbs)
S/I Coverage	28 cm (11 in)
R/L Coverage	24 cm (9.4 in)
Head-first or feet-first imaging	
Up to 14 elements in the FOV, when combined with the PA	

## Head Neck Unit with Open Face Adapter Specifications

Length	49.5 cm (19.5 in)
Width	38.8 cm (15.3 in)
Height	25.9 cm (10.2 in)
Weight of Open Face Adapter	1.3 kg (2.8 lbs)
S/I Coverage	28 cm (11 in)
R/L Coverage	24 cm (9.4 in)
Head-first or feet-first imaging	
Up to 12 elements in the FOV, when combined with the PA	








# Anterior Array

The GEM AA is a standard component of the GEM Suite that facilitates chest, abdomen, pelvis, and cardiac imaging with the GEM RF coil suite. The GEM AA is lightweight, flexible, thin and pre-formed to conform to the patient's size and shape. With 54 cm of S/I coverage, the coil permits upper abdominal and pelvic imaging without repositioning the patient. The 16 element electrical design supports parallel imaging in all 3 planes.

## Anterior Array Specifications

Length	55.6 cm (21.9 in)
Width	67.4 cm (26.5 in)
Height	3.3 cm (1.3 in)
Weight	2.8 kg (6.16 lb) resting on patient 3.9 kg (8.6 lb) with cable
S/I Coverage	54 cm (21.3 in)
R/L Coverage	Full 50 cm (19.7 in) FOV of the system
Head-first or feet-first imaging	
Up to 36 elements in the FOV, when combined with the PA	

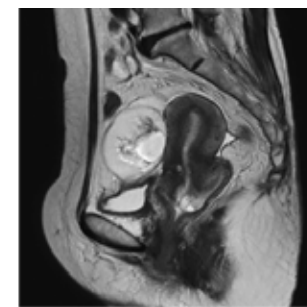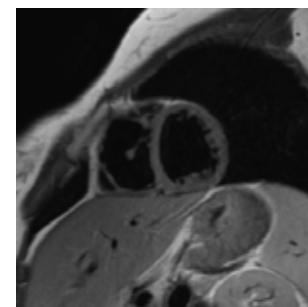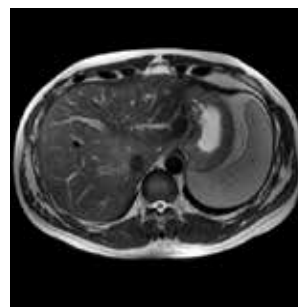



# Peripheral Vascular Array

The GEM PVA is an optional component of the GEM Suite that facilitates imaging of the thighs and lower legs. The high-density layout supports parallel imaging in all 3 planes. The coil incorporates an innovative hinge design between the upper & lower elements to simplify patient setup. In addition, to improve patient comfort, the lower leg section of the coil is fully supported by the GEM table and not the patient.

## Optional Peripheral Vascular/Lower Extremity Array

Length	105 cm (41.3 in)
Width	2nd station: 64.2 cm (25.3 in) 3rd station: 51.6 cm (20.3 in)
Height	24.8 cm (9.8 in)
Weight	9.1 kg (20.0 lbs)
S/I Coverage	104 cm (49.9 in) overall 2nd station: 52 cm (20.5 in) 3rd station: 52 cm (20.5 in)
R/L Coverage	Full 50 cm (19.7 in) FOV of the system
Head-first or feet-first imaging	
Up to 35 elements in the FOV, when combined with the PA	






\*Optional



# 16 channel Flex Coils



Knee and foot ankle positioner

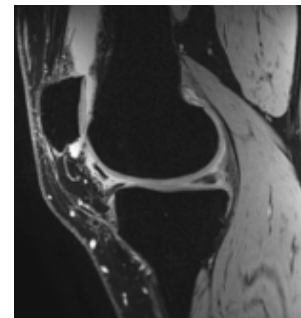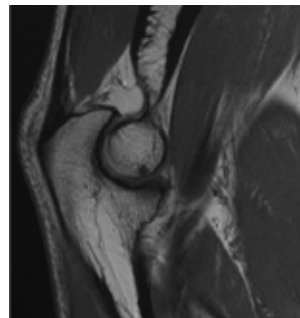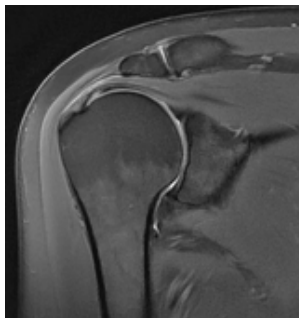

## GEM Flex Specifications

Coil	Dimensions (W x L x H)	Wrap Diameter	Elements	Weight
GEM Flex Large	23 cm x 71 cm x 5 cm	155 mm – 215 mm	16	1.2 kg
GEM Flex Medium	23 cm x 57 cm x 4 cm	115 mm – 155 mm	16	0.9 kg
GEM Flex Small	23 cm x 44 cm x 4 cm	90 mm – 125 mm	16	0.9 kg


# 18 channel T/R Knee Coil

## Benefits

- Transmit Receive 18 channel array design
- Large diameter to better accommodate anatomy
- High SNR for unique performance
- Parallel imaging compatible for speed



The 18 channel Transmit Receive Phased Array Knee coil is designed to acquire high SNR images of the knee. It is generously sized to effortlessly accommodate a wide range of the patient population. The two-part design provides a quick and efficient workflow. Offset imaging is fully supported with adjustable left-right coil positioning.

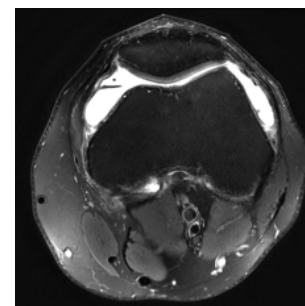

## Specifications

Coil	Approximate dimensions (W x L x H)	Approximate Diameter	Elements	Approximate Weight
18ch T/R Knee	544 mm x 500 mm x 282 mm	155 mm	18	6.8 kg





# 16 channel Shoulder Coil

## Benefits

- 16 channel phased array design
- Adjustable anterior plate for ease of positioning
- Parallel imaging compatible for speed



The Phased Array 16 channel Shoulder coil consists of a baseplate that supports a posterior hard shell connected to an adjustable anterior plate, designed to better accommodate the patient anatomy. The baseplate and customized pad provide easy right-left adjustment for off center positioning.

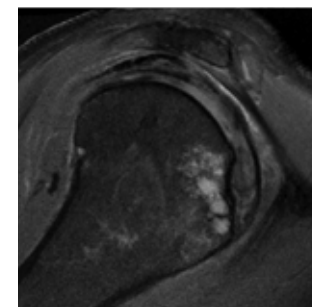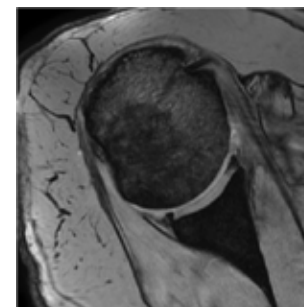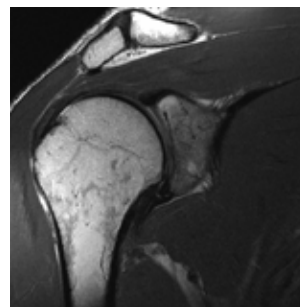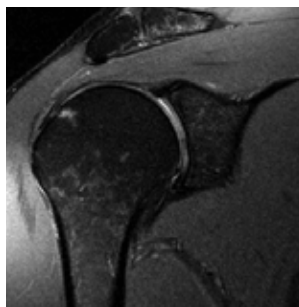

## Specifications

Coil	Approximate dimensions (W x L x H)	Elements	Approximate Weight
16ch Shoulder	280 mm x 280 mm x 311 mm	16	3.9 kg






# 8 channel Foot and Ankle Coil

## Benefits

- 8 channel dedicated foot and ankle phased array coil
- Optimized design to accommodate foot and ankle anatomy
- Slide and lock mechanism for easy positioning



The Phased Array 8 channel Foot and Ankle coil consists of a baseplate and a detachable hard shell coil that is designed for fast and easy positioning, comfortably accommodating the anatomy while providing proper immobilization.

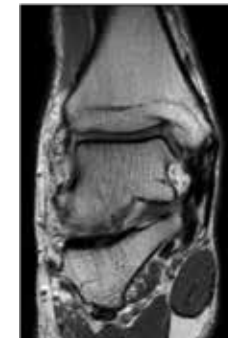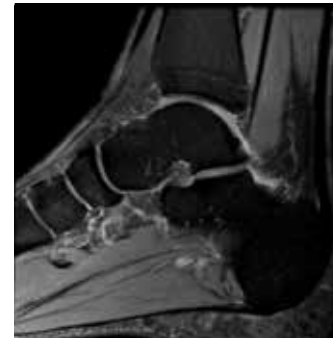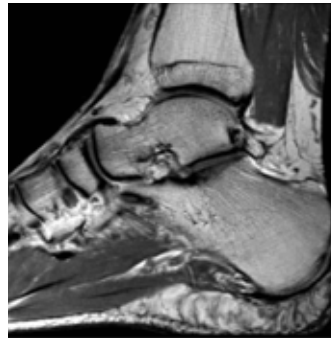

## Specifications

Coil	Approximate dimensions (W x L x H)	Elements	Approximate Weight
8ch Foot and Ankle	180 mm x 337 mm x 314 mm	8	3.1 kg
Baseplate	358 mm x 515 mm x 336 mm	–	3.8 kg




# 16 channel Sentinelle Breast Coil

## Benefits

- Variable Coil Geometry optimized for parallel imaging
- Ease access for breast biopsy procedures
- Visco-elastic padding for improved patient comfort
- High SNR, lower scan time and uniform fat suppression



The 16 channel Sentinelle Breast tabletop supports different coil configurations to provide optimal performance for high resolution bilateral, unilateral and biopsy mode breast imaging.

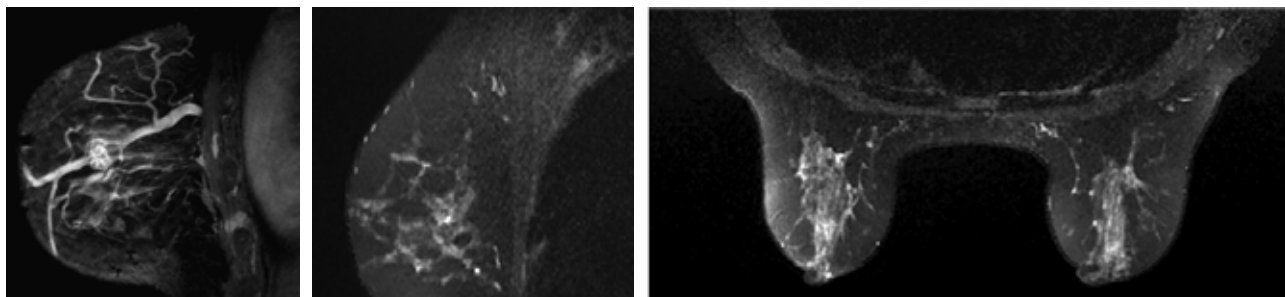

# RF Coils and Arrays

There are many optional receiver coils available to configure a Discovery™ MR750w 3.0T system to meet specific applications requirements. The coils listed below are commercially available at the time of printing and are optional with the system. Please contact your local GE sales representative for the most current list.



## 32-Channel Brain

- 32-channel array design for high signal-to-noise ratio (SNR) and excellent contrast to noise ratio (CNR)
- Optimized for fMRI studies
- Parallel imaging compatible in all three planes
- Sliding coil design for ease of patient positioning
- Coil dimensions
- Inner diameter: 24 cm (9.5 in)
- Coil size: 29cm OD x 30cm long (11.4 x 11.8 in) Base size: 41 cm x 48 cm x 19 in)



## HD Brain Array

- 8-channel, 8-element phased-array design
- Parallel imaging compatible
- Coil dimensions: 69 x 38 x 41 cm (27 x 15 x 16 in) 24 cm S/I coverage



## Small Anterior Array

- Up to 33 elements in the FOV when combined with PA
- Head-first or feet-first
- Optimized for parallel imaging
- Anterior coil dimensions (L x W x H): 45 cm (17.7 in) x 40.5cm (15.9 in) x 4.5 cm (1.8 in)
- Anterior coil weight: 2.95 kg (6.5 lbs)



# RF Coils and Arrays (continued)



## HD Breast Array

- 8-channel, 8-element phased-array design
- Optimized for uniformity, parallel imaging and VIBRANT
- Bilateral and unilateral breast imaging
- Biopsy plates available
- Coil dimensions: 50 x 54 x 25 cm (20 x 21 x 10 in)



## Shoulder Phased Array

- 3-channel phased-array coil
- Sleeve design
- Comprehensive shoulder imaging

Homogeneous penetration of the humeral head and neck, rotator cuff, glenoid labrum, acromium process, and glenohumeral articular surfaces



## HD Wrist Array Coil

- 8-channel, 8-element phased-array coil
- High SNR to enable high spatial resolution images
- Position overhead or at patient's side
- Coil dimensions: 34 x 23 x 14 cm (13 x 9 x 6 in) (includes base)
- 12 cm S/I coverage for wrist and hand



# RF Coils and Arrays (continued)



## Endo-Rectal

- 1 Channel, inflatable and auto-tuned
- Small FOV and high spatial resolution, sensitivity, and specificity for clear pictures that can help improve diagnosis and treatment planning
- Offers high resolution imaging of prostate, colon, rectum, cervix, and surrounding areas
- Provides visualization of prostate internal structures, including prostate capsule and neurovascular bundles
- Inflatable design conforms to prostate shape and size for immobilization and optimal placement
- Disposable RF Probe heads connect quickly to interface box to allow for rapid setup and disassembly



## 6-Channel Flex Array

- 6-channel, 6-element flex array
- High SNR and uniformity
- Can operate as a 6-channel coil with both coil halves or a 3-channel coil with one half
- Coil half dimensions (L x W x H): 30 cm (11.8 in) x 30 cm (11.8 in) x 2 cm (0.8 in)
- Maximum FOV: 24 cm (9.5 in)
- Cable length: 155 cm (61 in)!



# MR Enabled Therapy and Accessories

## Radiation Oncology Options

Combining the Discovery™ MR750w 3.0T with GEM advanced imaging capabilities with the Radiation Oncology Options offering helps minimize potential registration errors between MR and CT within radiation treatment plans, for improved confidence in tumor targeting and preservation of healthy tissue. Additionally, seamless integration with AdvantageSim MDtm simulation software and Integrated Registration on the GE AW workstation allows MR images to be easily incorporated into the Radiation Oncology workflow.

## MR Guided Focused Ultrasound

Your facility can offer a completely non-invasive treatment for uterine fibroids with the addition of an Exablate MR guided Focused Ultrasound therapy table to your MR system, which has been used in 6,500 procedures worldwide.

## Surgical Suite

The Surgical Suite offering is an effective solution for incorporating MR imaging into your surgery center. Through seamless integration with surgical navigation systems, surgeons can retrieve archived images and fuse them with newly acquired intra-operative MR images. This advanced technology can assist in real-time surgical decision making and improved tumor resections.





GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

express  
Workflow



# eXpress Workflow

eXpress Workflow is designed to standardize and accelerate workflow for patient setup, exam prescription, scanning and post processing. eXpress Workflow can begin before the patient enters the magnet room and exams can be completed within a few mouse clicks - delivering quality and consistency for all patients and from all technologists. At the same time, eXpress Workflow maintains the flexibility needed to rapidly adapt and optimize exams for patient specific situations.

## Exam Setup



## Modality Worklist

### Automated and standardized rapid setup

- Allows the MR protocol to be selected and linked to the patient record in advance of the patient's arrival
- For sites with full DICOM connectivity, select the patient from the Modality Worklist, start a new session and view the relevant exam details on the In-Room Operator Console
- Add critical patient information such as allergies, pre-medication, pregnancy status and history



## Protocol Tools

### Search, select and one click to share

- Protocol Libraries: GE Optimized (preloaded protocols) and Site Authored (customized and saved).
- Protocols can be saved based on patient demographics, anatomy, scan type, or identification number for rapid search
- Commonly used protocols can be flagged for quick selection from the Modality Worklist
- One-click to Share ProtoCopy - Enables a complete exam protocol to be shared with the click of a mouse and provides a process for managing protocols across multiple systems as well as saving protocols for back up
- Step-by-Step Protocol Notes - Guide the user through the entire clinical routine procedure via expert inputs
- Step-by-Step Protocol Notes - Can be edited by the user for site specific instructions
- Step-by-Step Video guides - Provide simplified video instructions on-console



# Patient Setup

## GEM eXpress Patient Table


### Safety, Comfort and Efficiency

- Reduce patient transfers – transfer outside the magnet room directly to the GEM Express table
- Accelerate emergency egress – can be undocked and removed by one user in under 30 seconds typically
- Automatic coil disconnect – in time sensitive situations the system coils are automatically disconnected
- Patient choice – feet-first or head-first positioning for all supported exams
- Reduce in-room patient setup and address privacy by fully preparing the patient and coils for an exam outside of the magnet room
- Integrate arm-boards and IV pole to support patient for transport
- Embedded GEM posterior array and multiple high density surface coil connectors
- IntelliTouch landmarking sensors
- Compatible second table, prepare the next patient outside the magnet room while scanning the current patient

GEM Express Patient Table	
Configuration	Detachable and mobile
Minimum & Maximum Height	70cm to 93cm continuous
Table Drive	Automated power-driven vertical Automated power-driven longitudinal
Longitudinal Speed	30 cm/sec (fast) and 0.5 cm/sec (slow)
Total Scannable Range	205 cm
Maximum Patient Weight for Scanning	227 kgs (500 lbs)
Maximum Patient Weight Detached and Mobile	227 kgs (500 lbs)
Maximum Lift Capacity	227 kgs (500 lbs)
Patient Transport Accessories	Self-storing non-ferrous IV pole Positioning Pads Immobilization Straps
Landmarking	Laser alignment with S/I and R/L alignment IntelliTouch touch sensors
Coil Connection Ports	2 high density, auto-sensing ports




# Patient Setup (continued)

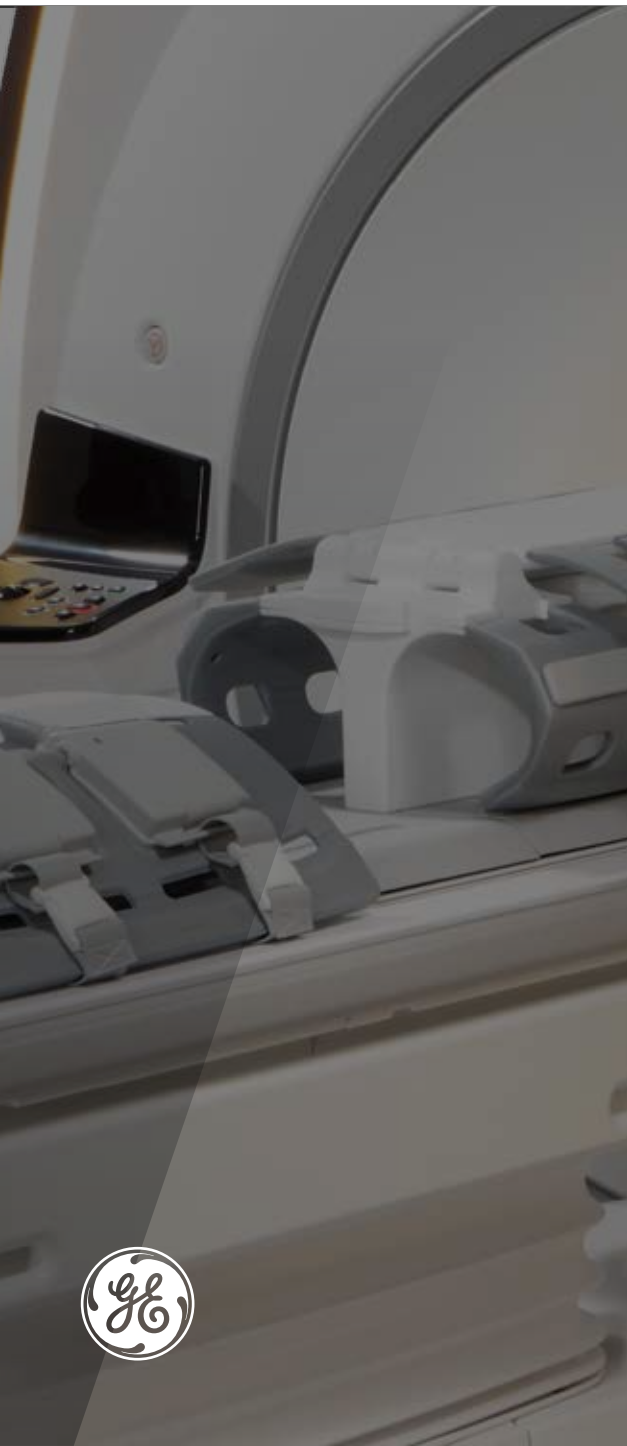

## GEM Coil Suite

### Multi-station coverage with one set-up

- Less coil lifting – GEM Posterior Array is embedded in the GEM Express Patient Table, eliminating the need to remove the coil
- Confirm coil connection – automatic information when coils are properly connected
- Head-first/Feet-First Scan – patients can be scanned feet-first or head-first for all supported exams

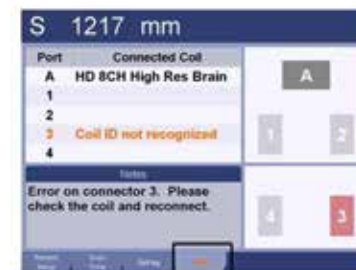

## IntelliTouch

### Touch to Landmark

- IntelliTouch sensors for simplified non-laser patient landmarking
- With IntelliTouch technology, the user can touch to complete
  - Patient landmarking
  - Localizing to the surface coil for auto-coil selection
  - Move patient to scan
  - Start scanning (with AutoStart activated)
  - Acquire, process and network images





# Patient Setup (continued)

## In-Room Operator Console and Control

### Full Control from table side

From the in-room operator console and controls, the user can:

- Position the table
- Return the table to home
- Stop the table movement
- Control multiple levels of in-bore ventilation and lighting
- Display of patient name, ID, study description
- Display and entry of patient weight
- Display and entry of patient orientation and patient position
- Cardiac waveform display and ECG/EKG lead confirmation
- Gating control for trigger select, invert and reset
- Respiratory waveform display
- IntelliTouch technology landmarking
- AutoStart to initiate scanning of the selected protocol
- Display connected coils and coil status
- Display of table location and scan time remaining
- Activate Screen Saver

The in-room display also allows for the integration of third-party visualization tools.





# Scanning

## Workflow Manager

### Linking and Auto Functions

AutoStart	Automatically initiates scanning of the selected protocol upon closure of the scan room door.
AutoCoil	Automatically determines the optimum coil elements to activate for scanning. If the prescribed field-of-view changes, AutoCoil automatically adjust the selection. The user has the option to review and edit the selection.
AutoScan	Automatically scans the prescribed series without user interaction. For series requiring a contrast injection, the Workflow Manager will pause and await user interaction.
Auto-calibration	For acquisitions that utilize ASSET parallel imaging or PURE surface coil intensity correction, Auto-Cal will prescribe and acquire a calibration scan based on the prescribed imaging volume.
AutoVoice	Delivers user selected, pre-recorded instructions to the patient at defined points in the acquisition to help ensure exam consistency. AutoVoice includes instructions in 14 languages and also allows the user to create and save unique instructions for specific local needs.
PB Navigators	Enable free-breathing body imaging for patients unable to breath-hold. The diaphragm tracker pulse automatically places and updates to streamline workflow and eliminate the set-up time associated with respiratory bellows. Auto Navigators can be used with a broad range of imaging techniques including dynamic contrast enhanced T1-weighted imaging.
ReadyBrain	Automates localizer acquisition, scan plane prescription, scanning, and post-processing for brain exams. ReadyBrain automatically calculates the mid-sagittal plane and determines the AC-PC line/OM line for 2D/3D prescription as well as corrects for extreme (>45 degree) rotation.
QuickSTEP	Automatically prescribes, acquires, and combines images from multiple stations. QuickSTEP acquires mask datasets and then secondary datasets from multiple stations (same locations), and automatically subtracts the mask datasets from the secondary datasets to create one subtracted series.
eXpress Prescan 2.0	Reduces pre-scan time for FSE-based techniques by up to 40% with a new calibration algorithm that reduces pre-scan time and consequently overall exam time.
Pause and Resume	Allows the user to pause a scan in progress, to respond to a patient need, and then resume mid-scan (without repeating scan).



# In-line Processing & In-line Viewing

## In-line Processing

### Automated post-processing

- Automated post-processing of specific applications
- Automatic opening and loading to advanced visualization tools when appropriate
- Automated in-line processing can be stored within the protocol

3D ASL series	Automatic compute and save
Diffusion weighted series	Automatic compute and save
Diffusion tensor series	Automatic compute and save
eDWI series	Automatic compute and save
Image filtering: A-E, SCIC, reFINE	Automatic compute and save
Maximum/Minimum Intensity Projection	Automatic compute and save
Reformat to orthogonal plane	Automatic compute and save
T2 map for cartilage evaluation	Automatic compute and save
3D Volume Viewer	Automatic load
BrainStat	Automatic load
FiberTrak	Automatic load
Image Fusion	Automatic load
Interactive Vascular Imaging	Automatic load
Pasting	Automatic load






# In-line Processing & In-line Viewing (continued)

## In-line Viewing

### Enhanced Visualization

In-line viewing allows the user to seamlessly and conveniently view, compare, and analyze images (during scan progress). The user simply selects the series, or multiple series, to view from the workflow manager, and the images are displayed along with the image display tools.

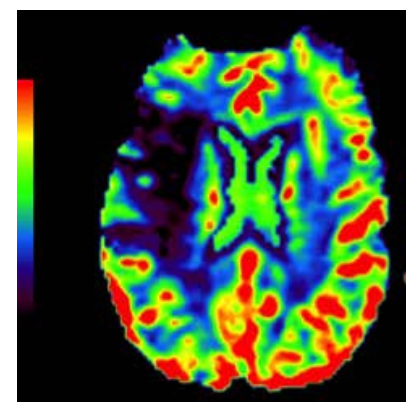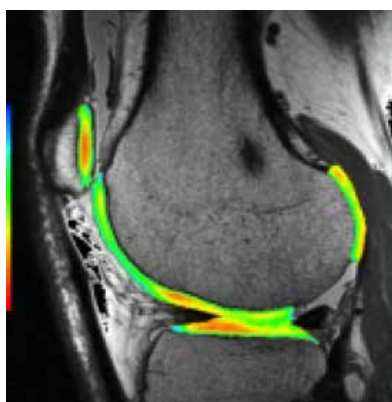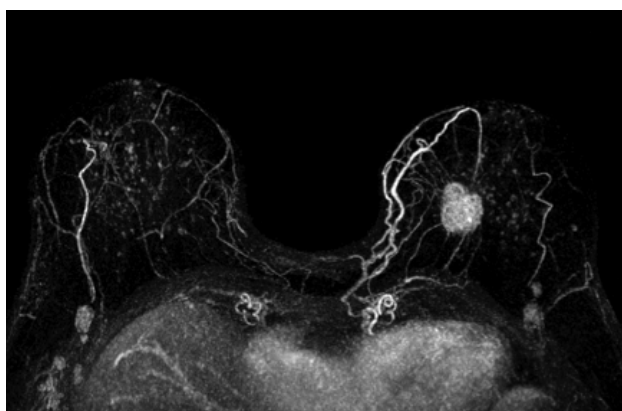

GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

Visualization





# READYView on MR Operator Console


## Benefits

- 3D ROI
- 3D Reformat MPR
- Auto-contour
- Distortion Correction
- Fusion & Registration
- MIP & HD MIP
- Motion Correction
- Multiparametric protocols
- Multiple graphics display
- Ratio AB/CD
- Reformat & Graphview
- Subtraction
- Volume Rendering
- Volume segmentation ROI

## Integrated Post Processing & Advanced Visualization

READYView is an image analysis software that allows the user to process dynamic or functional volumetric data and to generate maps that display changes in image intensity over time, echo time, b-value (Diffusion imaging), frequency (Spectroscopy). The combination of acquired images, reconstructed images, calculated parametric images, tissue segmentation, annotations and measurement performed by the clinician allows multi-parametric analysis and may provide clinically relevant information for diagnosis.

- Automatically selects the most relevant post processing protocol\*<sup>1</sup>
- Provides guided workflow and general assistance for the processing algorithms
- Multiparametric protocols selection for Brain, Breast, Liver, Knee and Pelvis studies when two or more functional series are present
- MR general review enables efficient reading of multi-contrast exams based on Smart Layout Technology
- One-click – to select and process functional data
- One-click – to save all generated parametric images
- One-click – to save and restore the state of processed images at any stage
- One ROI – display all multiparametric images and get all related functional values from a single ROI
- Export – display and export ROI statistics from the summary table
- Export graph values as csv file
- Customize workflows with adjustable layouts, personalized parameter settings, and custom review steps

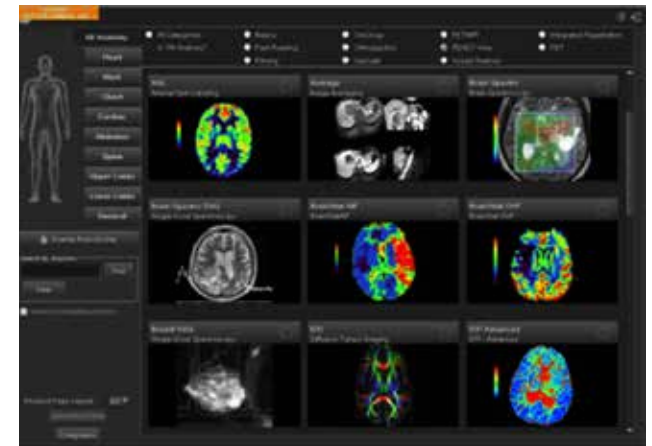

<sup>1</sup>When only one protocol is compatible with the selected data, the access is made through the One-Touch mode. If more than one protocol is compatible, the Protocol page opens for user selection.





## Standard Protocols

### READYView One-Touch

Protocols uses display intelligence with pulse sequence, image contrast and scan plane recognition to enable direct access between a unique post processing that is associated with the series selection.

### One-Touch ADC and eADC

Provide algorithms to process DWI1 images to generate ADC maps and eADC maps to eliminate T2 “shine through” in the isotropic (trace) DWI.

### One-Touch ASL

ASL READYView has algorithms that calculate Cerebral Blood Flow maps from a 3D ASL series. ASL acquisition is a non-invasive, one-click application that allows whole brain CBF measurements.

### One-Touch Brain & Prostate Spectroscopy


The READYView Brain and Prostate MR spectroscopy protocols are used to display functional maps for metabolites and metabolite ratios in the brain and prostate.

### One-Touch MR-Touch

READYView MR-Touch is a post process of an MR-Touch acquisition, which is a Phase Contrast (PC) application that generates an image contrast related to the shear stiffness of soft tissue. An algorithm is used to derive a relative stiffness map (Elastogram) and wave images from the phase images.

### One-Touch T2 MAP

The READYView T2 Map protocol post processes data sets acquired using the T2 Map (CartiGram) application. The T2 Map acquisition is displayed in READYView, where the T2 relaxation time color map is coded to capture T2 values from the TE range of the acquired images.

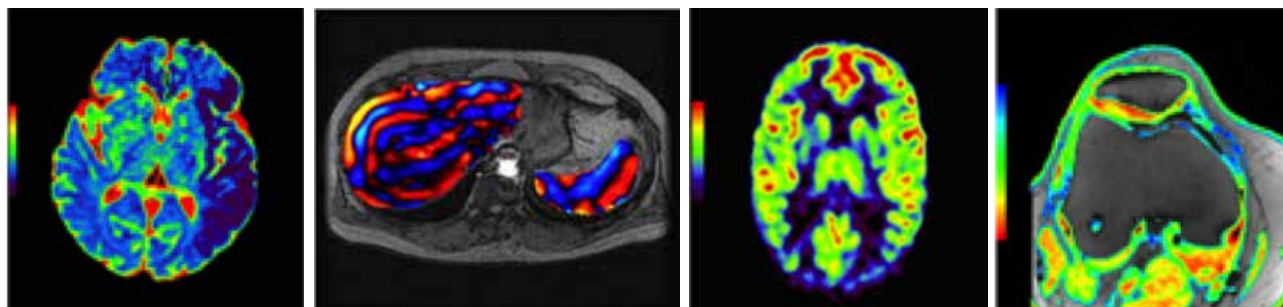

# READYView (continued)

## BrainStat

BrainStat is an MR Time Course imaging READYView protocol that provides accurate spatial resolution for brain tissue viability given by hemodynamic parameters: BV, BF, TTP, MTT (SVD), BAT, Tmax. These hemodynamic parameters can provide unique information on tissue changes and improve delineation of vascular-deficient or vascular-rich regions in normal and abnormal anatomy.

## MR Standard

MR Standard is a time course protocol. The READYView MR Standard is a time course protocol that can be used to create the following maps: enhancement integral (negative and positive), time to peak, mean time to enhance, maximum slope of increase, maximum slope of decrease.

## SER

SER is a time course protocol for analyzing T1-contrast changes. The READYView SER protocol can be used to create the following maps: Positive enhancement integral, signal enhancement ratio and maximum slope of increase.

## FiberTrak


FiberTrak is designed for the advanced analysis of MR images acquired with a DTI technique. It allows for processing of isotropic, ADC and FA maps among other options. The FiberTrak option augments this functionality to allow DTI processing to create: 2D color orientation maps, 2D color eigenvector maps and 3D tractography maps.

## fMRI

Functional imaging or BOLD provides fMRI analysis using the correlation coefficient algorithm to analyze an image set. Neuronal activity of either motor or cognitive functions can be mapped by fMRI through changes in signal intensity. The resulting functional maps can be used for mapping the motor cortex and higher cognitive regions of the brain.

## R2 Star

The R2 Star feature uses water proton transverse relaxation rates (R2) technique. It provides parametric maps for R2\* (Hz) and T2\* (ms). The R2\* values vary with tissue characteristics such as iron concentration.

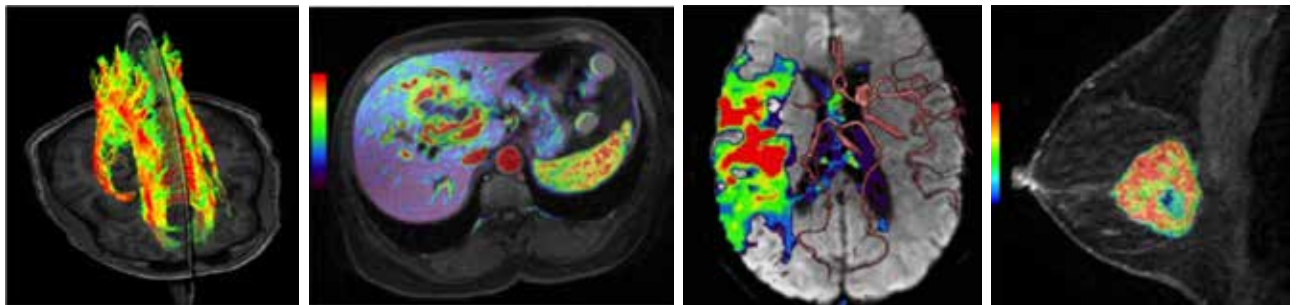


## Multiparametric Protocols: Visualization at a Glance


READYView multiparametric protocols provide a guided workflow to streamline post processing and analysis of multiparametric studies. All measurements can be obtained with one ROI and the user customizable workflow has the ability to display all processed maps in one screen.




**MR Brain\***  
Diffusion, Perfusion, Brain Spectroscopy,  
Brain SVQ DTI and ASL

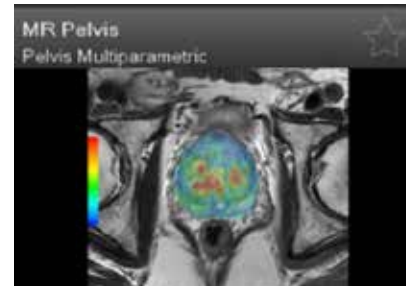

**MR Pelvis\***  
Diffusion, Perfusion, and  
Prostate Spectroscopy

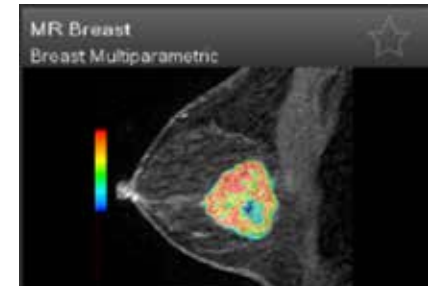

**MR Breast\***  
Diffusion, SER and Breast SVQ

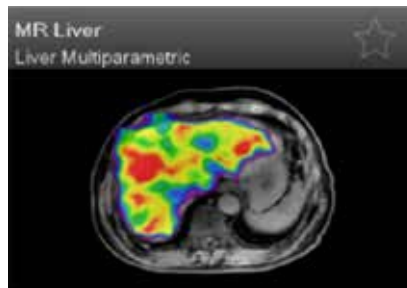

**MR Liver\***  
Diffusion, R2 Star or MR Touch



**MR Knee\***  
T2 Map

\*2 Requires two or more of the functional series selected.



GE Healthcare

# Discovery™ MR750w 3.0T Wide Bore MR with GEM Suite

Siting





# Siting and Other Specifications



## Typical Room Layouts

	System configuration minimum values
Magnet Room W x D	21.2 sq.m
Minimum Ceiling Height	2.5 m (8 ft. 2.4 in) min ceiling height
Equipment Room W x D	7.9 sq.m
Control Room W x D	3.2 sq.m
Minimum Total Area (M <sup>2</sup> )	


## Fringe Field

	Axial	Radial
0.5 mT (5 Gauss)	5.2 m	2.8 m
0.1 mT (1 Gauss)	7.8 m	4.8 m





# Siting and Other Specifications (continued)

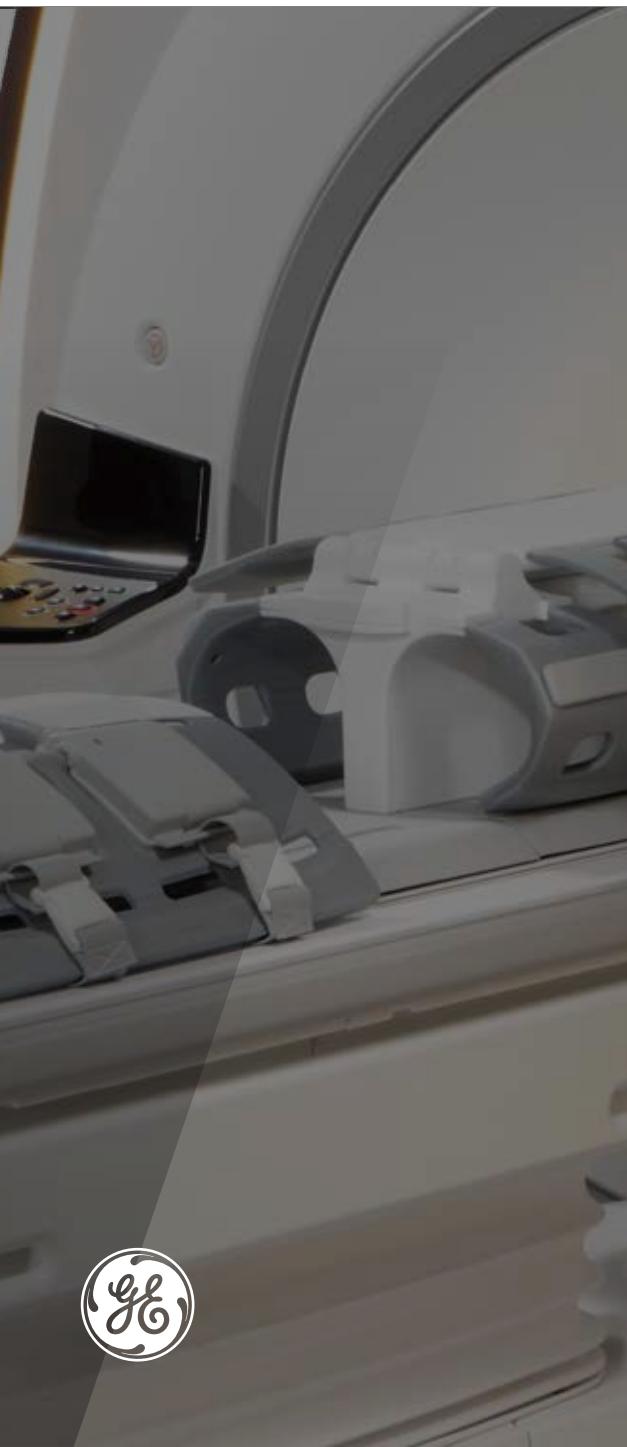

## Electrical Supply Requirements

Supply system recommended configuration:

- 3-phase grounded WYE with neutral and ground (5-wire system)

Note: Neutral must be terminated inside main disconnect control

Alternate configuration:

- 3-phase DELTA with ground (4-wire)
- Recommended grounded delta configuration
- Voltage: 480/415/400/380/Vrms

## Power Consumption

Power consumption depends on actual usage. They exclude consumption by the shield cooler compressor (9 kVA). The following values are approximate:

Standby (no scan)	<17 kVA
Typical Power†	41 KVA
Maximum continuous sustained power (> 5 secs)	99 kVA
Peak instantaneous power (< 5 secs)	123 kVA

## RF Shielding


100 db. for 10 - 100 MHz plane wave

## Workspace Monitor Positions

	Maximum field strength
LCD flat panel monitor	5 mT (50 Gauss)



# Siting and Other Specifications (continued)

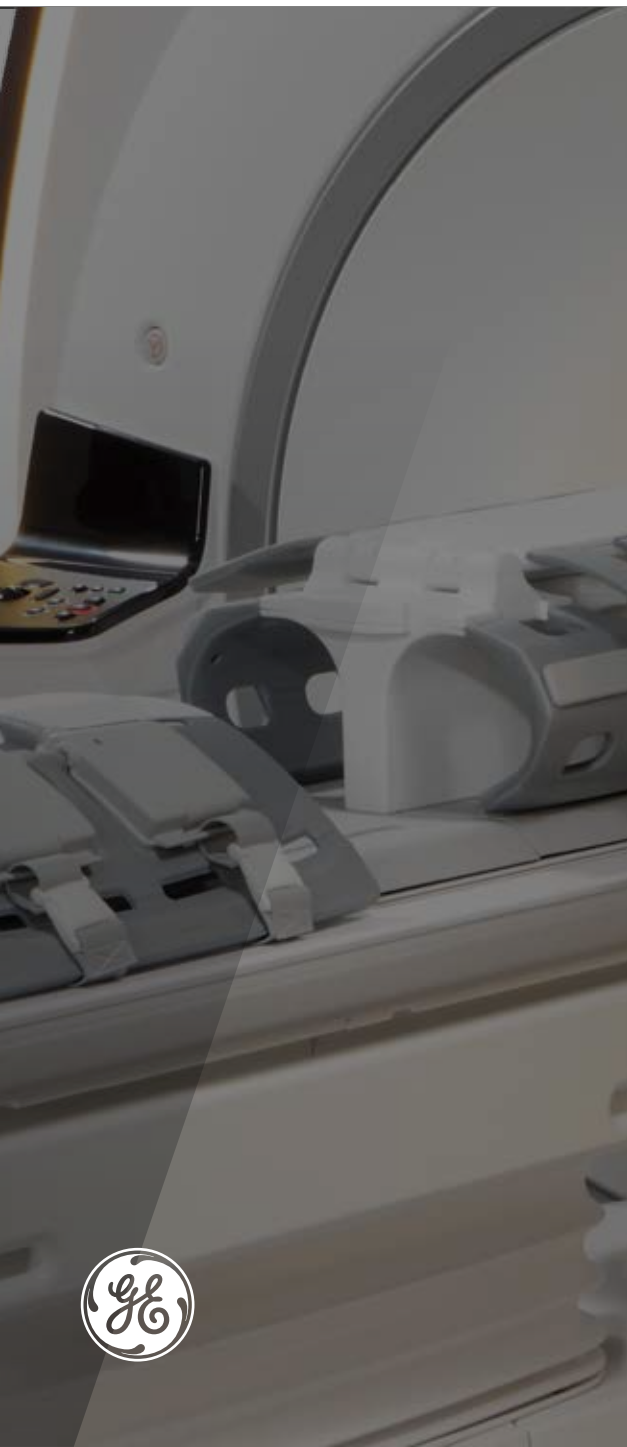

## Temperature and Humidity Requirements

	Magnet Room	Control Room	Equipment Room
Temperature	15 - 21° C	15 - 32° C	15 - 32° C
Max. temperature change rate	3° C/hour	3° C/hour	3° C/hour
Humidity (non-condensing)	30 - 60%	30 - 70%	30 - 70%
Max. humidity change rate	5% RH/hr	5% RH/hr	5% RH/hr

## Altitude Requirements

Upper limit	2600 m
Lower limit	-30 m



# Miscellaneous

## Alternative environments

Modular buildings may also be available (including air-conditioning, heating, chiller, RF shielding, additional magnetic shielding in walls). Contact your local GE representative for GE certified designs and vendors.

Please ask your local GE project manager for a comprehensive installation and siting manual.

## Filming considerations

Filming requires the Discovery™ MR750w 3.0T analog or digital filming.

Interface (purchased separately) unless DICOM print will be used exclusively for software filming to DICOM print peripheral devices. An Analog/VDB or Digital/LCAM camera interface is typically required for most installations.

## Accessory Package

- SPT phantom set with storage cart
- Customer diagnostic software
- Operator manuals
- Patient log books

## Emergency stop

Disconnects electrical power from RF and gradient components in the magnet room (duplicate control at the magnet).

## Warranty

The published GE warranty in effect on the date of shipment shall apply.

## InSite\* Remote Diagnostics

GE's unique remote service and applications support including magnet monitoring. Also allows downloading of applications software such as eFlexTrials program.

## Optional capabilities

Some features and capabilities listed in this data sheet are optional with a Discovery™ MR750w 3.0T and are subject to change without notice. Contact a GE representative for the most recent data.

## GE regulatory compliance

The Discovery™ MR750w 3.0T complies with all applicable safety standards, including but not limited to IEC60601-1 or ES60601-1 (Electromagnetic Compatibility).

Laser alignment devices contained within this system are appropriately labeled according to the requirements of the FDA's Center for Devices and Radiological Health (CDRH) and I EC 60825-1.





## Imagination at work

Product may not be available in all countries and regions.  
Contact a GE Healthcare Representative for more information.

GE Healthcare  
3200 N. Grandview Blvd.  
Waukesha, WI 53188  
USA

© 2016 General Electric Company – All rights reserved.

General Electric Company reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation.  
Contact your GE Representative for the most current information.

GE and GE Monogram are trademarks of General Electric Company.

GE Healthcare, a division of General Electric Company. \*Trademark of General Electric Company.  
Intel® and Xeon® are trademarks of Intel Corporation. NVidia® and Quadro® are trademarks of NVIDIA Corporation.