



Translate
7T research power
into clinical care
MAGNETOM Terra

MAGNETOM Terra – Translate 7T research power into clinical care

MAGNETOM Terra is designed to let you explore new territories in MRI by enabling powerful 7T research and enhancing clinical care. Uncover a whole new world of clinical insights with double SNR for more precision. Our advanced ultra-high-field (UHF) technology will keep you at the cutting edge of MRI, to attract the brightest minds to your facility, sharpen your competitive edge and strengthen your reputation. It delivers a fertile platform for unlocking research capabilities, publishing new insights first, and setting the pace in diagnostic imaging. Welcome to an exclusive research community. Welcome to a whole new world in MRI.

Welcome to clinical 7T.

Content

MAGNETOM Terra at a glance	04
Uncover a whole new world of clinical insights	06
with double SNR for more precision	08
Change the game in UHF business	26
with Siemens' 50% lighter 7T magnet	28
Unlock research beyond clinical limits	32
with 8-channel parallel transmit	34
Join the largest research community	50
with over 70% of all UHF users	52
Service and exchange	56
Technical specifications	60
Siemens Healthineers	62

MAGNETOM Terra – Translate 7T research power into clinical care

Unique Dual Mode functionality

- Uncover a whole new world of clinical knowledge with the flexibility to get more from your scanner
- Secure switch between research and clinical operation in less than 7 minutes
- Operating with 2 separate databases for clear research and clinical distinction

50% lighter 7T magnet technology

- First 7T scanner released for clinical use
- Lower weight and cold-shipment for easier integration in clinical environments
- Reduced operating costs thanks to Zero Helium boil-off

80/200 gradients, 8-channel pTX and up to 64 channels

- More power for greater diffusion MRI and functional MRI with 80/200 gradients
- Higher homogeneity for challenging body regions with 8-channel parallel transmit (in research mode)
- Higher acceleration factors with 64 receive channels

Double SNR for more precision

- 0.2 mm in-plane resolution to visualize previously unseen structures
- 0.14 cm³ voxel sizes for metabolic brain mapping¹ (in research mode)
- Submillimeter BOLD fMRI precision to visualize sub-cortical activations



> 70% global market share in 7T technology

- Over 70% of 7T and 100% of vendor-integrated > 7T MRI human scanners worldwide from Siemens
- 8 of 11 leading U.S. hospitals with a 7T, (2016–2017), trust Siemens when they decide for 7T investment³
- 74% of ISMRM UHF abstracts in 2017 were based on data from Siemens UHF systems³

The world's largest UHF community

- Largest installed base for exchanging ideas in a strong collaborative network
- An opportunity to enhance your reputation and competitiveness
- Incentive for the brightest minds in the MRI community to work with you



Uncover a whole new world of clinical insights

Discovering new ground in MRI can help you significantly enhance patient outcomes. Imaging at 7T offers more than double the SNR of 3T. This delivers potential for better lesion conspicuity, faster image acquisition to reduce motion artefacts, and earlier disease detection at submillimeter resolution. MAGNETOM Terra is the first 7T scanner released for clinical use. Its Dual Mode lets you switch between clinical and research tasks, unlocking new opportunities and providing a solid, well-founded platform for innovative results.

“Based on higher resolution, 7T provides new insights into gray and white matter disease in the brain, such as multiple sclerosis, focal cortical dysplasia, and hippocampal sclerosis. Furthermore, functional MR benefits from 7T based on a clinically relevant increase in functional sensitivity and specificity. In musculoskeletal imaging, 7T enhances the visualization of small joint structures and subtle pathologies, such as small meniscal tears, triangular fibrocartilage lesions, and early stages of cartilage degeneration.”⁴

Professor Siegfried Trattnig
Director of the MR Centre of Excellence,
MedUni Wien, Vienna, Austria

Uncover a whole new world of clinical insights – Double SNR for more precision with clinical applications in Dual Mode

Dual Mode flexibility

MAGNETOM Terra⁷ is the first 7T scanner released for clinical use. With selected neurological and musculoskeletal scan protocols, it has potential to uncover a whole new world of clinical care. Its unique Dual Mode functionality lets you switch between research and clinical operation, giving you flexibility to get more from your scanner.

Ultra-fine anatomical resolution

In brain and musculoskeletal MRI, 7T reveals details previously unseen at lower field strengths. For example, cerebral cortex imaging at 0.2 mm in-plane resolution may detect changes in cortical structure indicating early dementia. It also helps visualize

cortical microinfarcts and plaques in MS patients and delivers insight into the plaque-vessel relationship and iron accumulation.

Submillimeter fMRI

The BOLD contrast increases linearly with field strength. In clinical use, this could mean higher precision in oncology compared to 3T applications, for example, through smaller voxel sizes. Potentially, this can increase the accuracy of neurological pre-surgical evaluation of eloquent areas before tumor removal, while keeping scanning times viable².

Powerful image reconstruction

MAGNETOM Terra delivers improvements in workflow for easier operation and better patient handling. Leveraging the latest syngo MR E11 software platform, it lets you work in the same way as you do with cutting-edge 3T technology. What's more, it comes with the most powerful MaRS (Measurement and Reconstruction System) computing technology ever built.⁶



8 of 11 leading U.S. hospitals with a 7T, (2016–2017), trust Siemens when they decide for 7T investment³



Switch between
research and
clinical tasks with
Dual Mode

Ultra-fine
resolution
to visualize
details previously
unseen

Submillimeter
BOLD
fMRI precision for
pre-surgical
evaluation

Ultra-fast image
reconstruction⁶ and
syngo MR E11


Dual Mode offers the flexibility to switch from research to clinical tasks


Clinical Mode

- 1 transmit channel
- 11 kW RF power
- 2 coils (Head 32, Knee 28)
- Neuro and MSK optimized clinical applications


Research Mode

- single channel and 8-channel parallel transmit
- 8 x 2 kW RF power
- Wider range of RF coils
- Whole-Body WIP Applications

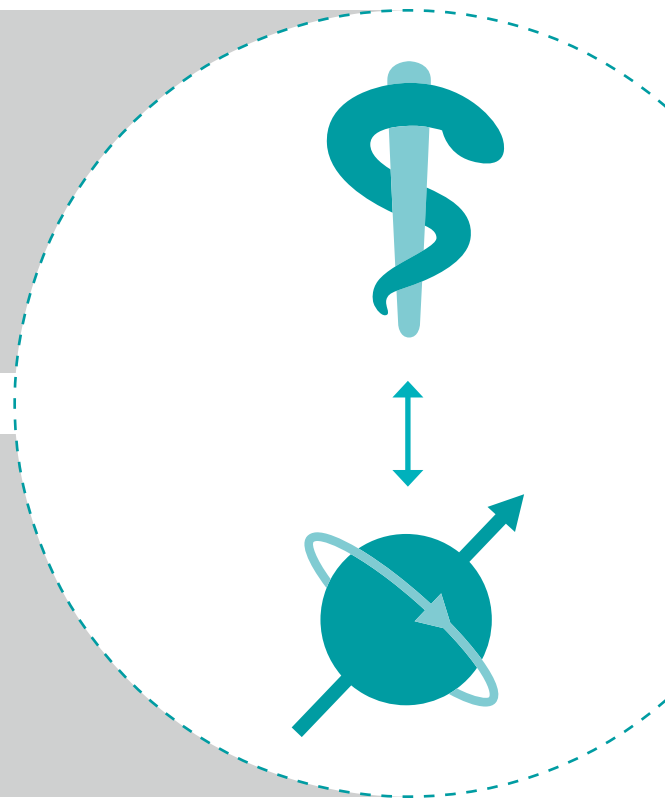

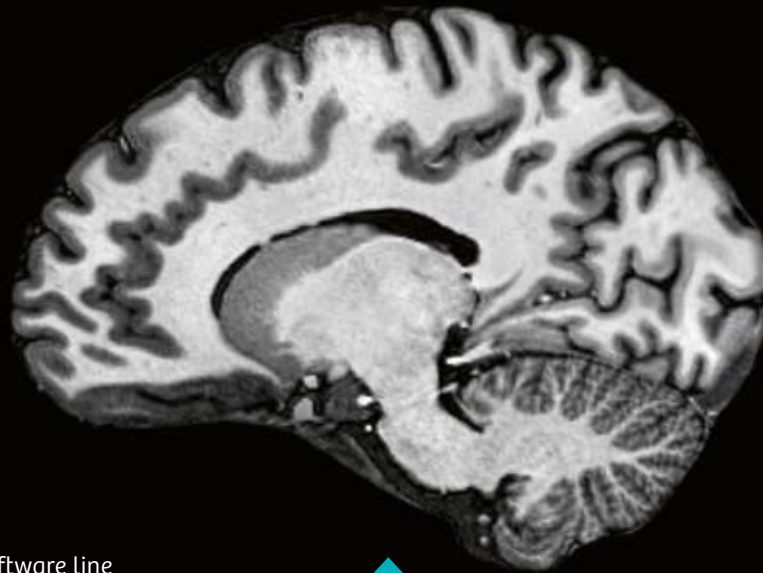

- syngo MR E11 software line
- XR Gradients 80/200
- Up to 64 receive channels²
- Latest MaRS computer
- 3rd order shims

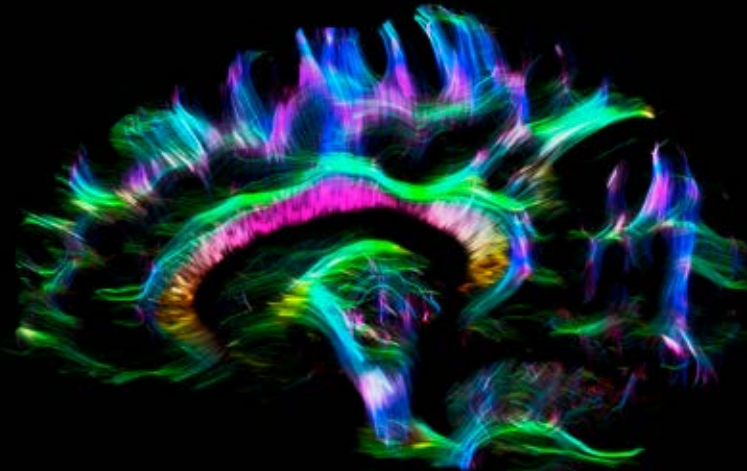


Clinical Mode – Multiple Sclerosis

Hyperintense MS lesion with hypointense center

DarkFluid TSE



0.3 x 0.3 x 3 mm³,
5:59 min

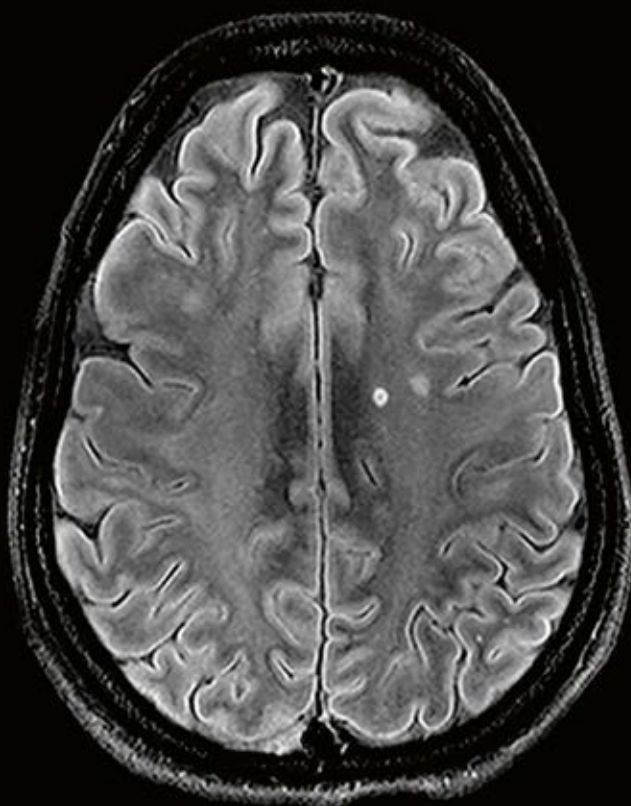

SWI minIP/phase

0.2 x 0.2 x 1.2 mm³,
5:38 min

Typical central vein and perivenular demyelination is visible.

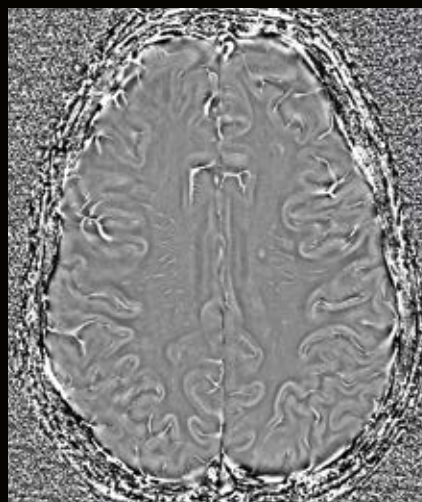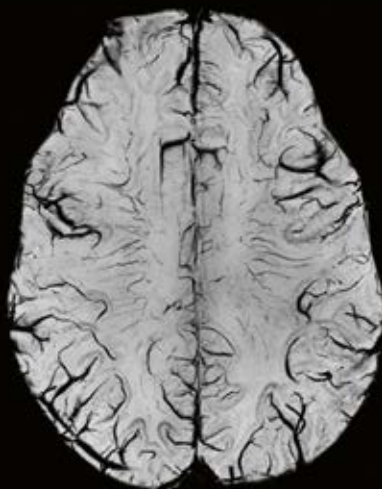

Multiple Sclerosis
with low lesion load



Multiple Sclerosis
with high lesion load



Dark Fluid TSE

0.3 x 0.3 x 3 mm³, 5:20 min

T2 TSE

0.2 x 0.2 x 3 mm³, 5:33 min

SWI

0.2 x 0.2 x 1.2 mm³,
5:38 min



Clinical Mode – Tumor


3D SWI of Glioblastoma

3D SWI minIP provides superior assessment of the microvasculature.

Erwin L. Hahn Institute for MRI, Essen, Germany

3 Tesla

0.85 x 0.72 x 2 mm³

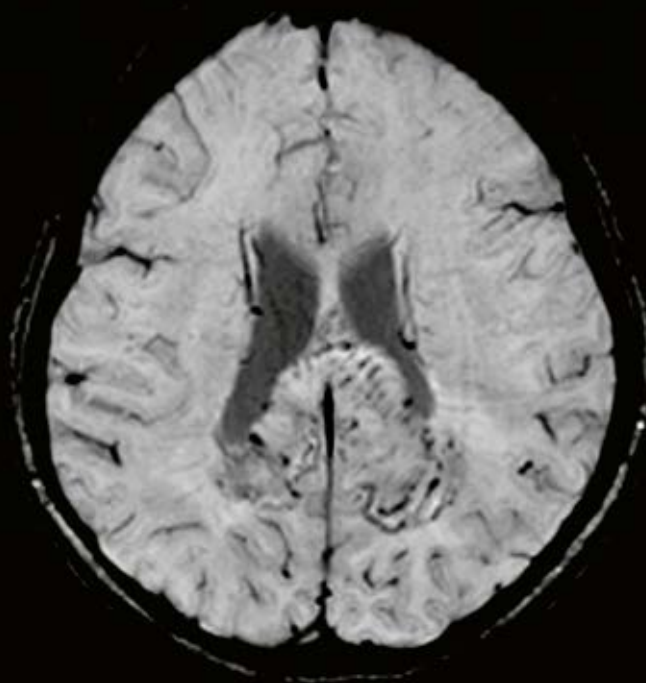

7 Tesla


0.25 x 0.25 x 1 mm³




Glioblastoma

Higher SNR for ultra-high 0.2 mm in-plane resolution for imaging tumor vascularization.

DKFZ, Heidelberg, Germany

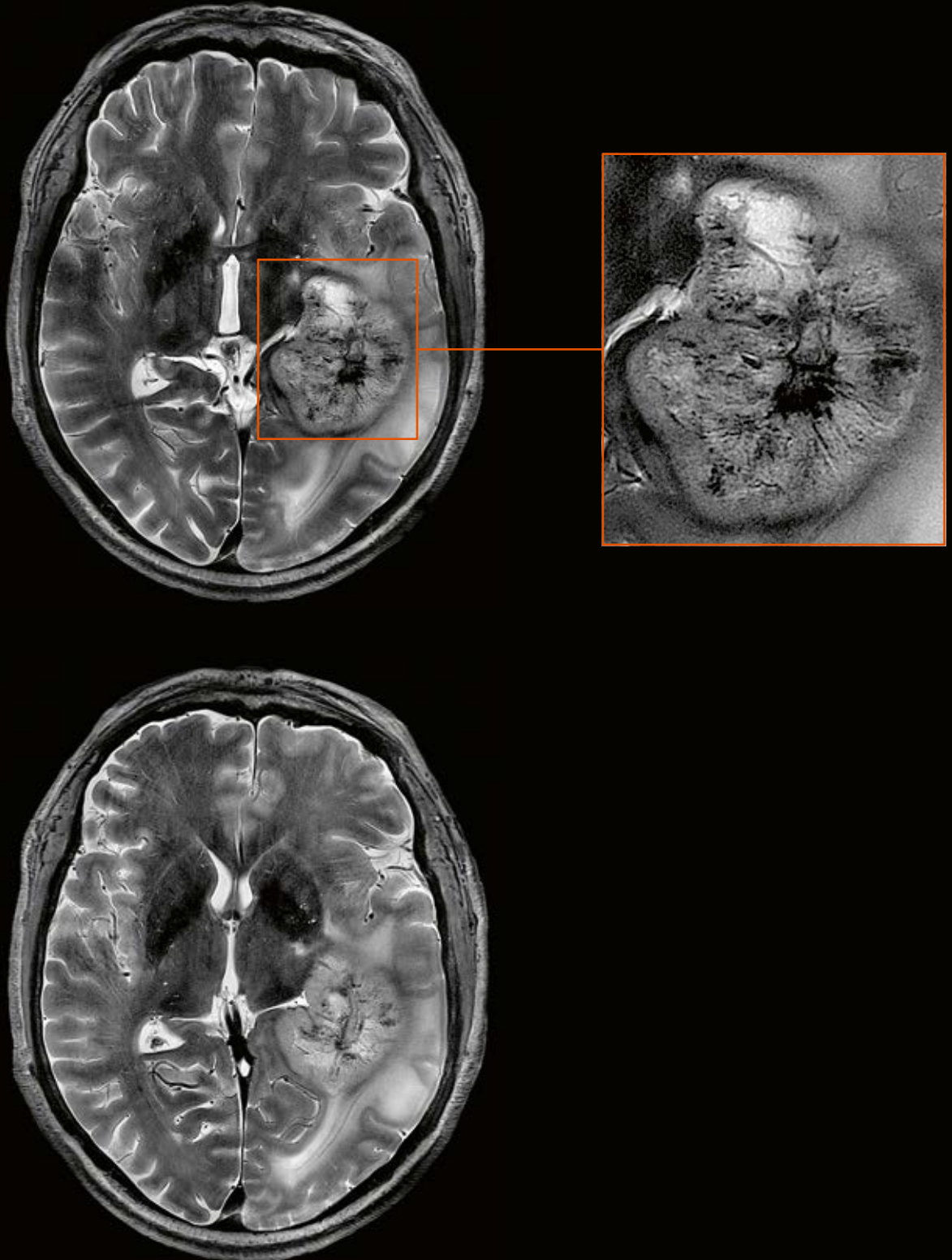

Clinical Mode – Stroke

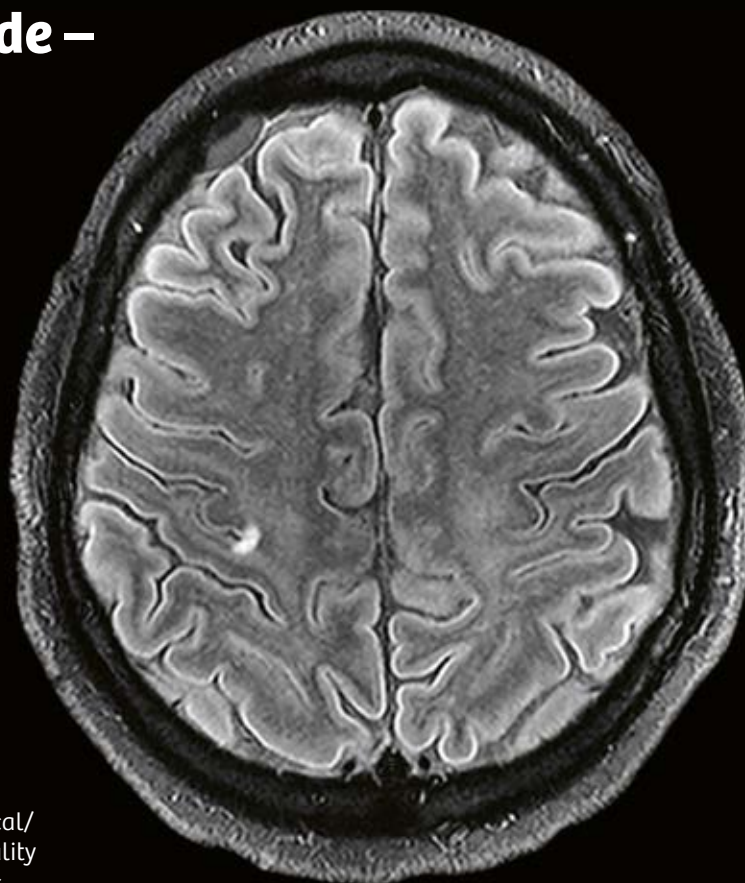


DarkFluid TSE

0.3 x 0.3 x 3 mm³,
5:59 min

Small focus of strong cortical/
subcortical signal abnormality
in the right precentral gyrus.

FAU, Erlangen, Germany

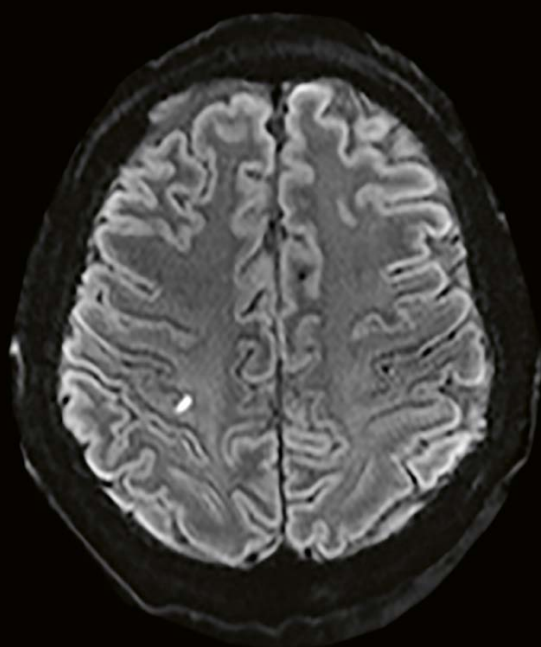


RESOLVE 1 x 1 x 3 mm³, 1:46 min



PD FS TSE 0.2 x 0.2 x 3 mm³, 5:14 min

Clinical Mode – Tumor

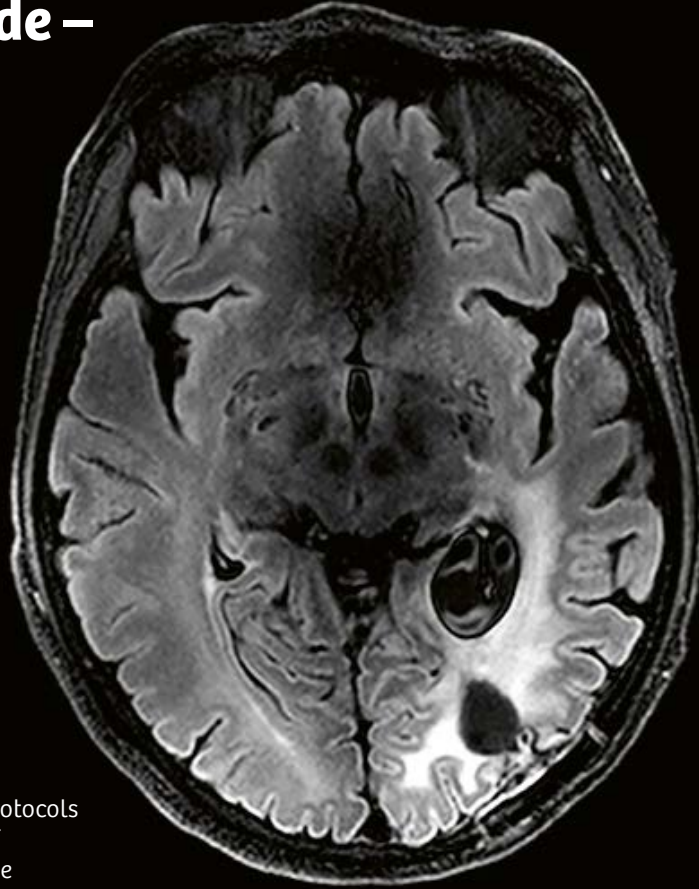

DarkFluid TSE

0.4 x 0.4 x 3 mm³,
4:22 min

High resolution standard protocols
for detailed visualisation of
pathologies, increased tissue
contrast and high resolution at 7T.

FAU, Erlangen, Germany




SWI minIP 0.2 x 0.2 x 3 mm³, 5:38 min



T2 TSE 0.2 x 0.2 x 3 mm³, 5:33 min

Clinical Mode – Healthy volunteer

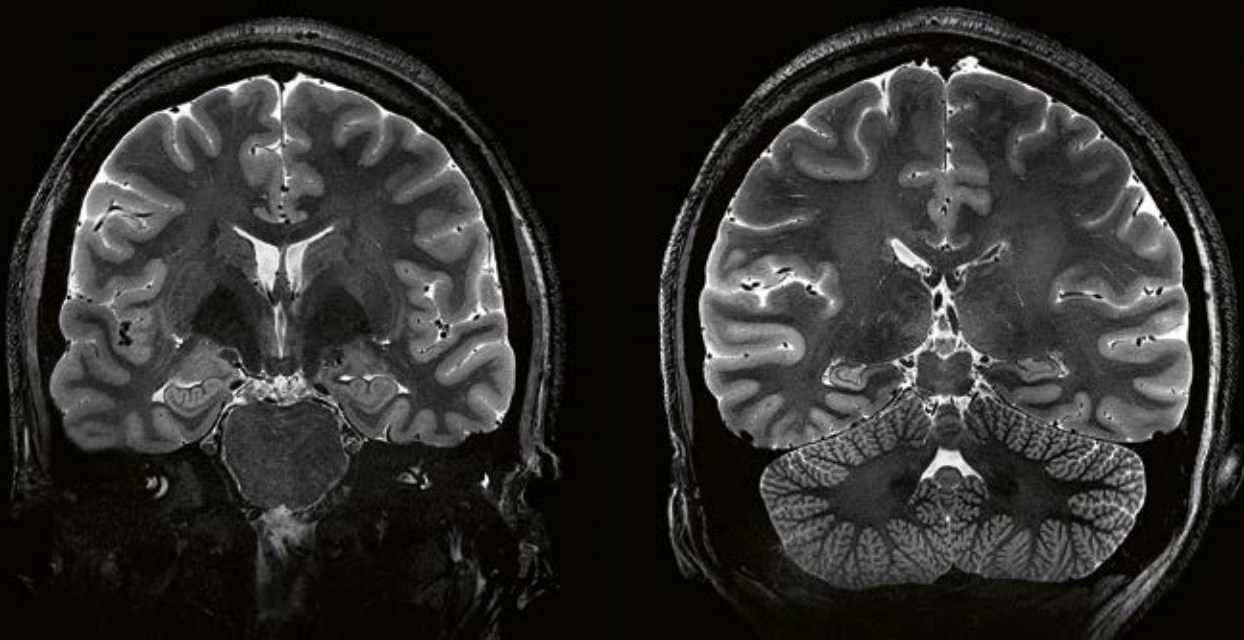

Hippocampus imaging


High-resolution imaging of the hippocampus at 0.25 mm in-plane resolution.

Scannexus, Maastricht, Netherlands


T2* weighted imaging

High-resolution imaging of the brainstem at 0.3 mm in-plane resolution.


MGH, Boston, USA

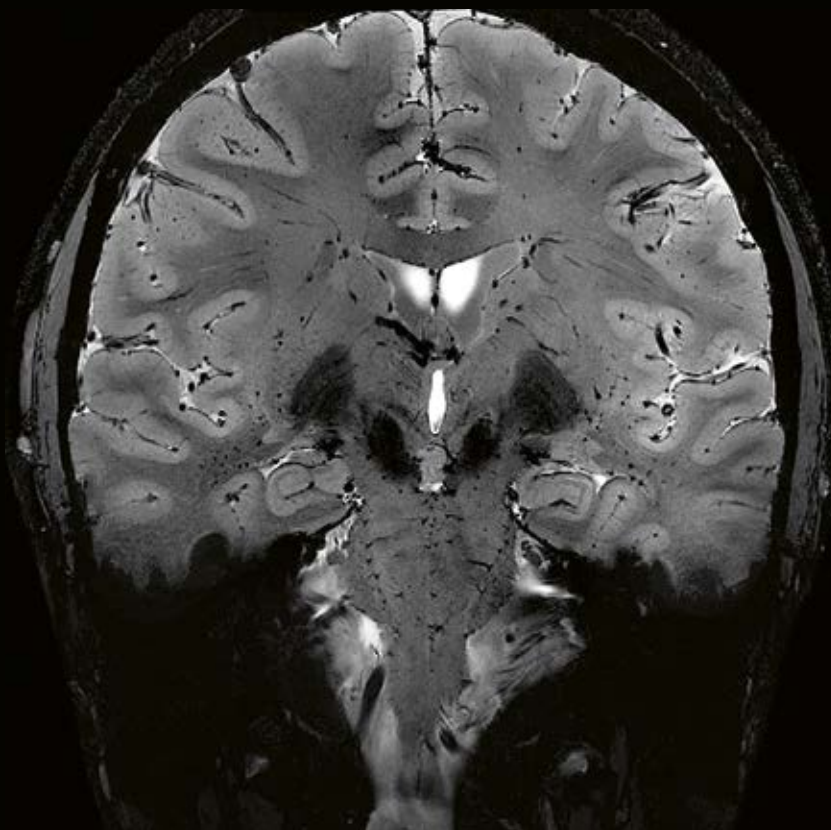

Clinical Mode – Diffuse axonal injury

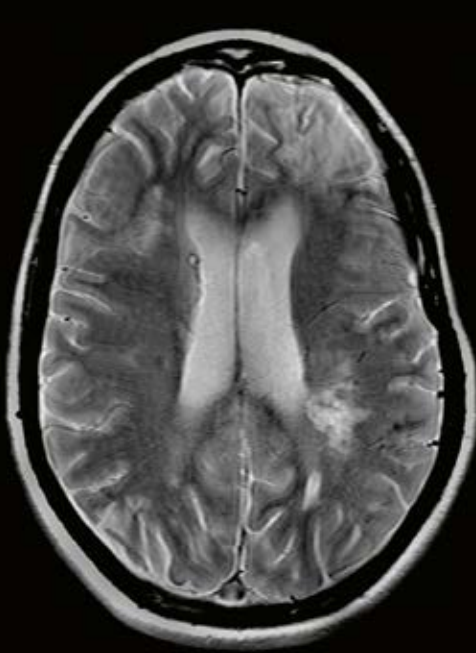

3 Tesla PD TSE, $0.7 \times 0.4 \times 5 \text{ mm}^3$, TA 2:38 min

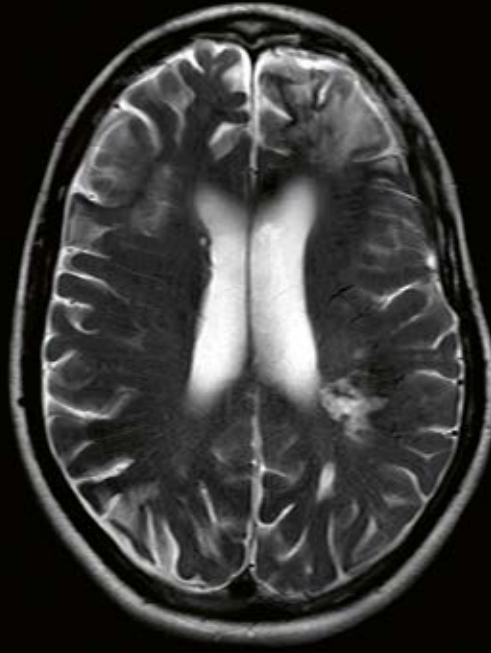

3 Tesla T2 TSE, $0.7 \times 0.4 \times 5 \text{ mm}^3$, TA 2:38 min



7 Tesla PD PD TSE, $0.2 \times 0.5 \times 3 \text{ mm}^3$, TA 3:09 min



7 Tesla T2 TSE, $0.2 \times 0.5 \times 3 \text{ mm}^3$, TA 3:09 min

Diffuse axonal injury

The higher sensitivity at 7T reveals hemosiderin from traumatic brain injury in PD images.

Erwin L. Hahn Institute for MRI, Essen, Germany


Clinical Mode – Healthy volunteer



Clear identification of anatomical structures with increased tissue contrast and high resolution at 7T.

0.6 x 0.6 x 0.6 mm³, 13:45 min


FAU, Erlangen, Germany

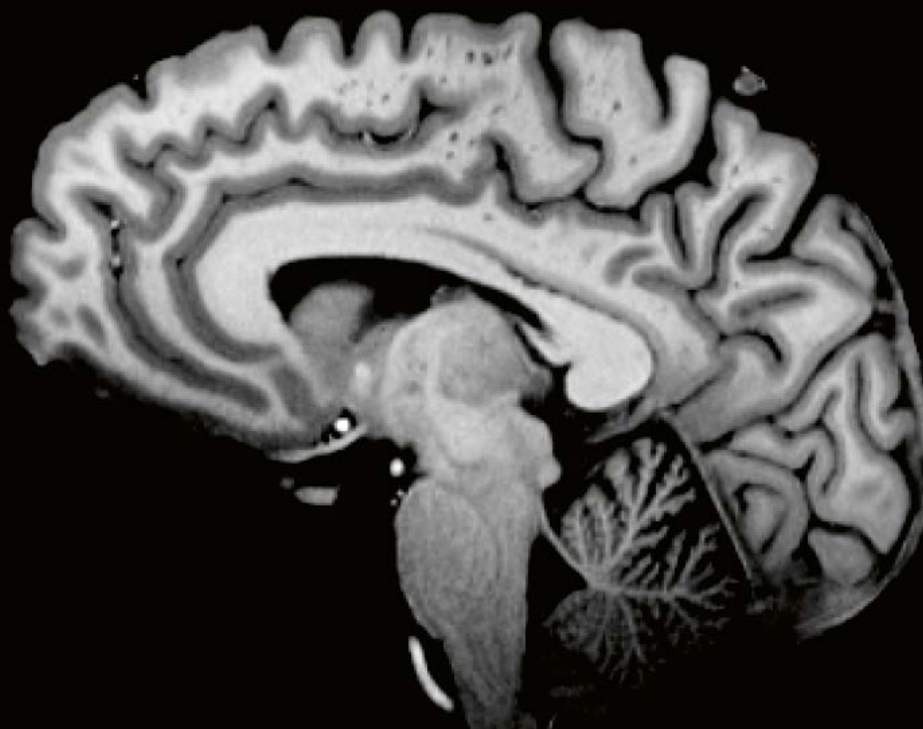

Clinical Mode – Stroke

Time of Flight (ToF) with 400 micron isotropic resolution reveals smallest vessels in the brain. The higher the signal and the longer the T1 at 7T are, the higher the quality of the Maximum Intensity Projection (MIP) gets.

Coronal

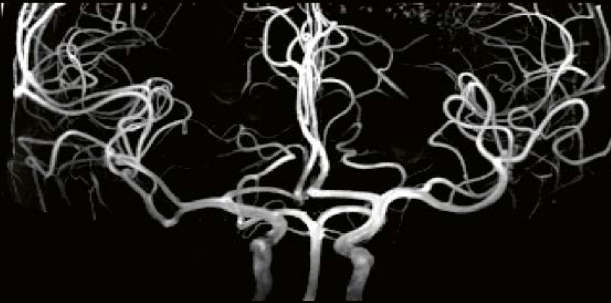

Sagittal




Axial

Visualize smallest vessels with
0.4 mm isotropic resolution.

0.4 x 0.4 x 0.4 mm³, 8:09 min

FAU, Erlangen, Germany

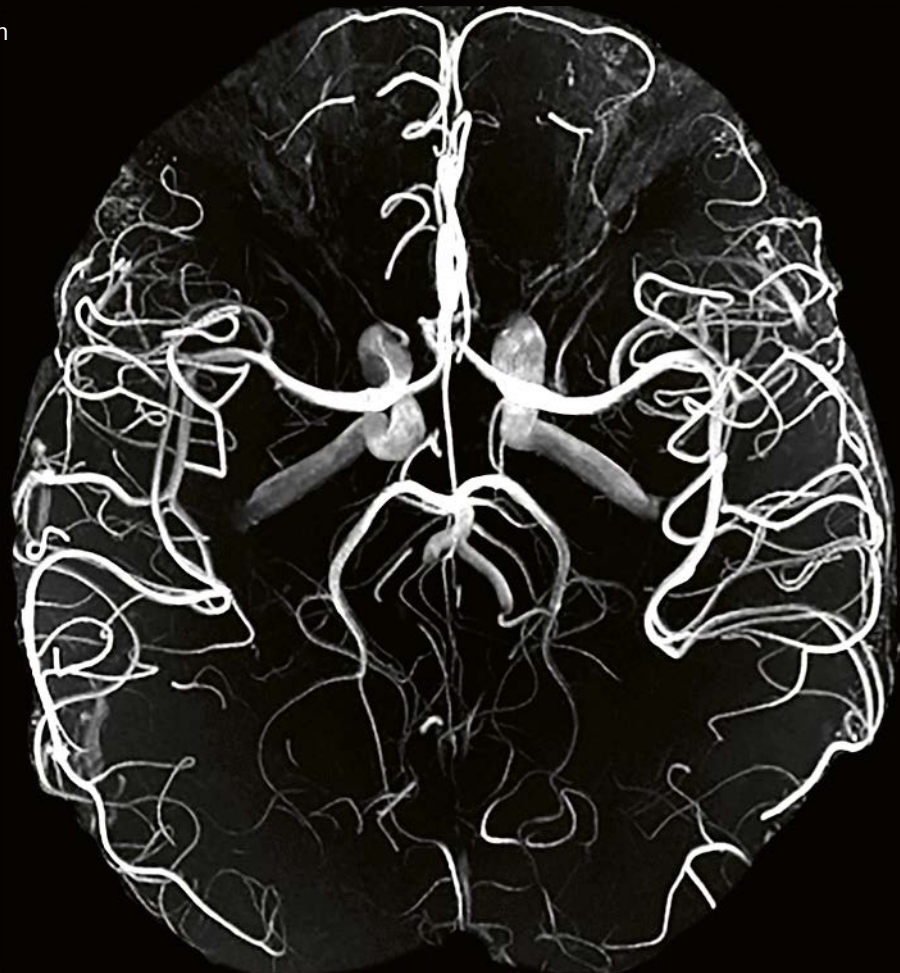


Clinical Mode – Healthy volunteer

0.2 mm in plane resolution

The basal ganglia Caudate, Putamen and Globus Pallidus can be differentiated.

Enlarged sections: cortical veins can be depicted.

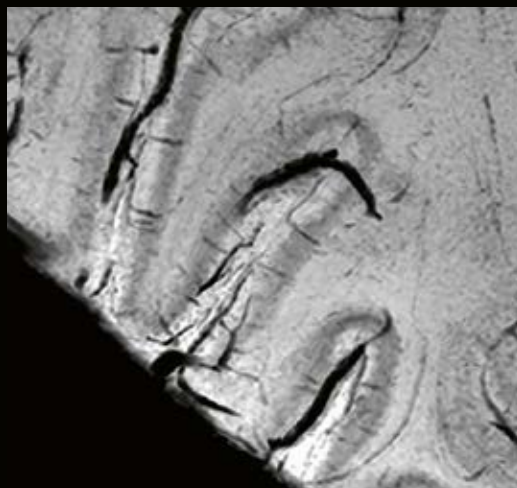

SWI 0.2 x 0.2 x 1 mm³, 10:59 min

Clinical Mode – Enchondroma


Fine structure visible in the lesion with different contrasts.



PD TSE FS $0.2 \times 0.2 \times 2.5 \text{ mm}^3$, 3:15 min



T2 TSE $0.3 \times 0.3 \times 2 \text{ mm}^3$, 3:24 min

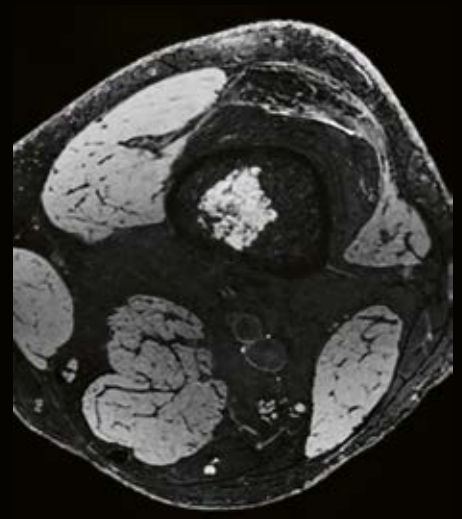

3D DESS $0.5 \times 0.5 \times 0.5 \text{ mm}^3$, 3:43 min

Clinical Mode – Healthy volunteer

Clear delimitation of anatomical structures, such as ligaments, vessels or cartilage.


T1 SE

0.2 x 0.2 x 2.5 mm³,
4:05 min



T1 qSE

0.3 x 0.3 x 2.5 mm³,
7:21 min

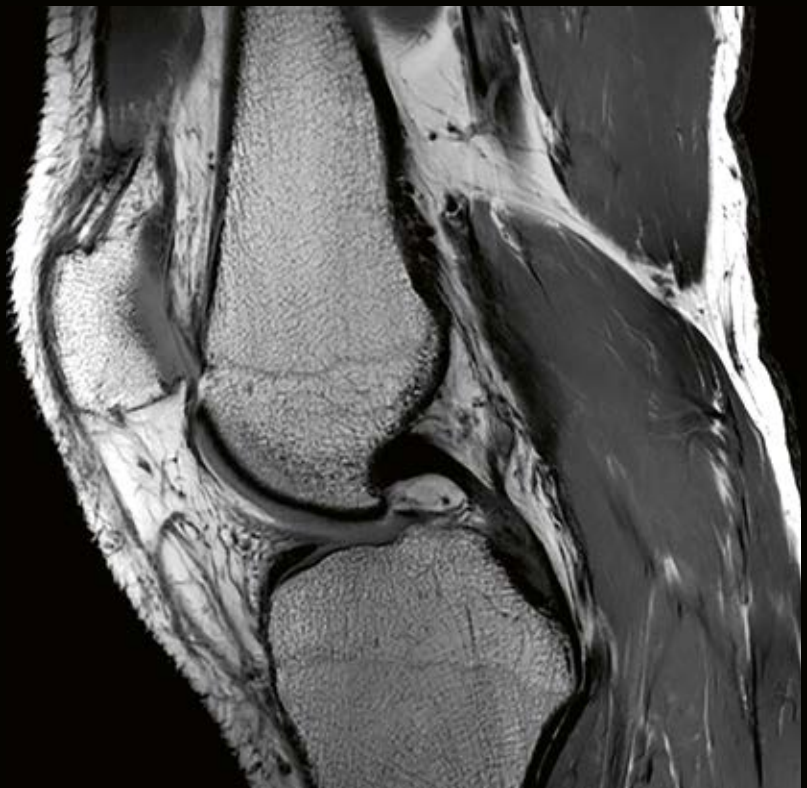

Clinical Mode – Healthy volunteer



PD qTSE FS

0.2 x 0.2 x 2.5 mm³,
3:15 min



T1 FL3D WE

0.5 x 0.5 x 0.5 mm³,
4:35 min



Change the game in UHF business

Medical research funding has stagnated in the last decade. Ensuring your high-end MRI endeavors have the right business impact is crucial in today's competitive environment. MAGNETOM Terra is the result of over 25 years of Siemens UHF innovations, culminating in the design of a brand-new, volume-produced 7T magnet. The magnet is 50% lighter than previous generations and supports easier integration into clinical environments. MAGNETOM Terra⁷ can help you become more competitive, while making a tangible difference to clinical care, research – and your business.

"When you talk to other people in the field, it is clear that Siemens has by far the greatest expertise in ultra-high-field imaging."⁴

Professor Rainer Goebel
University of Maastricht & scannexus,
Maastricht, The Netherlands

Change the game in UHF business with Siemens Healthineers' 50% lighter 7T magnet

Innovative magnet technology

Siemens Healthineers' 7T magnet is a milestone in MR magnet technology. Its unique design and thermally balanced materials minimize physical interactions between core components. The result is 50% lighter than previous generations, with a higher structural stability and a greater fundamental stress capacity. In addition, excellent homogeneity makes for enhanced image quality.

Easy clinical integration

Thanks to the lighter magnet, the scanner can be shipped cold via airfreight. What's more, you benefit from up to 50% faster installation time and ramp-up. Zero Helium boil-off translates into lower lifecycle costs and an improved eco-footprint. All this has the potential to enhance performance, lower resource consumption, improve sustainability, and reduce operating costs.

Increased competitiveness

MAGNETOM Terra⁷ can help you broaden research funding opportunities, making your institution stand out as a leader in life sciences. By being at the cutting edge of clinical care and research, you have the opportunity to increase competitiveness for grants, benefit from reduced complexity in clinical trials, and open up potential for clinical imaging reimbursements.

Forward-looking technology

An investment in MAGNETOM Terra⁷ is an investment in the future. Siemens is committed to serving the ultra-high-field community – today and tomorrow – with a host of outstanding innovations. From development and production, to service – all of MAGNETOM Terra's⁷ key components are delivered from a single reliable partner you can trust, for maximum peace of mind.



**First
7T scanner**
released for
clinical use

50%
lighter magnet
technology

**Lower weight and
cold-shipment**
for easy
integration

Zero
Helium boil-off

Proven innovations in the development and production of magnet technology

1980

The world's first superconducting whole-body MRI

1994

The world's first open MRI magnet

1989

The world's first 1.5T active-shielded magnet


1997

The world's first 3T active-shielded magnet

Award-winning development and production

Siemens Magnet Technology in Oxford, UK, has received seven Best Factory Awards and seven Queen's Awards for Enterprise in multiple categories, including for processes and design. The facility deploys leading-edge supply chain management methods, and prides itself on reliable, robust production and the highest standards of quality.



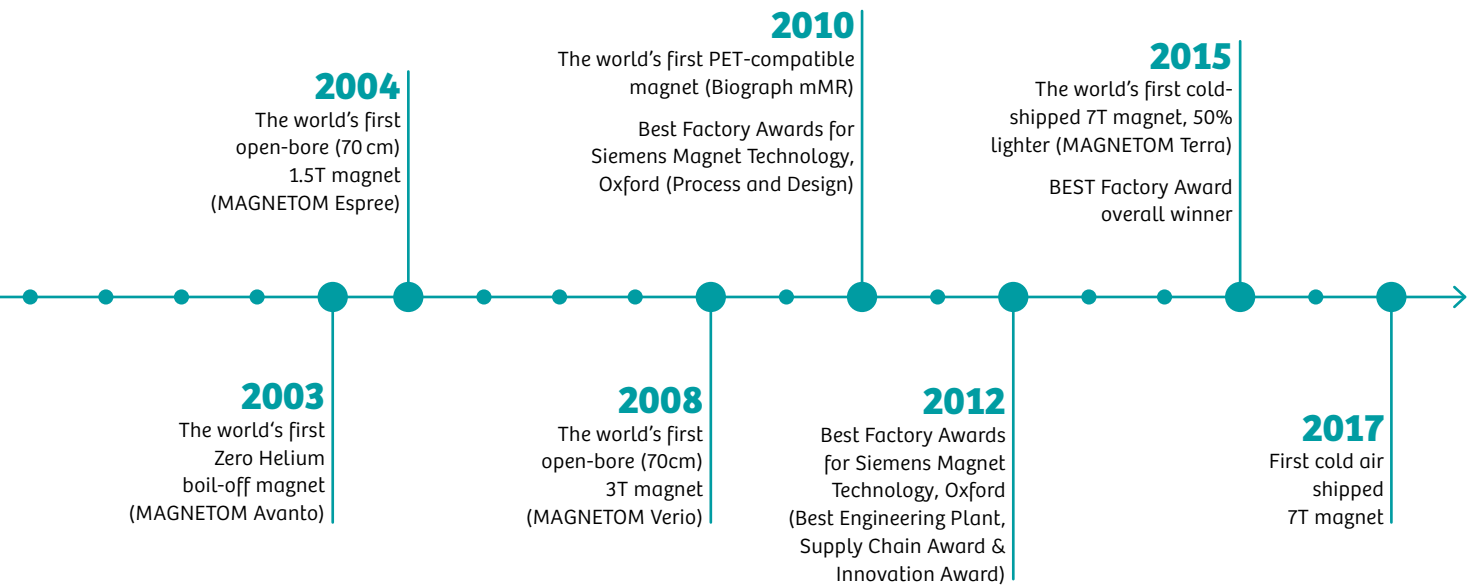

April 2015 – Installation of Siemens' first 7T magnet

During the 30 years that we have been producing 1.5T and 3T magnets, we have gained extensive engineering skills and well-founded process expertise. This knowledge and experience has led to the development and production of our own 7T magnet.



Universitätsklinikum Erlangen, Germany

"We are extremely proud at Siemens Magnet Technology to have developed the 7T magnet at the heart of the MAGNETOM Terra. Once again our expert design and process teams have demonstrated how their innovative thinking has led to a product that has pushed forward the boundaries of magnet technology. It is wonderful to see how seamlessly the manufacture of this flagship product has already been integrated into our award winning facility."

Ralph Seidler
Managing Director,
Siemens Magnet Technology



Unlock research beyond clinical limits

When it comes to research, the freedom to push the boundaries is imperative for gaining a competitive edge. For neuroscience and clinical research applications, MAGNETOM Terra delivers ultra-high SNR and 8-channel pTX for imaging challenging body regions. In addition, it has up to 64 receive channels for higher acceleration factors and 80/200 gradients for maximum flexibility. Moreover, this powerful, reliable scanner supports basic research by helping you develop groundbreaking technologies, set new trends and translate your outcomes into clinical routine.

“The increased spatial resolution offered by 7T MRI enables us to study fine-grained activation patterns within cortical areas and investigate detailed functional topography of the cerebral cortex in individual human subjects. This will provide us with a deeper understanding of the human brain and its connectomics in healthy and diseased populations.”⁴

Professor Kamil Ugurbil
Director of the Center for Magnetic Resonance Research (CMRR),
Minneapolis, Minnesota, USA

Unlock research beyond clinical limits with 8-channel parallel transmit

Enhanced images with pTX

Image quality and speed are key, but inhomogeneities may present challenges, for example, in body MRI. MAGNETOM Terra's 8-channel pTX technology helps you overcome these issues and generate images of excellent quality. This particularly promising technology has the potential to support your own hardware developments.

Ultra-high resolution spectroscopy


Proton magnetic resonance spectroscopy at 7T not only delivers metabolic information, but also gives accurate anatomical insight. Ultra-high 0.14 cm³ resolution has the potential to reveal valuable new diagnostic information for clinical applications – including patients with tumors, epilepsy, multiple sclerosis or other neurodegenerative diseases.

More power for your research

MAGNETOM Terra offers a host of cutting-edge research functionalities, providing access to works in progress packages and powerful hardware configurations. 80/200 gradients and up to 64 receive channels deliver enhanced capabilities for your studies. What's more, the scanner gives you the freedom to explore and develop new clinical applications only possible at ultra-high-field strengths.

Open platform architecture

MAGNETOM Terra provides a flexible, fertile ground for your own UHF hardware and software developments. For example, Siemens collaboration partners benefit from technical support and direct access to the sequence, the Image Calculation Environment (ICE), and imaging protocols.

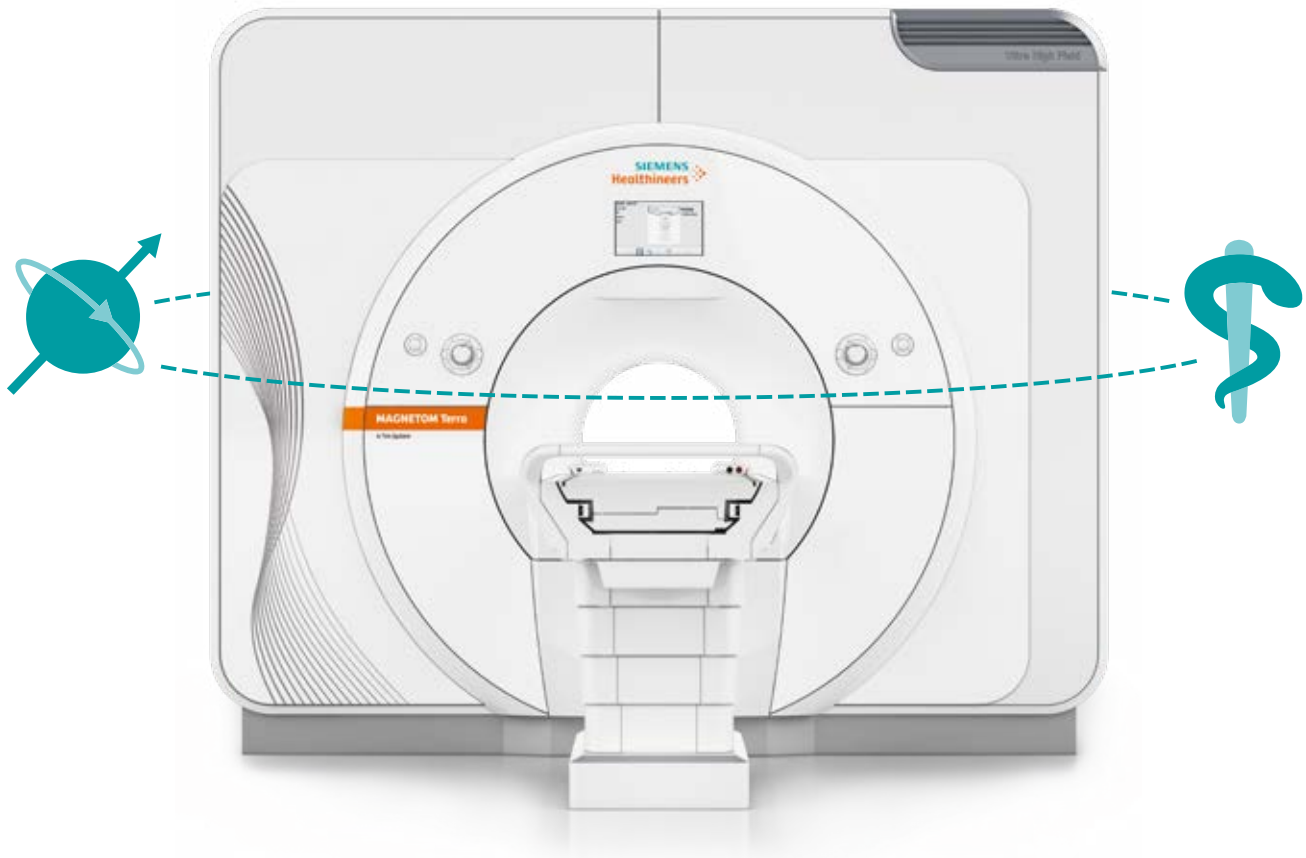

8-channel
pTX for higher
homogeneity

80 / 200
gradients and
64 receive channels
for more research
power

Open
platform
architecture
for own
developments


0.14 cm³
voxel size in
spectroscopy for
metabolic brain
mapping¹



Transparent fibers

Tracks calculated with spherical deconvolution based on diffusion-weighted EPI acquisitions with 1mm isotropic resolution covering the whole brain. The high SNR provided by 7T allows resolving crossing fibers in many brain sub-regions.

Max Planck Institute, Leipzig, Germany

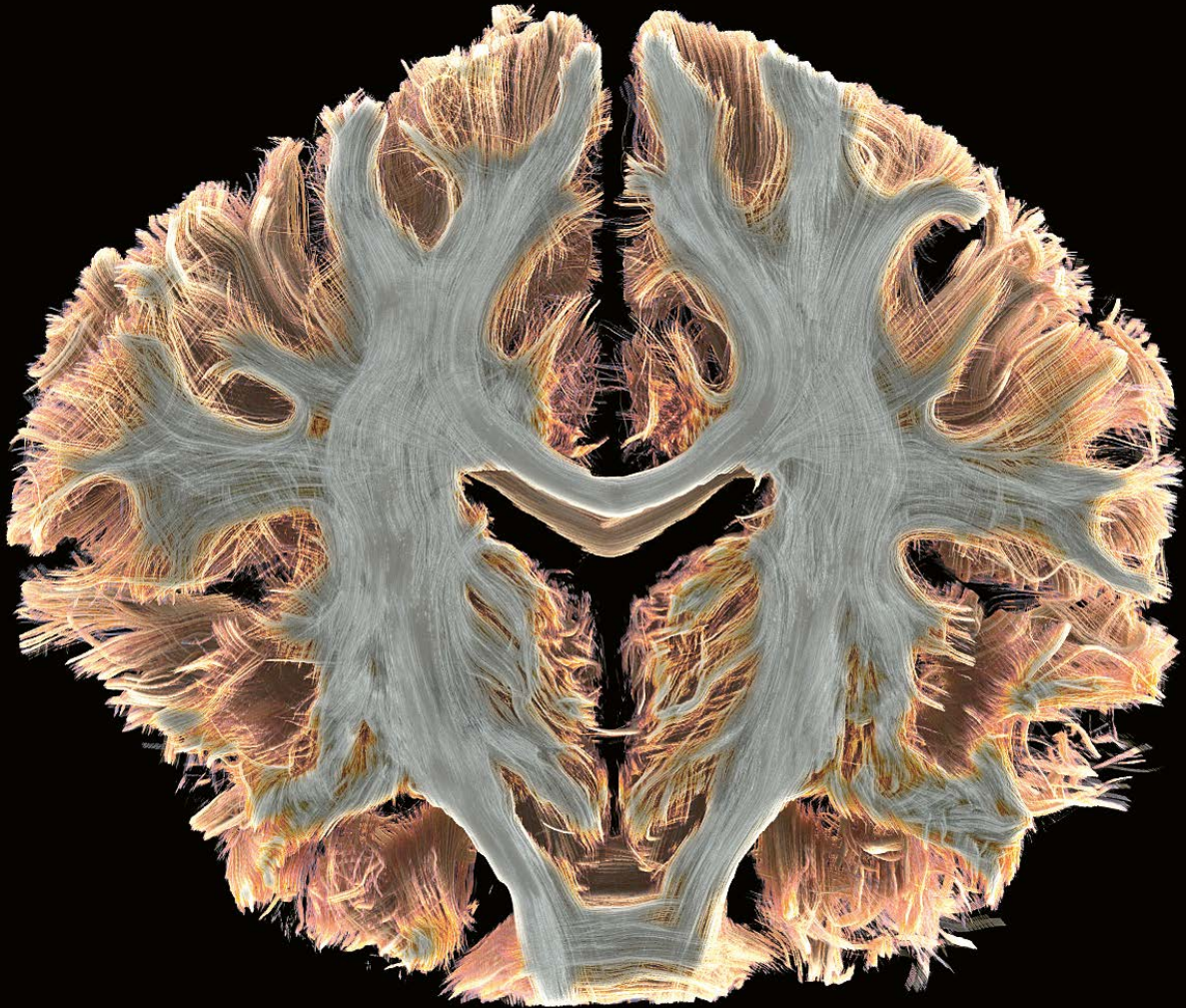

In-vivo histology


syngo.via Frontier, the research extension of *syngo.via*, helps bridge the gap in post-processing translational research. Cinematic rendered images based on MR data sets may be used for patient counseling, surgery planning, or teaching purposes.⁵



Transparent fibers


Tracks calculated with spherical deconvolution based on diffusion-weighted EPI acquisitions with 1 mm isotropic resolution covering the whole brain. The high SNR provided by 7T allows resolving crossing fibers in many brain sub-regions.

Max Planck Institute, Leipzig, Germany

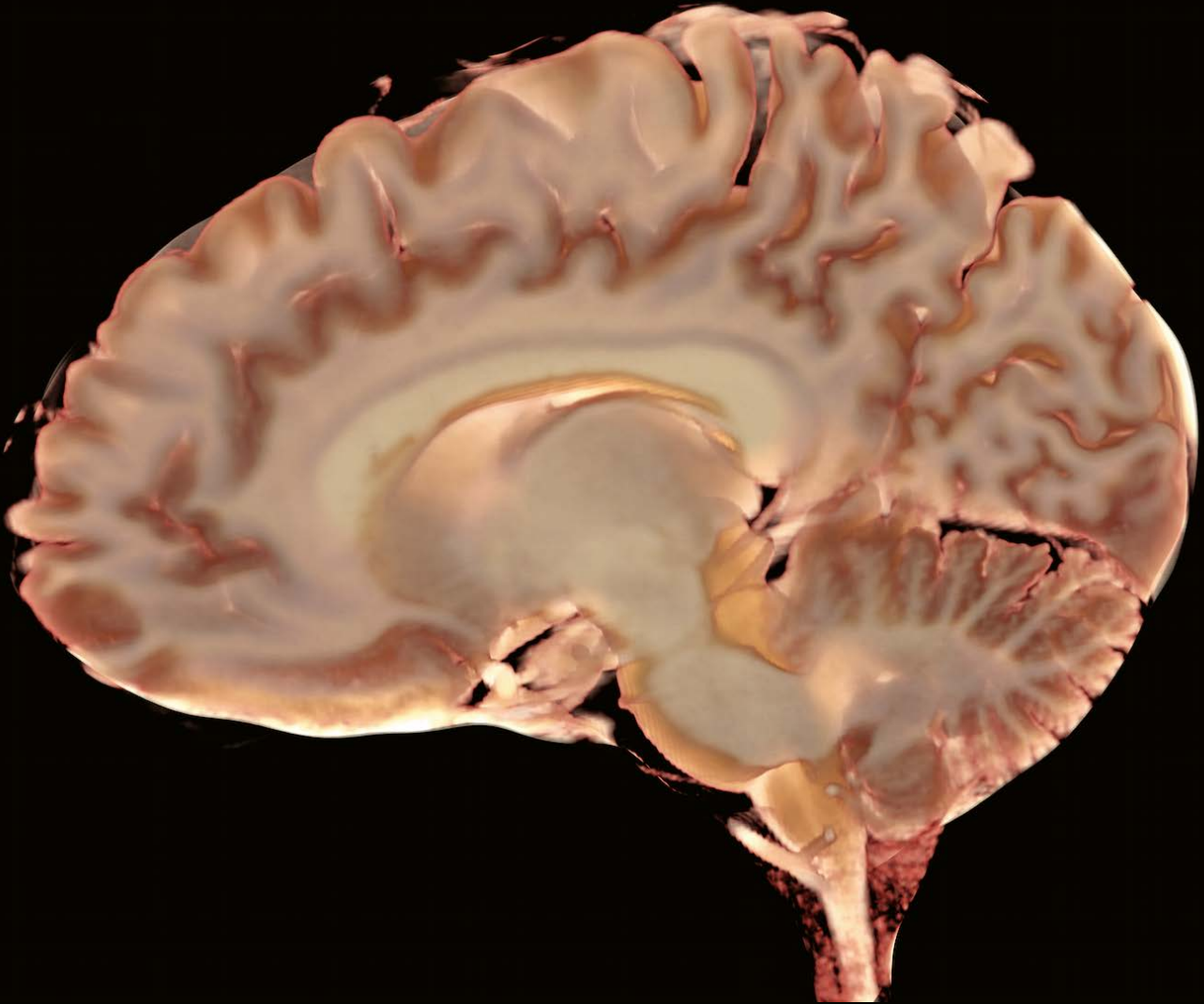


In-vivo histology

syngo.via Frontier, the research extension of *syngo.via*, helps bridge the gap in post-processing translational research. Cinematic rendered images based on MR data sets may be used for patient counseling, surgery planning, or teaching purposes.⁵

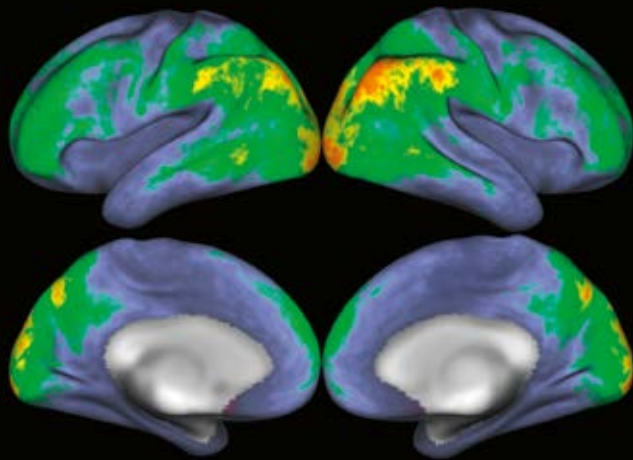

3 Tesla

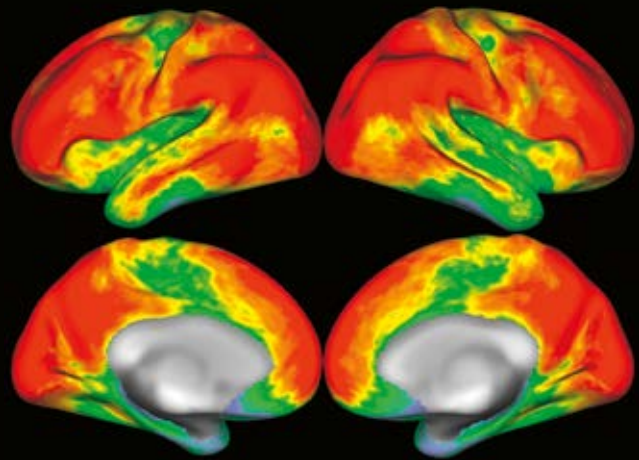


7 Tesla



Functional MRI

Contrast-to-noise ratio maps in resting state fMRI.

*Consortium The Human Connectome Project.
CMRR, Minnesota, USA
Washington University St. Louis, USA
Oxford University, UK*

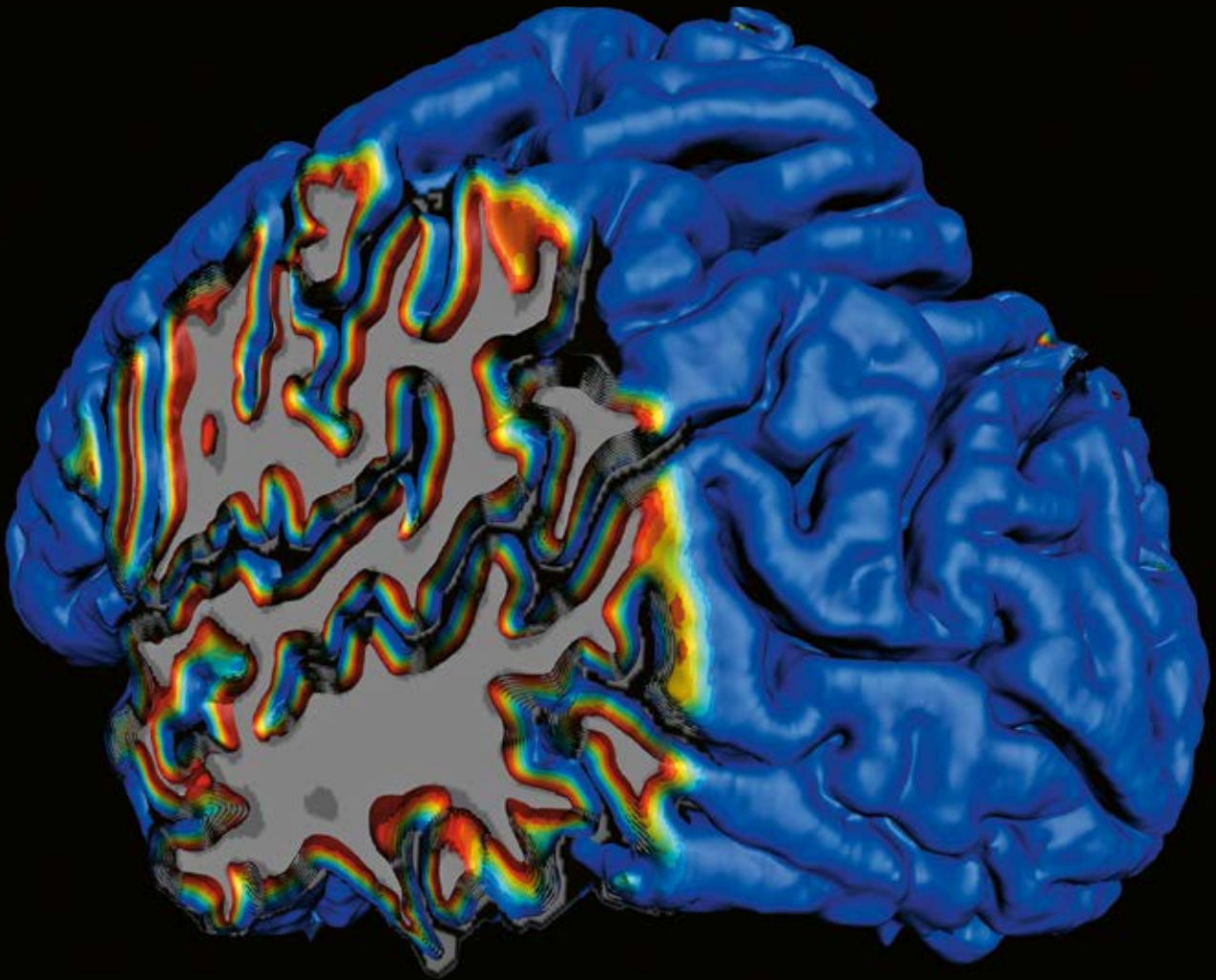

Depiction of cortical layers

Post-processed high-resolution anatomical MR data reveals reconstructed surfaces at different cortical depth levels. The inner red surface runs along the white/gray matter boundary. The outer blue surface runs along the outer (pial) boundary of the cortex.


Scannexus, Maastricht, The Netherlands



Tissue segmentation

Delineation of subcortical nuclei in the thalamus and brainstem at 1mm isotropic resolution. White matter nulled MPAGE (top TA 8:52 min) and gray matter nulled MPAGE (bottom TA 10:38 min).

MGH, Boston, USA

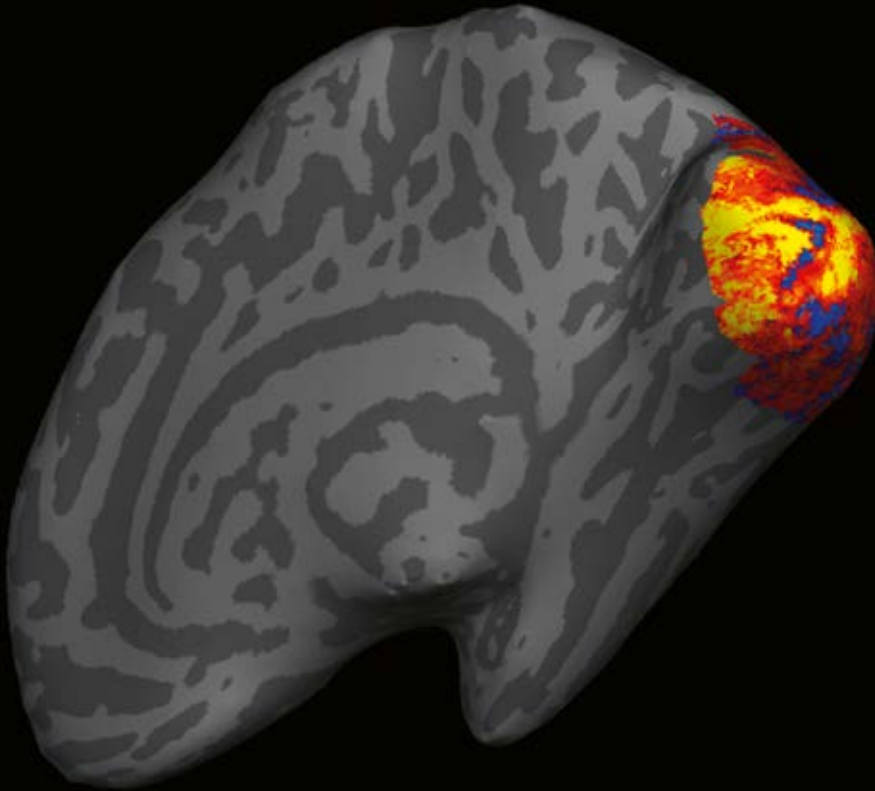

High-resolution fMRI

Cortical-layer-specific activation with fMRI at 1 mm isotropic resolution, inflated view. The fMRI visual stimulus was designed to activate a pattern in the shape of the number "7" using the known retinotopic mapping in the human visual cortex.

MGH, Boston, USA

*"Siemens provides the best possible open hardware and software environment to explore these new transmit and receive concepts, all of which have proved essential to allow ultra-high fields to fulfill their potential for the benefit of human health."*⁴

Professor Lawrence L. Wald, Director
MGH NMR Core at Martinos Center, Department of Radiology,
Boston, Massachusetts, USA



Spine imaging

Ultra-high resolution of the cervical spine using a custom-built spine coil.

MGH, Boston, USA




3 Tesla 0.5 mm in-plane



7 Tesla 0.5 mm in-plane




7 Tesla 0.3 mm in-plane

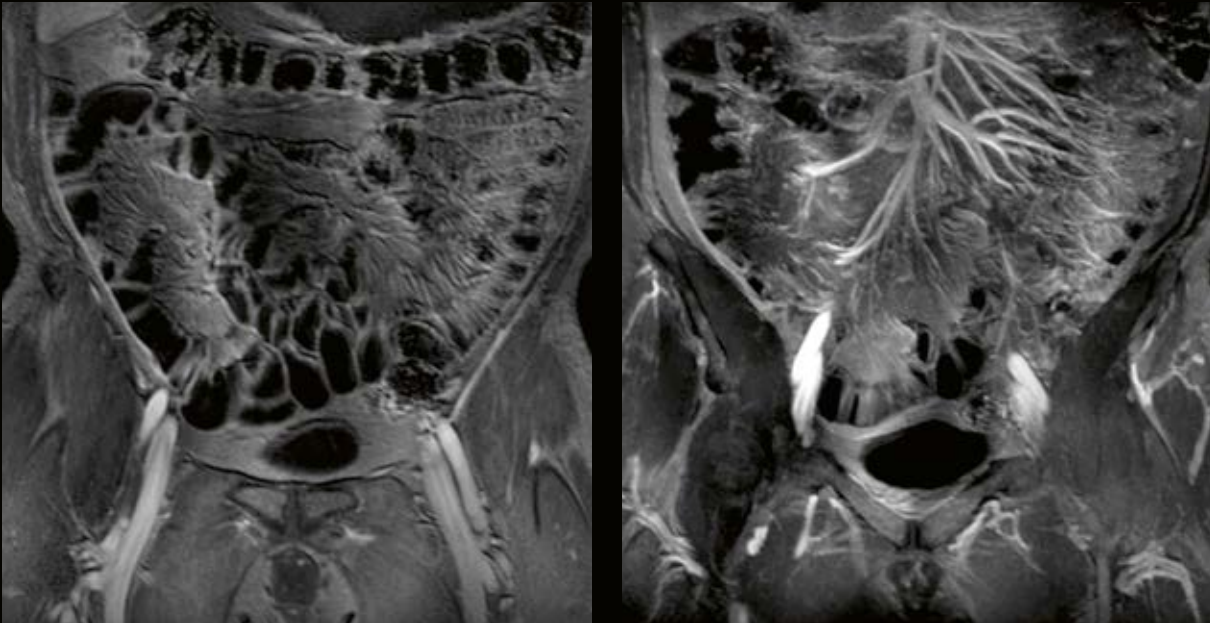

Body imaging

Left image, 3D VIBE FatSat, right image, thin MIP from the 3D VIBE FatSat. Images acquired using pTX and custom-built coils.

Erwin L. Hahn Institute for MRI, Essen, Germany

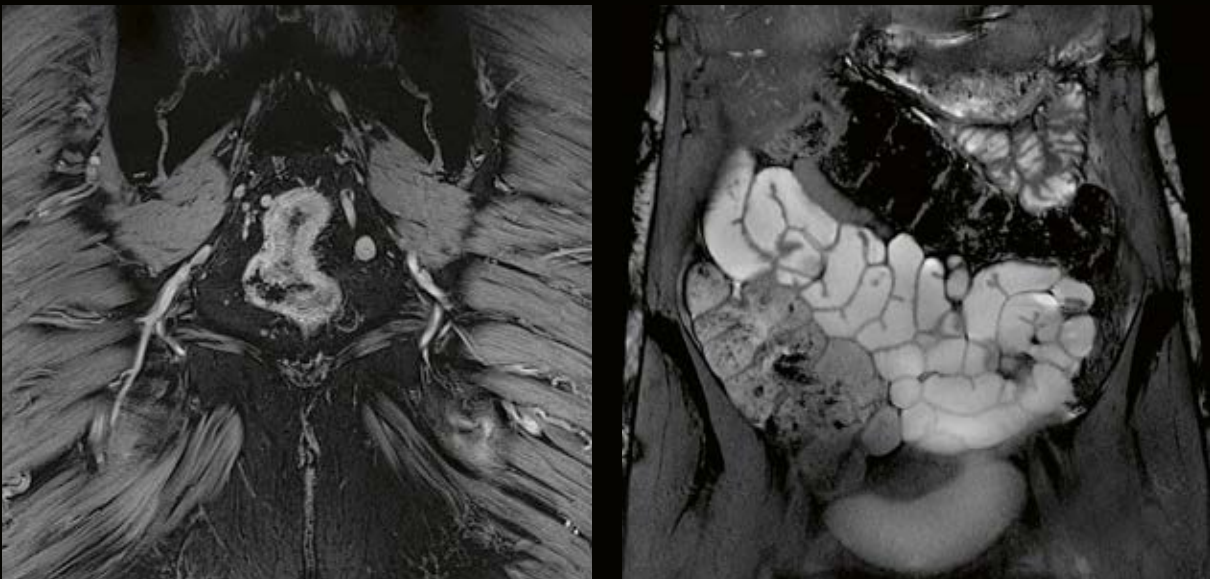


Body imaging

Image left, rectum carcinoma imaged with a ce-FLASH ($0.3 \times 0.6 \times 2 \text{ mm}^3$, TA 2:14 min). Images acquired using pTX and custom-built coils. Image right, abdominal small bowel imaging with TrueFISP ($0.4 \times 0.8 \times 2 \text{ mm}^3$, TA 26 s).


Erwin L. Hahn Institute for MRI, Essen, Germany




Short axis view FLASH cine retro
(1 x 1 x 4 mm³, GRAPPA 2, TA 18 s).

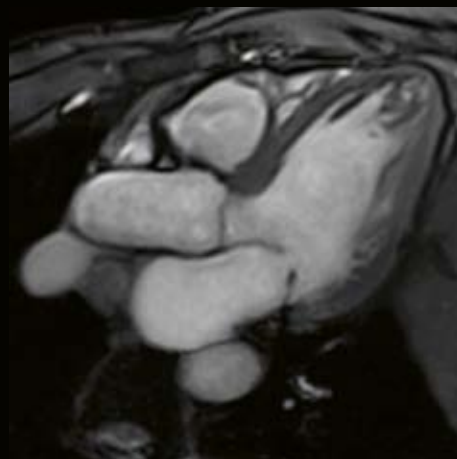

Right ventricular output tract FLASH cine retro
(1.3 x 1 x 4 mm³, SENSE 2, TA 13 s).

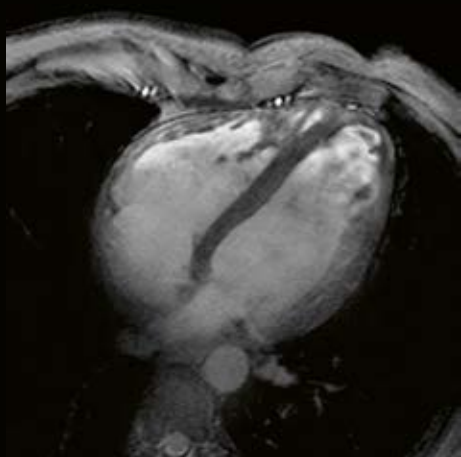

Four chamber view FLASH cine retro
(1 x 1 x 4 mm³, GRAPPA 2, TA 16 s).


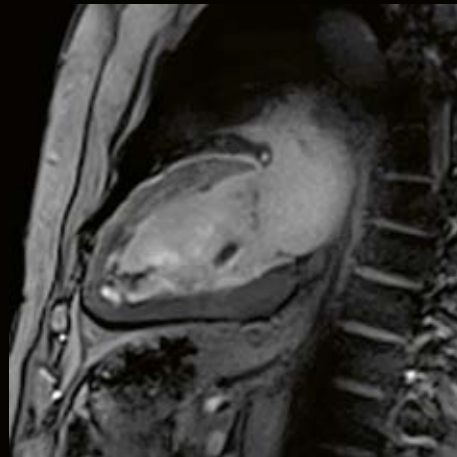


Two chamber view FLASH cine retro
(1.2 x 1 x 4 mm³, SENSE 2, TA 13 s).

Cardiac imaging

Accelerated T1-weighted FLASH
acquisitions.

Berlin Ultrahigh Field Facility, Berlin, Germany

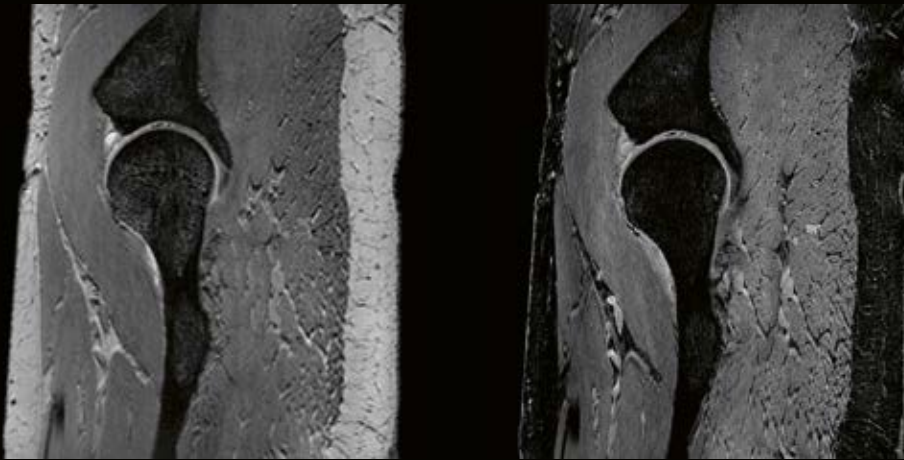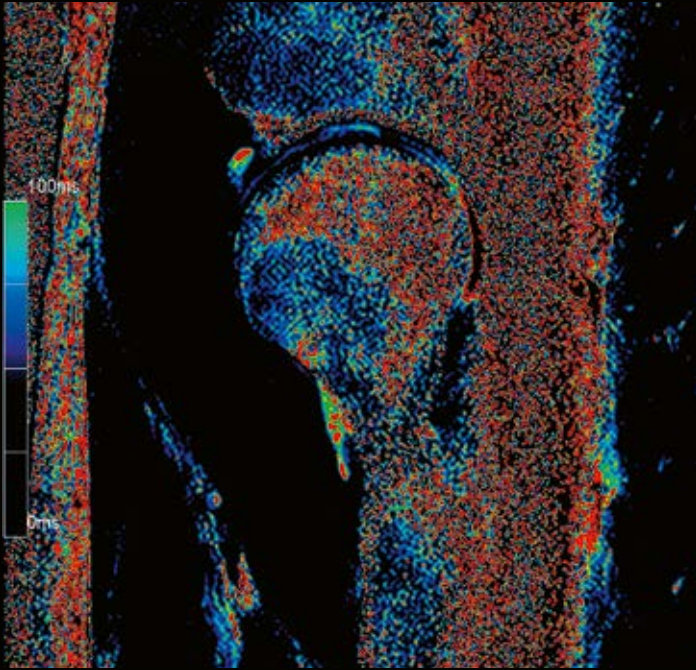


Hip cartilage transplant

Cartilage transplant visible in 3D DESS ($0.37 \times 0.74 \times 0.74 \text{ mm}^3$, TA 5:11), 3D VIBE SPAIR ($0.19 \times 0.39 \times 0.8 \text{ mm}^3$, TA 5:58) and T2 MapIt ($0.25 \times 0.5 \times 2.5 \text{ mm}^3$, TA 4:47) 14 months after Autologous Chondrocyte Transplantation (ACT).


Erwin L. Hahn Institute for MRI, Essen, Germany




Morphologic PD TSE F ($0.4 \times 0.4 \times 2 \text{ mm}^3$)

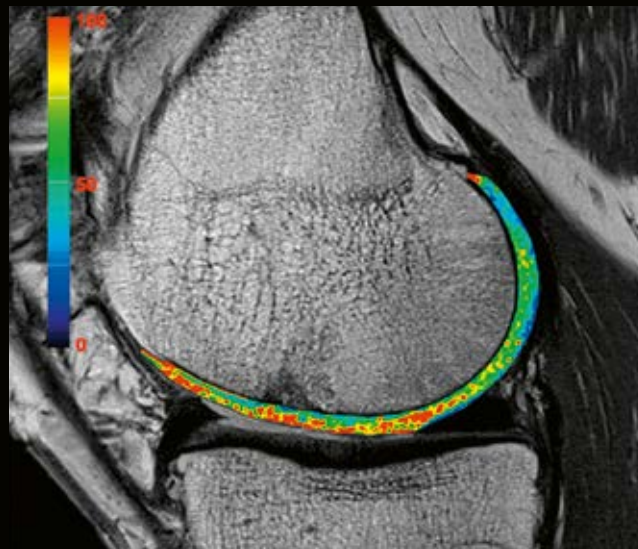

T2 map ($0.6 \times 0.6 \times 1 \text{ mm}^3$), T2 in ms
More water, disturbed collagen architecture visible

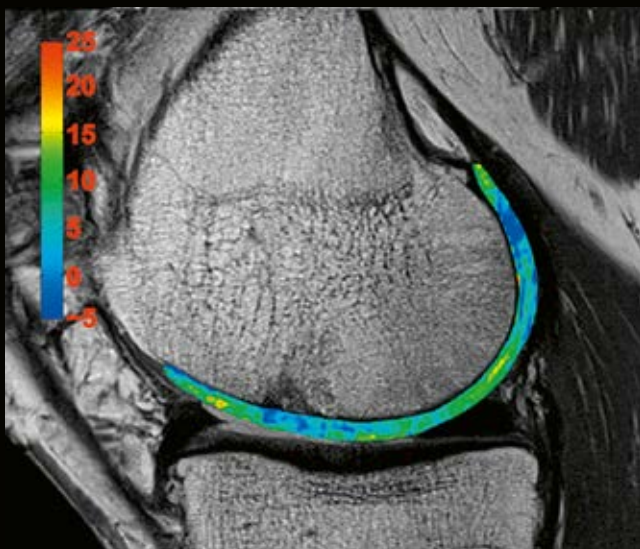

gagCEST image ($0.7 \times 0.7 \times 3 \text{ mm}^3$)
gagCEST asymmetries in [%] lower values, less PG content

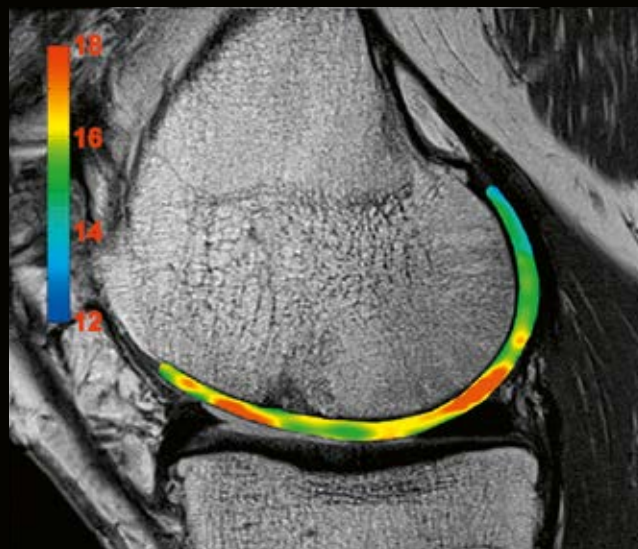


Sodium image
Sodium SNR lower values, less PG content


Biochemical imaging using CEST

Male patient nine years after Autologous Osteochondral Transplantation (AOT) in the medial femoral condyle.


MedUni Wien, Vienna, Austria

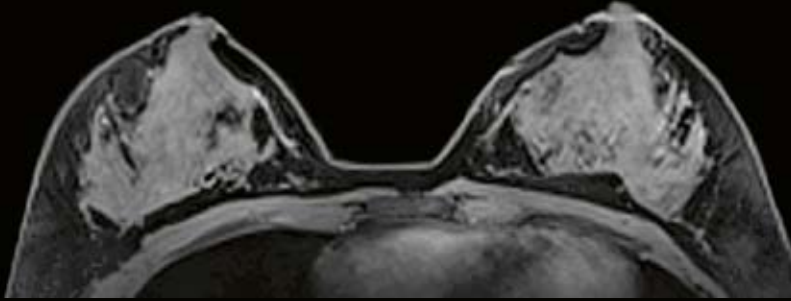

Pre-contrast

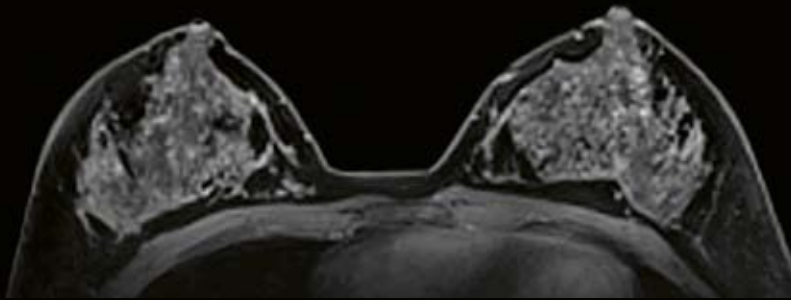

Post-contrast

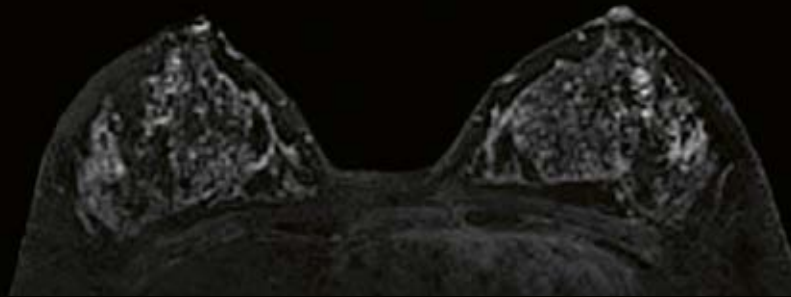

Subtraction

Breast imaging

High-resolution bilateral breast imaging in short acquisition time. Pre- and post-contrast 3D FLASH SPAIR, 1.4 mm isotropic resolution, TA 1:48 min/series.

NYU Langone Medical Center, New York, USA



Join the largest research community

Your reputation plays a pivotal role in your institution's success. MAGNETOM Terra has the power to let you go deeper than ever before, making your research and patient outcomes stand out from the rest. What's more, this leading-edge technology can help you attract the brightest minds to your facility, further enhancing your capabilities. MAGNETOM Terra has the potential to put your organization firmly on the map, offering access to an exclusive network of expertise and broad scope for collaboration and exchange.

"When we were in a position to order a 7T system, Siemens was the logical choice."⁴

Professor Peter Jezzard
Professor of Neuroimaging,
University of Oxford, Oxford, UK

Join the largest research community with over 70% of all UHF users

Enhance your reputation


MAGNETOM Terra helps you achieve your research goals, giving you the opportunity to publish first and become a true opinion leader. This advanced technology has the potential to strengthen your position by attracting the brightest brains to your facility. The scanner lets you deliver previously unseen insights that could improve patient outcomes and further enhance your reputation.

Expand your network


Even if you are taking your first steps in ultra-high-field imaging, you will never be alone. Siemens has proven expertise in UHF MRI and cultivates links with an extensive network of users. As a result, you benefit from the experience of others and can share your own ideas. Siemens is the global leader in 7T – with a market share of over 70% and more than 25 years of experience in this field.

Exchange your ideas with peers

When you become part of the Siemens UHF community, you join an exclusive group of outstanding MRI experts. Through collaboration and exchange with other leaders in your field, you can extend your own knowledge and gain deeper insights. Siemens' regular user meetings and an online discussion board are the ideal platforms to interact with your peers.

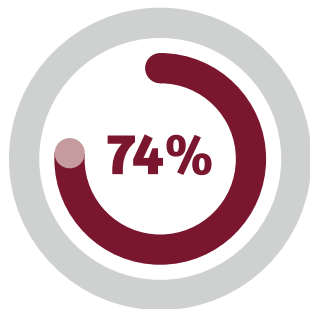


74% of ISMRM UHF abstracts
in 2017 were based on data
from Siemens UHF systems.²

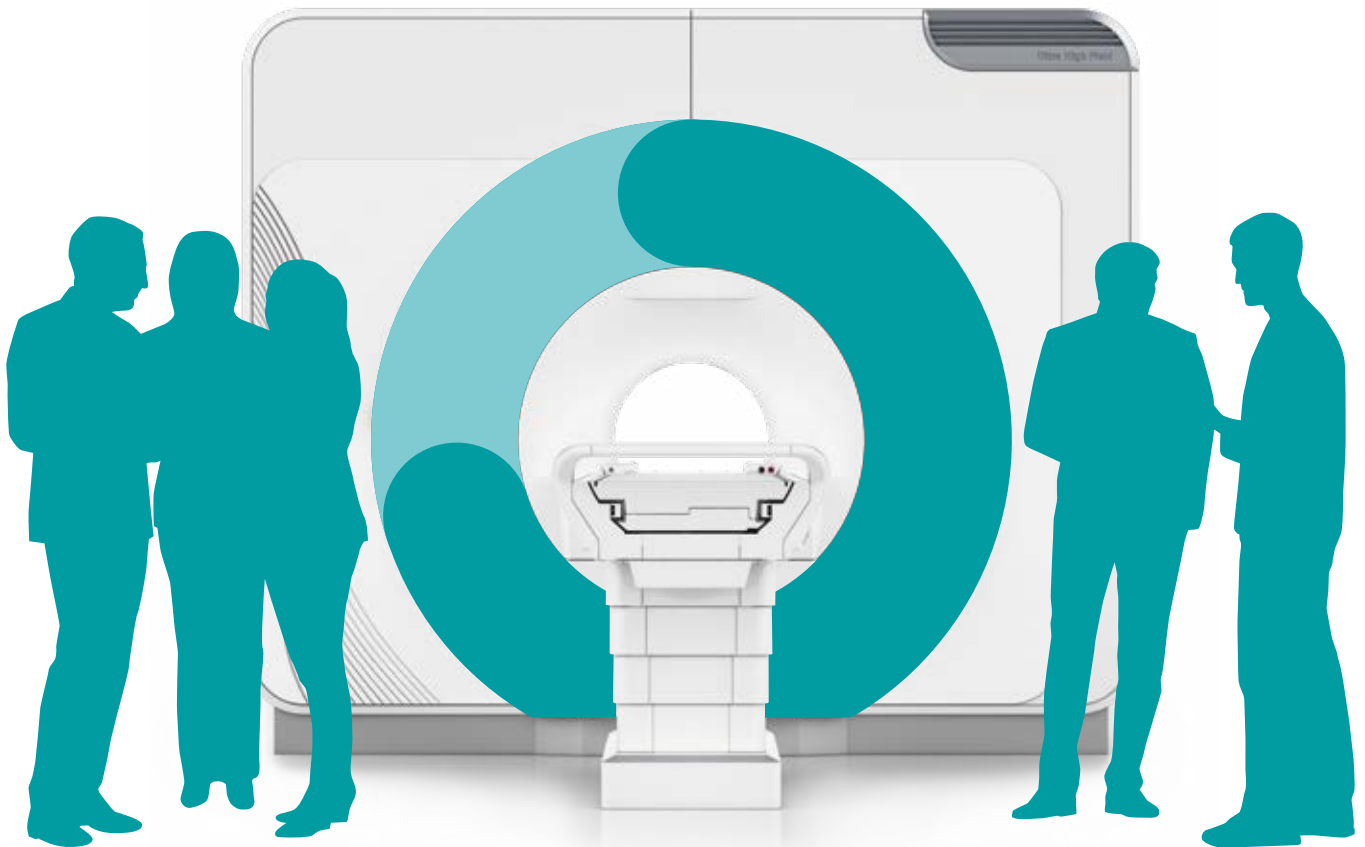

Over **70%**
of 7T scanners
deployed
worldwide are
from Siemens

Attract and
retain the
brightest minds and
**publish
first**

Strong
network
for collaboration
and peer-to-peer
exchange

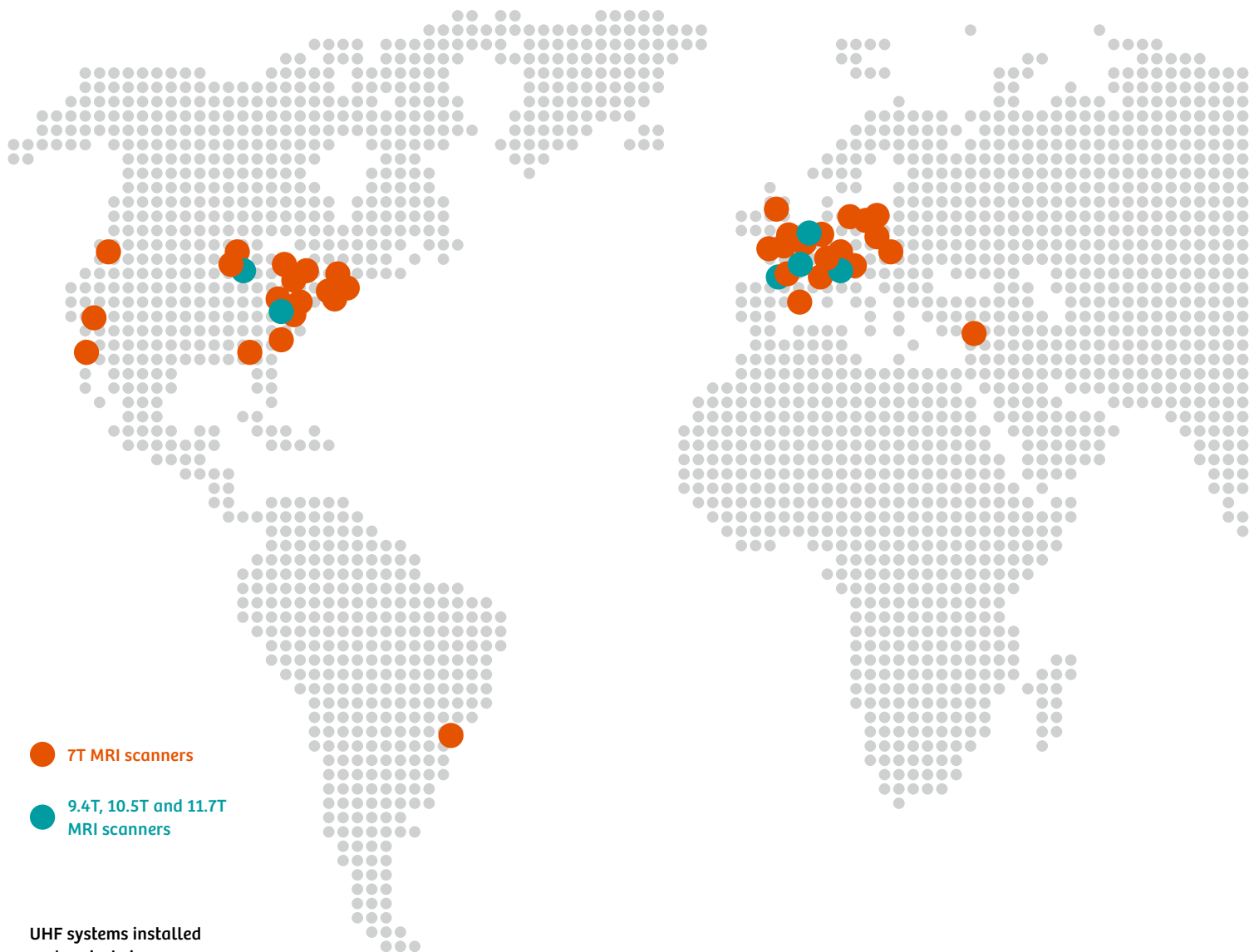


● 7T MRI scanners

● 9.4T, 10.5T and 11.7T MRI scanners

UHF systems installed and projects in progress

- 1 Athinoula A. Martinos Center for Biomedical Imaging of MGH, Boston, Massachusetts, USA
- 2 Leibniz Institute for Neurobiology (LIN), Magdeburg, Germany
- 3 Bernard and Irene Schwartz Center for Biomedical Imaging (CBI) of New York University Langone Medical Center, New York City, New York, USA
- 4 Center for MR Research (CMRR), University of Minnesota, Minneapolis, Minnesota, USA
- 5 Neuroscience Research Institute (NRI) of Gachon University of Medicine and Science, Incheon, South Korea
- 6 Advanced Imaging Research Center (AIRC), Oregon Health & Science University, Portland, Oregon, USA
- 7 Erwin L. Hahn Institute for Magnetic Resonance Imaging (ELH), Essen, Germany
- 8 Center for Imaging in Biomedicine (CIBM), École polytechnique fédérale de Lausanne (EPFL), Lausanne, Switzerland
- 9 Max Planck Institute for Biological Cybernetics (MPI KYB), Tübingen, Germany (9.4T)
- 10 NeuroSpin, French Alternative Energies and Atomic Energy Commission (CEA), Saclay, France
- 11 NeuroSpin, French Alternative Energies and Atomic Energy Commission (CEA), Saclay, France (11.7T)
- 12 Magnetic Resonance Research Center (MRRC), University of Pittsburgh Medical Center (UPMC), Pittsburgh, Pennsylvania, USA
- 13 Max Planck Institute for Human Cognitive and Brain Sciences (MPI), Leipzig, Germany
- 14 Excellence Center for Highfield MR, Medical University of Vienna (MUW), Vienna, Austria
- 15 German Cancer Research Center (DKFZ), Heidelberg, Germany
- 16 Institute of Neuroscience and Medicine (INM), Research Centre Jülich, Jülich, Germany (9.4T)
- 17 Center For Magnetic Resonance And Optical Imaging (MMRCC), University of Pennsylvania Health System (HUP), Philadelphia, Pennsylvania, USA
- 18 Berlin Ultrahigh Field Facility (B.U.F.F.), Experimental and Clinical Research Center (ECRC), Berlin, Germany
- 19 State Key Laboratory of Brain and Cognitive Science, Institute of Biophysics, Chinese Academy of Sciences (CAS), Beijing, China
- 20 Oxford Centre for Functional MRI of the Brain (FMRIB), University of Oxford, Oxford, UK

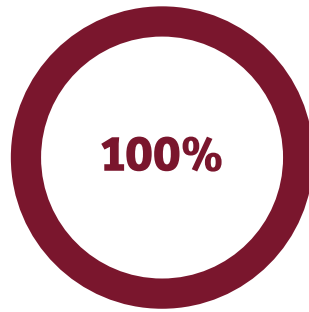

100% of vendor-integrated human MRI scanners with a field strength higher than 7T are from Siemens.




Over 70% of 7T scanners deployed worldwide are from Siemens.

- 21** Magnetic Resonance Imaging Research Center, Auburn University, Auburn, Alabama, USA
- 22** Center for MR Research (CMRR), University of Minnesota, Minneapolis, Minnesota, USA
- 23** Functional MRI Facility (FMRIF), National Institute of Mental Health and Neurological Disorders and Stroke, National Institutes of Health (NIH-NIMH & NINDS), Bethesda, Maryland, USA
- 24** National Institute of Neurological Disorders and Stroke, National Institutes of Health (NIH-NINDS), Bethesda, Maryland, USA (11.7T)
- 25** National Institute of Information and Communication Technology (NiCT) / Center for Information and Neural Networks (CiNET), Osaka, Japan

- 26** Center for MR Research (CMRR), University of Minnesota, Minneapolis, Minnesota, USA (10.5T)
- 27** Center for Imaging of Neurodegenerative Diseases (CIND), San Francisco VA Medical Center, UCSF, San Francisco, California, USA
- 28** German Center for Neurodegenerative Diseases (DZNE), Bonn, Germany
- 29** Biomedical Research Imaging Center (BRIC), University of North Carolina (UNC), Chapel Hill, North Carolina, USA
- 30** Maastricht Brain Imaging Centre (M-BIC), Maastricht University, Maastricht, The Netherlands
- 31** Maastricht Brain Imaging Centre (M-BIC), Maastricht University, Maastricht, The Netherlands (9.4T)
- 32** Mt Sinai School of Medicine, New York City, New York, USA
- 33** Cleveland Clinic, Cleveland, Ohio, USA
- 34** Centre for Advanced Imaging, University of Queensland, Brisbane, Queensland, Australia
- 35** Royal Melbourne Hospital, University of Melbourne, Victoria, Australia
- 36** University of Sao Paulo (USP), Sao Paulo, Brazil
- 37** Centre d'Exploration Métabolique par Résonance Magnétique (CIMEREM), Marseille, France
- 38** Centre for Functional and Metabolic Mapping, Robarts Research Institute, London, Ontario, Canada
- 39** National Institute for Physiological Sciences (NIPS), Okazaki, Japan

- 40** Kyoto University, Kyoto, Japan
- 41** Zhejiang University, Hangzhou, China
- 42** Brigham and Women's Hospital (BWH), Boston, USA
- 43** University of Southern California (USC), Los Angeles, California, USA
- 44** Cardiff University Brain Research Imaging Centre (CUBRIC) Cardiff, UK
- 45** Wolfson Brain Imaging Centre (WBIC), University of Cambridge, Cambridge, UK
- 46** Imaging Centre of Excellence (ICE), South Glasgow University Hospital, Glasgow, UK
- 47** Magnetic Resonance Research Center (MRR), Yale University, New Haven, Connecticut, USA
- 48** Comprehensive Heart Failure Center (CHFC), Würzburg University Hospital, Würzburg, Germany
- 49** Weizmann Institute of Science, Tel Aviv, Israel
- 50** Mayo Clinic, Rochester, USA
- 51** Toronto Western Hospital (TWH), University Health Network (UHN), Toronto, Canada
- 52** National Institute of Health, National Institute on Drug Abuse (NIH-NIDA), Bethesda, Maryland, USA
- 53** Forschungszentrum Jülich, Jülich, Germany
- 54** CRC, University of Liege, Liege, Belgium
- 55** Houston Methodist, Houston, Texas, US
- 56** Athinoula A. Martinos Center for Biomedical Imaging of MGH, Boston, Massachusetts, USA

Service and exchange – Comprehensive services

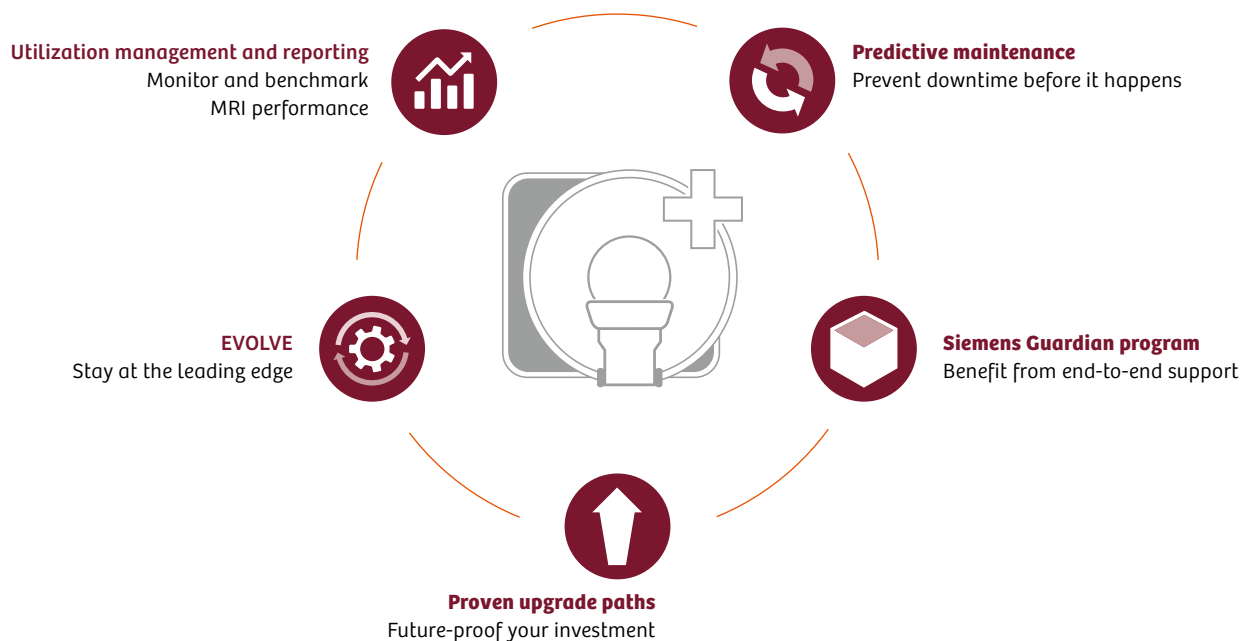

Siemens' end-to-end services ensure you stay at the leading edge of MRI technology throughout the entire system lifecycle – from installation, to operation, to upgrades, to ongoing support. Moreover, our diverse communication platforms and communities keep you up to speed on the world of MRI and enable you to share your ideas and experiences with your peers.



Utilization management and reporting

This powerful solution gives you more from your MRI scanner. It allows you to monitor KPIs and benchmark your system against other Siemens MRI machines at any facility or organization. So you can keep track of your MRI performance, and reap the maximum reward from your scanner.

Predictive maintenance

When systems go down, it impacts both your ability to care for your patients and your bottom line. Siemens provides a predictive maintenance service to help you minimize lost time. It informs you when a part of your MRI system is likely to fail, enabling you to plan repairs and prevent downtime before it happens.

EVOLVE

Keep your hardware and software up to date at all times – a key factor in enhancing performance and diagnostic quality. You receive all applicable upgrades for software and the *syngo* OS, plus at least one workstation hardware upgrade within the first six years.


Siemens Guardian program

This program provides the latest service technology so you can better manage your MRI system. It combines many features in a single package – offering real-time system monitoring, expert advice to improve workflow efficiency, proactive maintenance, and support. Moreover, it guarantees defined repair times, giving you complete peace of mind.

Proven upgrade paths

With MAGNETOM scanners, taking your MRI system to the next level is simplicity itself, thanks to clearly defined upgrade paths. In fact, Siemens has built an entire organization (CDV) to help customers truly maximize their system life – and increase their return on investment.

Service and exchange – Peer-to-peer information

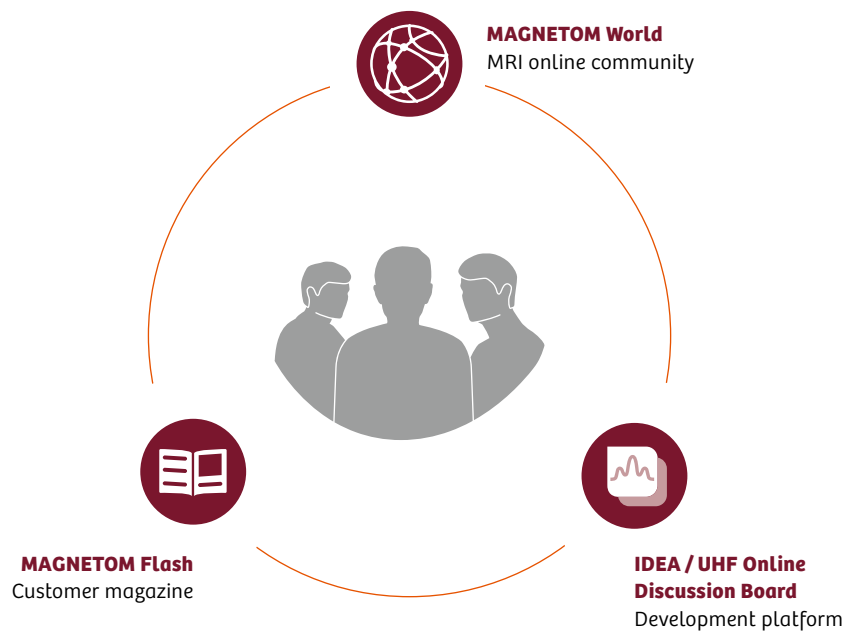

On MAGNETOM Flash:

“An excellent and useful combination of technological and clinical articles that both keep one up to date with advances in MRI and provide practical assistance for day-to-day practice – good and interesting learning material.”⁴

Mark Lourensz
St Vincent's Hospital, Fitzroy, Victoria, Australia



MAGNETOM World

Siemens Healthineers' global MRI community offers peer-to-peer support and information. Radiologists, cardiologists, technologists, and physicists have all contributed with publications, presentations, training documents, case studies, and more – all freely available to you via this unique network. Plus, the bi-annual MAGNETOM World Summit is the ideal opportunity to share and exchange ideas.

MAGNETOM Flash

Published quarterly, the MR customer magazine features up-to-date clinical case studies, application tips and technical and product information relevant to you. All content is carefully compiled by experts to meet the needs of today's MRI users in both clinical and research scenarios. In fact, 98.5% of readers report that MAGNETOM Flash is clinically relevant.

IDEA / UHF Online Discussion Board

IDEA⁷ is an open development platform for the largest and most active 3T and UHF research communities in the world. It unites users from across the globe and fosters innovation in the field of MRI. Members collaborate online at www.mr-idea.com and at an annual meeting. IDEA includes an exclusive area, the UHF Online Discussion Board, to help users focus on topics of interest, as well as find and communicate with the right peers.

Visit MAGNETOM World

siemens.com/magnetom-world

Technical specifications

MAGNETOM Terra Technical specifications

Field strength	7 Tesla
Bore size	60 cm
Magnet length	270 cm
System length ⁸	297 cm
System weight (in operation) ⁸	< 25 tons
Minimum room size ⁸	65 m ² / (w/o pTX and 3rd order shim option)
Dual Mode functionality	Clinical Mode Research Mode
RF transmit	TimTX-1, TimTx-8 (in research mode)
RF receive	Tim [32 x 32] [64 x 64]
Gradient strength	XR gradients (80 mT/m @ 200 T/m/s)
Helium consumption	Zero Helium boil-off technology
Quiet Suite/DotGO	Available in research mode
Local coils in clinical mode	1TX/32RX head coil 1TX/28RX knee coil



Why Siemens Healthineers?

At Siemens Healthineers, our focus is to help healthcare providers succeed in today's dynamic environment.

Healthcare providers around the world have long relied upon our engineering excellence – leading-edge, high-quality medical technologies across a broad portfolio. Our technologies touch an estimated 5 million patients globally every day.¹⁰ At the same time, they help hospital departments to continuously improve their clinical, operational, and financial outcomes.

We now consolidate this unprecedented volume of data and insights and turn them into pioneering enterprise and digital health services. With those, we maximize opportunities and share risk for the success of your entire health system.

Partnerships are built on people. With Siemens Healthineers, there is no team more committed and more connected than we are to realize your success together.

Siemens Healthineers.
Engineering success. Pioneering healthcare. Together.

Not for distribution in the USA. On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products included in this brochure are available through the Siemens sales organization worldwide. Availability and packaging may vary by country and are subject to change without prior notice. Some/All of the features and products described herein may not be available in the United States. Some products are still under development and not commercially available yet. Their future availability cannot be ensured. The information in this document contains general technical descriptions of specifications and optional features which do not always have to be present in individual cases. Siemens reserves the right to modify the design, packaging, specifications, and options described herein without prior notice. Please contact your local Siemens sales representative for the most current information.

Note: Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced.

For accessories, please visit:
www.siemens.com/medical-accessories

- 1 Scheenen et al., Magn Reson Mater Phy (2008) 21:95–101
- 2 S. Trattnig, et al. Clinical applications at ultrahigh field (7 T). Where does it make the difference? NMR Biomed. 2015.
- 3 Data on file.
- 4 The statements by Siemens' customers described herein are based on results that were achieved in the customer's unique setting. Since there is no "typical" hospital and many variables exist (e.g., hospital size, case mix, level of IT adoption) there can be no guarantee that other customers will achieve the same results.
- 5 Rendered with a Siemens internal cinematic rendering prototype.
- 6 Compared to previous
- 7 This website provided by Siemens AG may be used solely in accordance with the general terms and conditions of use, available prior to registration / login on the website itself.
- 8 Minimum total space requirement for magnet, electronics, and console room
- 9 Siemens AG, "Sustainable healthcare strategy – Indicators in fiscal 2014", pages 3–4.

.....

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen, Germany
Phone: +49 9131 84-0
siemens-healthineers.com