

### APPLICATION

Wrist SPEEDER (MJAJ-162A/S1) is designed for optimal imaging of the hand and wrist using SPEEDER technology (Toshiba's advanced parallel imaging technology).

### APPLICABLE SYSTEM

This coil is applicable to the following system.

- Vantage Titan™ 3T

### Requirements for use

System software V1.20 or later

### COMPOSITION

Wrist SPEEDER.....	1 set
• Wrist SPEEDER Coil.....	1
• Coil base .....	1
• Pads.....	1 set
Operation manual.....	1

### PERFORMANCE SPECIFICATIONS

This array coil permits up to 6 elements to be selected, providing an excellent signal-to-noise ratio.

Hand mode is used for whole-hand imaging with 6 elements. Wrist mode is used for wrist imaging with 6 smaller elements, which results in a higher SNR.

Wrist SPEEDER can be used in both the side position and the overhead position. The overhead position is permitted by setting Wrist SPEEDER on the detached coil base.



### APPLICABLE CODES

IEC60601-1 (1988)	+ Amendment 1 (1991)
	+ Amendment 2 (1995)
IEC60601-1-2 (2001)	+ Amendment 1 (2004)
IEC60601-2-33 (2002)	+ Amendment 1 (2005)
	+ Amendment 2 (2007)

### IEC CLASSIFICATION





Class II

Type BF Applied Part

### DIMENSIONS AND MASS

Unit	Dimensions W × D × H mm	Mass kg (Approx.)
Wrist coil	152 × 254 × 207	2.6
Coil base	300 × 500 × 20	3.6
Other		1.5

## OUTLINE DRAWINGS

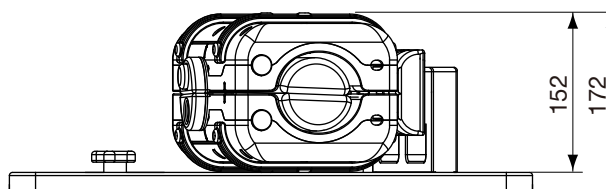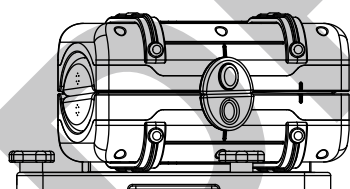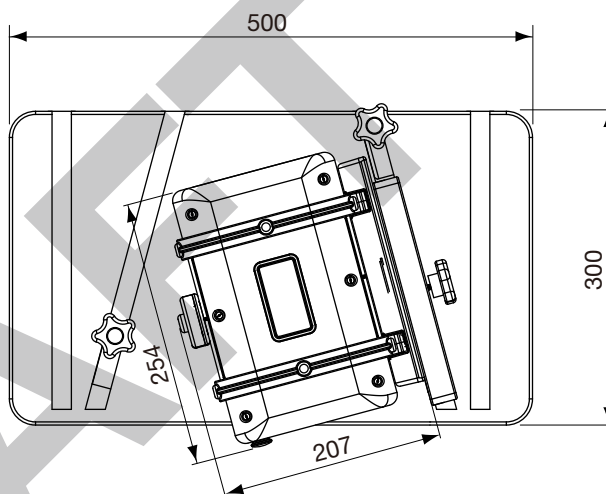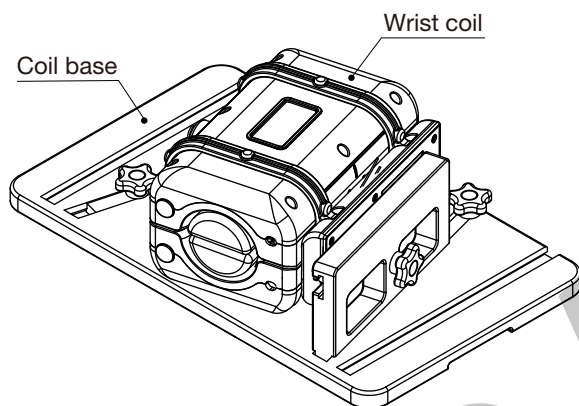

Unit: mm



**TOSHIBA MEDICAL SYSTEMS CORPORATION**

<http://www.toshibamedicalsystems.com>

©Toshiba Medical Systems Corporation 2010 all rights reserved.  
Design and specifications subject to change without notice.  
2010-7 TME/KI

Toshiba Medical Systems Corporation meets internationally recognized standards for Quality Management System ISO 9001, ISO 13485.

Toshiba Medical Systems Corporation Nasu Operations meets the Environmental Management System standard, ISO 14001.

Made for Life and Vantage Titan are trademarks of Toshiba Medical Systems Corporation.

Produced in Japan